



Applications of Environmental TEM for Catalysis Research

Wagner, Jakob Birkedal; Damsgaard, Christian Danvad; Cavalca, Filippo Carlo; Fiordaliso, Elisabetta Maria; Gardini, Diego; Hansen, Thomas Willum

Published in:
Proceedings of CISCEM 2014

Publication date:
2014

Document Version
Publisher's PDF, also known as Version of record

[Link back to DTU Orbit](#)

Citation (APA):
Wagner, J. B., Damsgaard, C. D., Cavalca, F. C., Fiordaliso, E. M., Gardini, D., & Hansen, T. W. (2014). Applications of Environmental TEM for Catalysis Research. In *Proceedings of CISCEM 2014*

General rights

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from the public portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain
- You may freely distribute the URL identifying the publication in the public portal

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.

Applications of Environmental TEM for Catalysis Research

Jakob B. Wagner*, Christian D. Damsgaard, Filippo Cavalca**, Elisabetta M. Fiordaliso, Diego Gardini, Thomas W. Hansen

Center for Electron Nanoscopy, Technical University of Denmark, Fysikvej building 307, DK-2800 Kgs. Lyngby, Denmark

*jakob.wagner@cen.dtu.dk

**present address: Haldor Topsøe A/S, 2800 Kgs. Lyngby, Denmark

In the quest for a better understanding of the dynamics of heterogeneous catalysts under working conditions, *in situ* characterization plays a key role. Catalyst efficiency and sustainability are strongly linked to the dynamic morphology and atomic arrangements of the active entities of these complex materials. Here, we provide examples of how environmental transmission electron microscopy (ETEM) can be applied in catalysis research.

ETEM is a unique tool combining imaging and spectroscopy capabilities with a spatial resolution approaching 1 Å at elevated temperatures and controlled gaseous environments. As with any other characterization technique, ETEM has drawbacks. Limitations in sample geometry and field of view call for use of complementary *in situ* characterization techniques. Such techniques could be *in situ* X-ray based techniques such as XRD, XAS and EXAFS. The combination of techniques will paint a more complete picture of the sample under a reactive environment.

In order to take full advantage of the complementary *in situ* techniques, transfer under reaction conditions is essential. Here we introduce the *in situ* transfer concept by use of a dedicated TEM transfer holder capable of enclosing the sample in a gaseous environment at temperatures up to approx. 900°C. The holder is compatible with other *in situ* technique set-ups [1].

Figure 1 shows the concept of the transfer holder. In this experiment, cuprous oxide has been initially characterized in the ETEM at room temperature in vacuum, then transferred to an *in situ* XRD set-up, reduced in hydrogen, and finally transferred back to the microscope in 10⁵ Pa of H₂ at 220°C.

In general, the catalytic behaviour of nanoparticle systems can be greatly influenced by the synthesis procedure used. The different stages of the formation process such as drying, calcination, and reduction can all be optimized by tweaking parameters such as temperature, time of treatment and gaseous environment. Following the process of identified volume subsets of the catalyst material by means of electron beam based characterization at the different stages of synthesis gives insight in the formation mechanisms of individual nanoparticles from precursor to active catalyst.

Converting solar energy into chemical bonds and thereby increase the storage capability of light energy harvesting requires the use of photocatalysts. In order to characterize and study such catalysts under relevant conditions we have developed a TEM holder capable of exposing a sample to visible light inside the microscope. This way photo induced reactions and phenomena such as Cu₂O degradation (Figure 2) are studied *in situ* under different gaseous environments [2].

[1] C. D. Damsgaard, H. Zandbergen, T. W. Hansen, I. Chorkendorff, J. B. Wagner, *Microsc. and Microanal.* 20 1038 (2014).

[2] F. Cavalca, A. B. Laursen, J. B. Wagner, C. D. Damsgaard, I. Chorkendorff, T. W. Hansen, *ChemCatChem* 5, 2667 (2013).

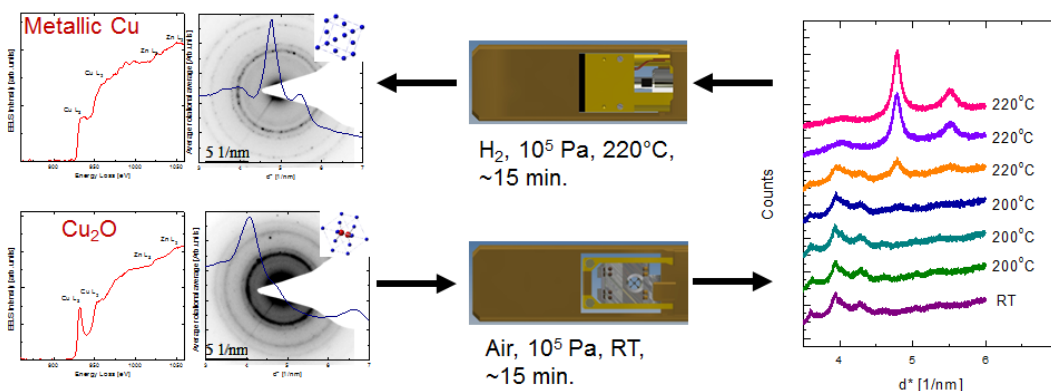


Figure 1: In situ transfer of copper containing nanoparticles between ETEM and in situ XRD setup. The cuprous oxide particles are reduced in the XRD and transferred under elevated temperature in hydrogen to the ETEM keeping the copper metallic.

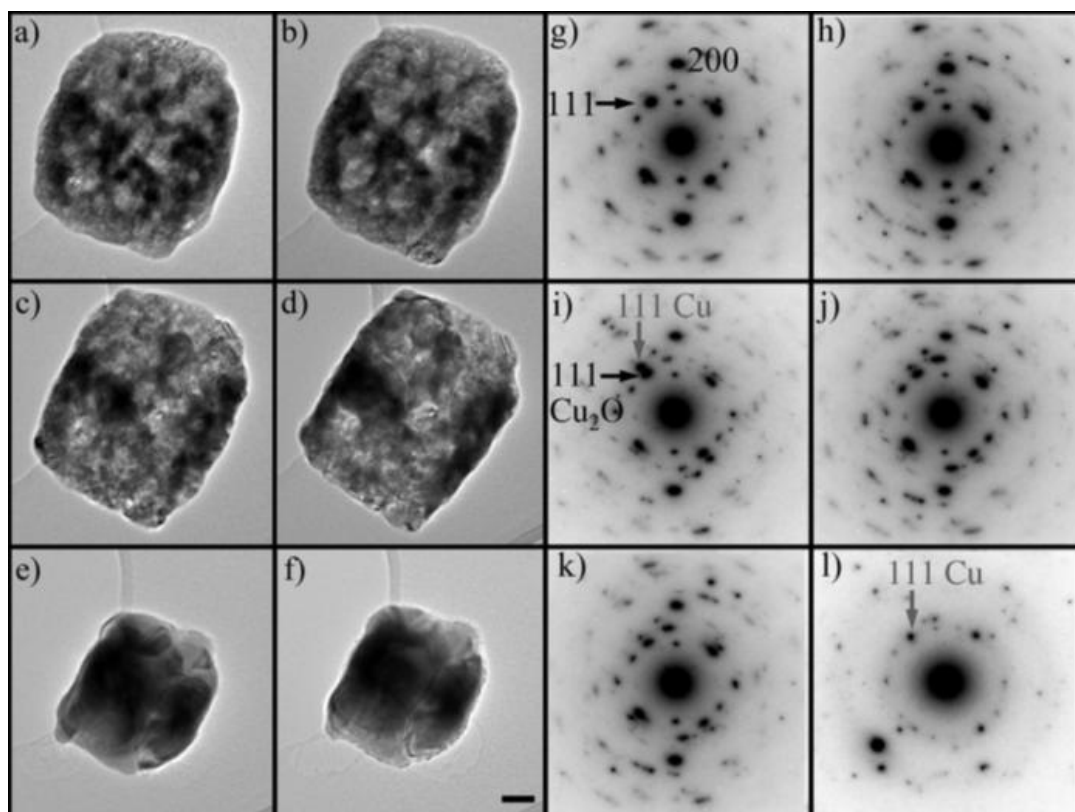


Figure 2: Dynamic study of degradation of Cu_2O to metallic Cu in 500 Pa H_2O under irradiation of visible light in situ in the ETEM. a)-f) TEM imaging. Scale bar 50nm. g)-l) corresponding SAD. Each frame separated by 15 min.