



## The 19th EURL-AR Proficiency Test Salmonella, Campylobacter and genotypic characterisation 2015

Karlslose Pedersen, Susanne; Cavaco, Lina; Hendriksen, Rene S.; Aarestrup, Frank Møller

*Publication date:*  
2016

*Document Version*  
Publisher's PDF, also known as Version of record

[Link back to DTU Orbit](#)

*Citation (APA):*  
Karlslose Pedersen, S., Cavaco, L., Hendriksen, R. S., & Aarestrup, F. M. (2016). *The 19th EURL-AR Proficiency Test Salmonella, Campylobacter and genotypic characterisation 2015*. National Food Institute, Technical University of Denmark. [http://eurl-ar.eu/data/images/reports/salm-camp-gen-report\\_2015\\_isbn\\_inclapp.pdf](http://eurl-ar.eu/data/images/reports/salm-camp-gen-report_2015_isbn_inclapp.pdf)

---


### General rights

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from the public portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain
- You may freely distribute the URL identifying the publication in the public portal

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.

# The 19th EURL-AR Proficiency Test *Salmonella, Campylobacter* and genotypic characterisation 2015

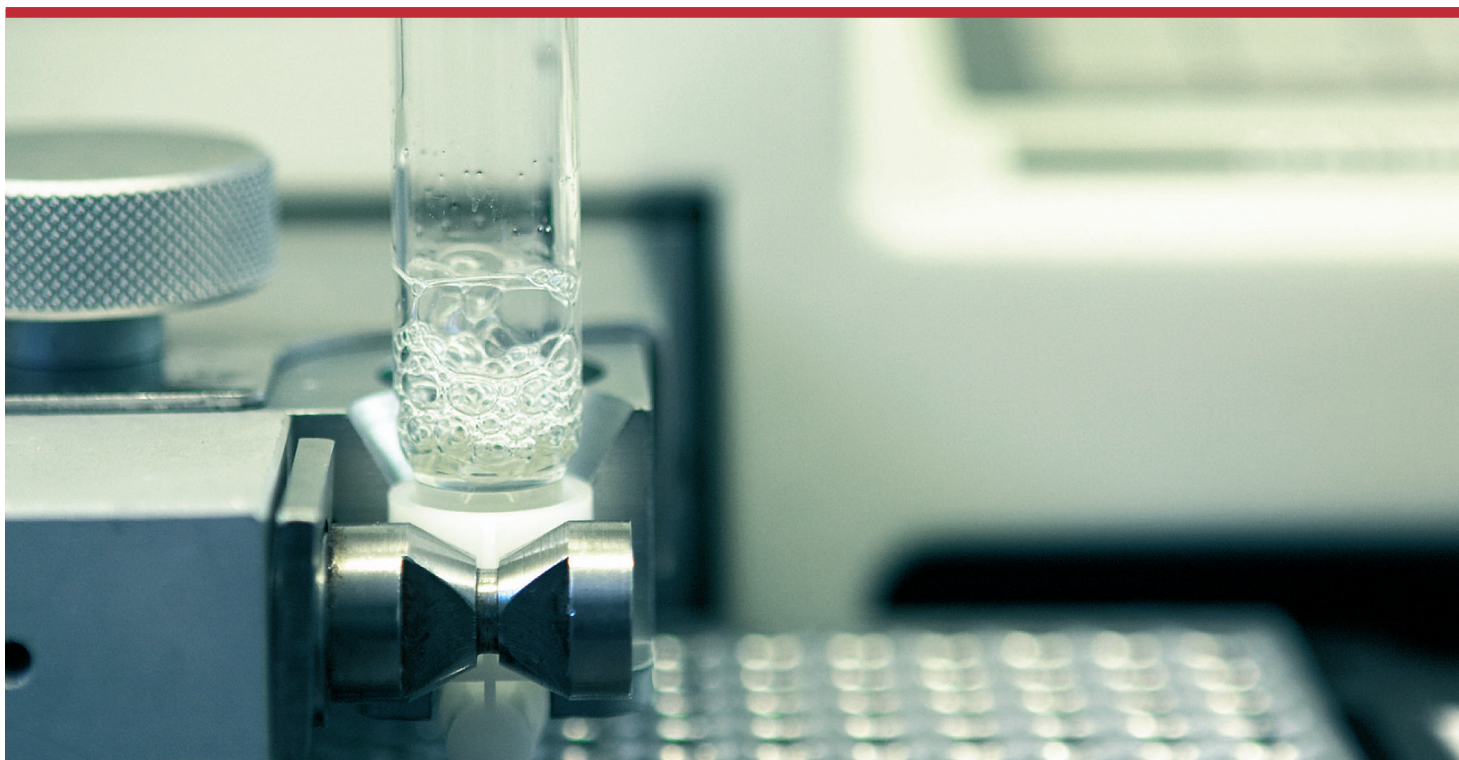

Susanne Karlsmose  
Lina M. Cavaco  
Rene S. Hendriksen  
Frank M. Aarestrup



DTU Food  
National Food Institute

**The 19th EURL-AR Proficiency Test *Salmonella*, *Campylobacter* and genotypic characterisation 2015**

1. edition, June 2016

Copyright: National Food Institute, Technical University of Denmark

Photo: Mikkel Adsbøl

ISBN: 978-87-93109-86-5

The report is available at

[www.food.dtu.dk](http://www.food.dtu.dk)

National Food Institute

Technical University of Denmark

Mørkhøj Bygade 19

2860 Søborg

# Index

1. Introduction.....	4
2. Materials and Methods .....	5
2.1 Participants in EQAS 2015.....	5
2.2 Strains .....	6
2.3 Antimicrobials .....	6
2.4 Distribution .....	7
2.5 Procedure.....	7
3. Results .....	9
3.1 Data omitted from the report.....	9
3.2 Methods .....	9
3.3 Deviations, overall .....	9
3.4 Deviations by laboratory .....	11
3.5 Deviations by reference strains .....	13
3.6 Genotypic characterisation .....	14
4. Discussion .....	14
4.1 <i>Salmonella</i> trial.....	15
4.2 <i>Campylobacter</i> trial .....	15
4.3 Genotypic characterisation .....	16
5. Conclusions .....	16
6. References .....	17
Appendix 1	Pre-notification
Appendix 2	Participant list
Appendix 3a	Reference values (MIC) – <i>Salmonella</i>
Appendix 3b	Reference values (MIC) – <i>Campylobacter</i>
Appendix 4a	Welcome letter
Appendix 4b	Protocol, text
Appendix 4c	Protocol, test forms
Appendix 4d	Instructions for opening and reviving lyophilised cultures
Appendix 4e	Subculture and maintenance of Quality Control strains
Appendix 5	Quality Control ranges for ATCC reference strains
Appendix 6a	Reference strain results – <i>E. coli</i> ATCC 25922
Appendix 6b	Reference strain results – <i>C. jejuni</i> ATCC 33560
Appendix 7a	Expected and obtained results – <i>Salmonella</i>
Appendix 7b	Expected and obtained results – <i>Campylobacter</i>
Appendix 8a	Deviations, <i>Salmonella</i>
Appendix 8b	Deviations, <i>Campylobacter</i>
Appendix 9	Genotypic characterisation (optional)

# 1. Introduction

This report describes and summarises results from the nineteenth proficiency test trial conducted by the National Food Institute (DTU Food) as the EU Reference Laboratory for Antimicrobial Resistance (EURL-AR). This proficiency test focuses on antimicrobial susceptibility testing (AST) of *Salmonella* and *Campylobacter* and is the ninth External Quality Assurance System (EQAS) conducted for these microorganisms (the first was EQAS 2006). In addition, the proficiency test includes categorization of the relevant *Salmonella* strains as presumptive AmpC-, ESBL-, or carbapenemase producing organisms, and identification of the *Campylobacter* species as either *C. jejuni* or *C. coli*,

In addition, for the seventh time, an optional element was included, consisting of genotypic characterization by PCR/sequencing of antimicrobial resistance genes. This optional component included characterization of genes related to production of AmpC, ESBL- and carbapenemases in the *Salmonella* test strains.

This EQAS aims to: i) monitor the quality of AST results produced by National Reference Laboratories (NRL-AR), ii) identify laboratories which may need assistance to improve their performance in AST, and iii) determine possible topics for further research or collaboration.

In reading this report, the following important considerations should be taken into account:

1) Expected results were generated by performing Minimum Inhibitory Concentration (MIC) determinations for all test strains in two different occasions at the Technical University of Denmark, National Food Institute (DTU Food). These results were then verified by the United States Food and Drug Administration (FDA), Centre for Veterinary Medicine. Finally, a fourth MIC determination was performed at DTU Food after preparation of the agar stab culture for shipment to participants to confirm

that the vials contained the correct strains with the expected MIC values.

2) Evaluation is based on interpretations of AST values determined by the participants. This is in agreement with the method used by Member States (MS) to report AST data to the European Food Safety Authority (EFSA), and complies with the main objective of this EQAS, i.e. to evaluate and improve the comparability of surveillance data on antimicrobial susceptibility of *Salmonella* and *Campylobacter* reported to EFSA by different laboratories, as stated in the protocol.

3) The EURL-AR network agreed on setting the accepted deviation level for laboratory performance on AST to 5%. For the optional genotypic characterisation, no specific acceptance level has been set.

Evaluation of a result as “deviating from the expected interpretation” should be carefully analyzed in a self-evaluation procedure performed by the participant including also considerations related to any corrective actions introduced in the laboratory. Note, that since methods used for MIC determination have limitations, it is not considered a mistake to obtain a one-fold dilution difference in the MIC of a specific antimicrobial when testing the same strains. If, however, the expected MIC is close to the breakpoint value for categorizing the strain as susceptible or resistant, a one-fold dilution difference - which is acceptable - may result in two different interpretations, i.e. the same strain can be categorized as susceptible and resistant. This result will be evaluated as correct in one case but incorrect when the evaluation is based on interpretation of MIC values. This report is based on evaluation of AST interpretations, therefore some participants may find their results classified as incorrect even though the actual MIC they reported is only a one-fold dilution different from the

expected MIC. In these cases, the participants should be confident about the good quality of their performance of AST by MIC. In the organization of the EQAS, we try to avoid these situations by choosing test strains with MIC values distant from the breakpoints for resistance, which is not always feasible for all strains and all antimicrobials. Therefore, the EURL-AR network unanimously established in 2008 that if there are less than 75% correct results for a specific strain/antimicrobial combination, the reasons for this situation must be further examined and, on selected occasions explained in details case by case, these results may subsequently be omitted from the evaluation report.

This report is approved in its final version by a technical advisory group composed by

competent representatives from all NRL-ARs. This group meets annually at the EURL-AR workshop.

All conclusions presented in this report are publically available. Participating laboratories are identified by codes and each code is known only by the corresponding laboratory. The full list of laboratory codes is confidential information known only by relevant representatives of the EURL-AR and the EU Commission.

The EURL-AR is accredited by DANAK as provider of proficiency testing (accreditation no. 516); working with zoonotic pathogens and indicator organisms as bacterial isolates (identification, serotyping and antimicrobial susceptibility testing).

## 2. Materials and Methods

### 2.1 Participants in EQAS 2015

A pre-notification (App. 1) to announce the EURL-AR EQAS on AST of *Salmonella* and *Campylobacter* was distributed on the 23<sup>th</sup> June 2015 by e-mail to the 43 NRLs in the EURL-AR-network including all EU countries and Iceland, Norway, Serbia, Switzerland and Turkey. All EU MS as well as Iceland, Norway, and Switzerland were represented as participants for both *Salmonella* and *Campylobacter*. In addition to the AST of *Salmonella* and *Campylobacter*, an optional genotypic characterization by PCR/sequencing of antimicrobial resistance genes of the AmpC-, ESBL- and carbapenemase-producing *Salmonella* test strains was offered.


Appendix 2 shows that 29 of the 34 participating NRLs were appointed by the individual Member States' Competent Authority. Five additional laboratories were included; one from each of the following countries: Iceland, the Netherlands, Norway, Spain, and

Switzerland. These were invited to take part in the EQAS 2015 on the basis of their participation in previous EQAS iterations and/or affiliation to the EU network. These laboratories were charged a fee for their participation in the EQAS, whereas the NRLs from EU Member States participated free of charge.

Figure 1 illustrates that of the 31 participating countries, all tested both *Salmonella* and *Campylobacter*. Nine laboratories participated in the optional genotypic characterisation of the ESC-producing *Salmonella* test strains (not illustrated in Figure 1; see Appendix 2).

The results from the NRLs designated by the MS are presented and evaluated in this report in addition to national reference laboratories in affiliated non-MS; i.e. results from 31 countries consisting of 31 laboratories submitting *Salmonella* results and 31 laboratories submitting *Campylobacter* results. Results from the two laboratories not designated by the MS but enrolled on equal terms as these are not further presented or evaluated in this report.



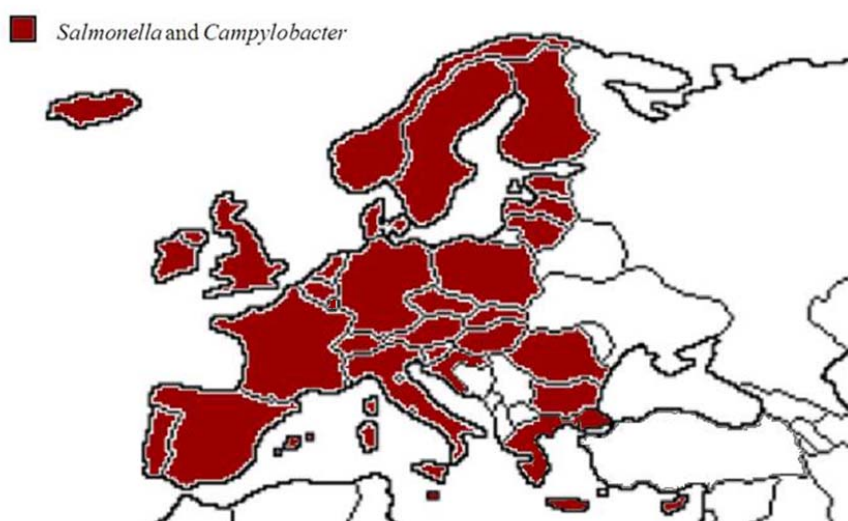

**Figure 1:** Participating countries that performed antimicrobial susceptibility testing of *Salmonella* and *Campylobacter*.

## 2.2 Strains

Eight *Salmonella* strains and eight *Campylobacter* strains were selected for this trial among isolates from the strain collection at DTU Food on the basis of antimicrobial resistance profiles and MIC values. For quality assurance purposes, one strain per bacterial species has been included in all EQAS iterations performed to date, representing an internal control.

Prior to distribution of the strains, AST was performed on the *Salmonella* and *Campylobacter* strains at DTU Food and verified by the US Food and Drug Administration (FDA). When MIC-values were not in agreement but varied +/- one MIC-step, the value obtained by DTU Food was selected as the reference value. The obtained MIC values served as reference for the test strains (App. 3a and 3b). Results from the following antimicrobials were not verified by FDA: cefepime, cefotaxime, cefotaxime/clavulanic acid, ceftazidime, ceftazidime/clavulanic acid, colistin, ertapenem, imipenem, meropenem, temocillin, tigecycline and trimethoprim for *Salmonella* and furthermore, streptomycin for

*Campylobacter*.

Reference strains *Escherichia coli* CCM 3954 (ATCC 25922) and *Campylobacter jejuni* CCM 6214 (ATCC 33560) were provided to new participating laboratories with instructions to store and maintain them for quality assurance purposes and future EQAS trials.

## 2.3 Antimicrobials

The antimicrobials tested in this EQAS are listed in the protocol (App. 4b).

The antimicrobials tested correspond to the panel of antimicrobials listed in Decision 2013/652/EU.

Guidelines for performing AST were set according to the Clinical and Laboratory Standards Institute (CLSI) document; M7-A10 (2015) "Methods for Dilution Antimicrobial Susceptibility Tests for Bacteria That Grow Aerobically; Approved Standard - Tenth Edition"; M100-S25 (2015) "Performance Standards for Antimicrobial Susceptibility Testing" (Twenty-Fifth Informational Supplement) and document VET01-A4 (2013) "Performance Standards for Antimicrobial Disk and Dilution Susceptibility Tests for Bacterial



Isolated From Animals” (Approved Standard – Fourth Edition).

MIC results were interpreted by using the interpretative criteria listed in Decision 2013/652/EU. Where values were not available, the list of interpretative criteria was supplemented with CLSI-interpretative criteria as described and indicated in the protocol (App. 4). No interpretative criteria were available to determine the interpretation of MIC-values from testing of azithromycin, cefepime and temocillin. Results of ESC detection tests were interpreted according to the at deadline most recent EFSA recommendations (EFSA Journal 2012; 10(6):2742).

The selection of antimicrobials used in the trial for *Salmonella* were: ampicillin (AMP), azithromycin (AZI), cefepime (FEP), cefotaxime (FOT), cefotaxime/clavulanic acid (FOT/CI), ceftazidime (TAZ), ceftazidime/clavulanic acid (TAZ/CI), chloramphenicol (CHL), ciprofloxacin (CIP), colistin (COL), ertapenem (ERT), gentamicin (GEN), imipenem (IMI), meropenem (MER), nalidixic acid (NAL), sulfonamides (sulfamethoxazole) (SMX), tetracycline (TET), tigecycline (TGC), temocillin (TRM) and trimethoprim (TMP).

Minimum Inhibitory Concentration (MIC) determination of the *Salmonella* test strains was performed using the Sensititre system from Trek Diagnostic Systems Ltd, UK. For ESC confirmatory test, the analysis included MIC determination by microbroth dilution.

For *Campylobacter* the following antimicrobials were included: ciprofloxacin (CIP), erythromycin (ERY), gentamicin (GEN), nalidixic acid (NAL), streptomycin (STR), and tetracycline (TET). MIC determination for the *Campylobacter* testing, was performed using the Sensititre systems from Trek Diagnostic Systems Ltd, UK, according to guidelines from the CLSI document M45-A2 (2010) “Methods for Antimicrobial Dilution and Disk Susceptibility

Testing of Infrequently Isolated or Fastidious Bacteria” (Approved Guideline – Second Edition) and VET01-S2 (2013) “Performance Standards for Antimicrobial Disk and Dilution Susceptibility Tests for Bacteria Isolated From Animals” (Second Informational Supplement). Participants of the *Campylobacter* EQAS were additionally requested to identify the species of the *Campylobacter* spp. as either *C. jejuni* or *C. coli*.

## 2.4 Distribution

On the 14<sup>th</sup> October 2015, bacterial strains in agar stab cultures (*Salmonella* spp.) or charcoal swabs in transport media (Stuarts) (*Campylobacter* spp.) together with a welcome letter (App. 4a) were dispatched in double pack containers (class UN 6.2) to the participating laboratories. The shipment (UN3373, biological substances category B) was sent according to International Air Transport Association (IATA) regulations.

## 2.5 Procedure

Protocols and all relevant information were uploaded on the EURL-AR website (<http://www.eurl-ar.eu>), thereby EQAS participants could access necessary information at any time.

Participants were instructed to subculture charcoal swabs immediately, store the agar stabs at 4°C (dark) and the freeze-dried strains cool and dark until performance of AST. Information related to the handling of the test strains and reference strains (App. 4b, c, d, e) was made available. Participants receiving an ATCC reference strain were requested to save and maintain this strain for future proficiency tests.

The participants were instructed to apply the interpretative criteria listed in the protocol (App. 4). Instructions for interpretation of AST results allowed for categorization of results as resistant or susceptible. Categorisations as ‘intermediate’ were not accepted.



The EURL-AR is aware that there are two different types of interpretative criteria of results, clinical breakpoints and epidemiological cut-off values. The terms 'susceptible', 'intermediate' and 'resistant' should be reserved for classifications made in relation to the therapeutic application of antimicrobial agents. When reporting data using epidemiological cut-off values, bacteria should be reported as 'wild-type' or 'non-wild-type' (Schwarz et al., 2010). Due to the different methods of AST used by the participants and also to simplify the interpretation of results, throughout this report, we will still maintain the terms susceptible and resistant, even in cases where we are referring to wild-type and non-wild-type strains.

As regards the method for performing the antimicrobial susceptibility testing, the protocol referred to Decision 2013/652/EU and instructed participants to perform a dilution method, i.e. microbroth dilution or agar dilution. Results obtained by methods not complying with the description in Decision 2013/652/EU were not included in the present analysis.

A mandatory part of the proficiency test was to detect ESC-producing strains and interpret results according to the most recent EFSA recommendations (EFSA Journal 2012; 10(6):2742) as described in the protocol.

Results from QC reference strains would consist of MIC values for the reference strains *E. coli* (ATCC 25922) and *C. jejuni* (ATCC 33560). The results were evaluated towards the quality control ranges according to the relevant guidelines; i.e. the CLSI documents VET01-S2 (2013) or M100-S25 (2015) (App. 5).

For the optional genotypic characterisation of the ESC-producing *Salmonella* test strains, participating laboratories were requested to report the genes conferring resistance to extended-spectrum beta lactam antimicrobials. The organizers, however, decided to include none-ESC TEM-genes resulting in *bla*<sub>TEM-1</sub> registered as an expected gene, also. The

genes listed in the table in the protocol (App. 4b) were included in the test. Identification of additional genes not listed in the protocol was not evaluated by the database. The results were evaluated based on the actual genes and variants identified.

The participating laboratories were encouraged to use their own laboratory's method(s) for the genotypic characterisation. The expected results for this component of the EQAS were obtained by whole-genome-sequencing and subsequent analysis using the ResFinder 2.1 platform available at <http://cge.cbs.dtu.dk/services/ResFinder/>. The positive identification of genes was not verified elsewhere.

All participating laboratories were invited to enter the obtained results into an electronic record sheet at the EURL-AR web-based database through a secured individual login and password. The record sheet contained space for reporting the results obtained for the QC reference strains. Alternatively, it was offered the possibility to fill-in a record sheet (provided with the protocol) and to send it to the EURL-AR by fax, mail or email.

In addition, participants were encouraged to complete an evaluation form available at the EURL-AR database with the aim to improve future EQAS trials.

The database was finally closed and evaluations were made available to participants on the 9<sup>th</sup> December 2015. After this date, the participants were invited to login to retrieve an individual, database-generated report which contained an evaluation of the submitted results including possible deviations from the expected interpretations. Deviations in the interpretation as resistant or susceptible were categorised as 'incorrect', as were also deviations concerning confirmation of an isolate as extended spectrum beta-lactamase- (ESBL-), ampC- or carbapenemase-producer.

## 3. Results

The participants were asked to report results, i.e. MIC values and the categorisation as resistant or susceptible. Only the categorisation was evaluated, whereas the MIC values were used as supplementary information.

### 3.1 Data omitted from the report

As mentioned in the introduction, the EURL-AR network established that data should be examined and possibly omitted from the general analysis if there are less than 75% correct results based on strain/antimicrobial combination (see Appendix 8 for an overview of correct/incorrect results). In the present EQAS this occurred in two cases which have been examined and consequently omitted from the analysis; 1) S10.4/imipenem (expected interpretation was 'resistant', however, 30% (9 laboratories) found the strain susceptible to imipenem. All but one of the deviating interpretations were based on MIC values two steps from the expected, one was three steps from the expected; 2) C10.4/tetracycline (expected interpretation was 'susceptible', however, 26% of participants found the strain 'resistant' to tetracycline. All deviating interpretations were based on MIC values only one step from the expected. Both these combinations were subsequently omitted from further analysis.

### 3.2 Methods

The agar dilution method and MIC determination were evaluated together as they are both quantitative methods giving results corresponding to the MIC of the bacterial strain tested.

In the *Salmonella* as well as the *Campylobacter* trial, 30 laboratories performed microbroth dilution and one performed agar dilution.

With the aim to conclude on the strain's presumptive ESBL, AmpC or carbapenemase

phenotype, two panels of antimicrobials were included in the testing of the *Salmonella* strains. The test strains found resistant to cefotaxime, ceftazidime or meropenem on the first panel were additionally tested on the second panel according to the protocol indications.

### 3.3 Deviations, overall

The list of deviations is shown in Appendix 8a and 8b. Figure 2 shows the total percentage of deviations from the expected results of AST performed by participating laboratories. The internal control strains mainly followed the trend in deviation level of the different EQAS trials (Figure 2). The deviation level in 2015 is acceptable for both the *Salmonella* and the *Campylobacter* trials. For the current EQAS, for both microorganisms, it appears that there has been a decrease in the level of deviations, to 0.7% for *Salmonella* and 1.6% for *Campylobacter* in 2015 from levels at 2.4% and 4.0%, respectively in 2014.


#### 3.3.1 *Salmonella* trial

For the *Salmonella* strains, 99.3% of the AST's were interpreted correctly. The number of AST's performed and the percentage of correct results for the individual strains in the EQAS, are listed in Table 1. Variations of obtained correct results ranged from 98.5-99.8% between the *Salmonella* strains. Table 2 illustrates the percentage of correct AST per antimicrobial by bacterial species. The level of correct AST was at 97.8% (tigecycline) or above, for all the *Salmonella* test strains.

#### ESC-producing *Salmonella* test strains

Confirmation of beta-lactamase production is a mandatory component of this EQAS.

According to the protocol, which was based on the EFSA recommendations, the confirmatory test for ESC-production requires use of both cefotaxime (FOT) and ceftazidime (TAZ) alone

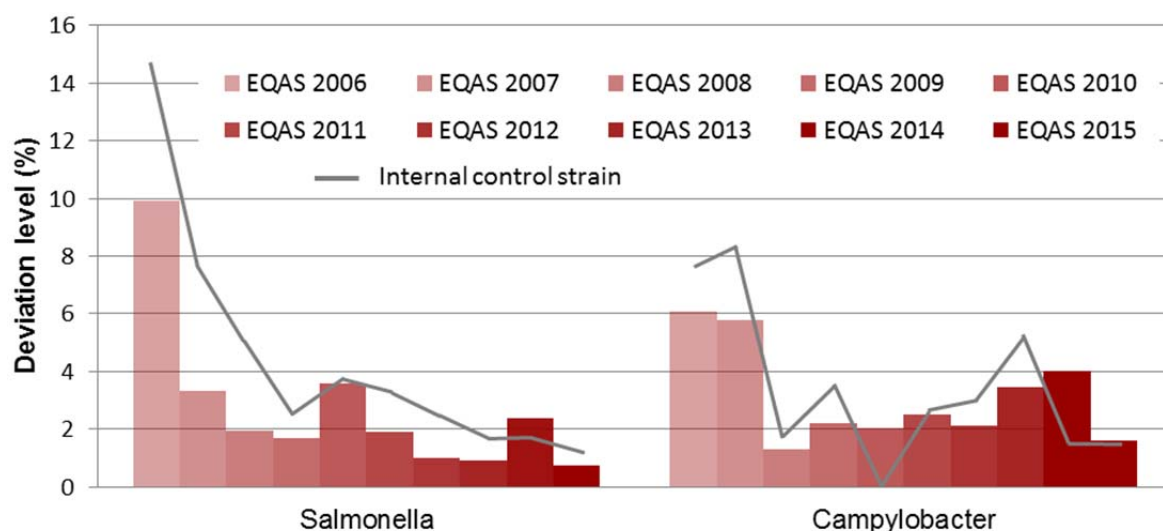

**Figure 2:** A comparison between the EURL-AR EQAS's since 2006, showing the total percentage of deviations for antimicrobial susceptibility testing performed by participating laboratories.

and in combination with a  $\beta$ -lactamase inhibitor. The MIC value for either antimicrobial agent (FOT or TAZ) tested in combination with clavulanic acid should be compared to the corresponding MIC when tested alone. Synergy is indicated if a three dilution steps difference is observed between the two MIC values (i.e. if the FOT:CTX/CI or TAZ:TAZ/CI ratio  $\geq 8$ ) (CLSI M100 Table 2A; Enterobacteriaceae). Participants were instructed to test strains presenting resistance to cefotaxime (FOT), ceftazidime (TAZ or meropenem (MERO) on the second panel of antimicrobials.

The classification of the phenotypic results was based on the most recent EFSA recommendations (EFSA 2012), indicating:

- Presumptive ESBL-phenotype: strains with positive synergy test, susceptible to ceftazidime and resistant to ceftazidime
- Presumptive ESBL+pAmpC-phenotype: strains with positive or negative synergy test, resistant to ceftazidime and resistant to ceftazidime

- Presumptive pAmpC phenotype: strains with negative synergy test resistant to ceftazidime and susceptible to ceftazidime
- Presumptive carbapenemase phenotype: strain resistant to meropenem
- Unusual phenotype: any other combinations

In this EQAS, all laboratories uploaded results for the strains exhibiting resistance to the cephalosporins tested.

The strains S-10.3, S-10.7 and S-10.8 were ESBL-producers, and S-10.4 was a carbapenemase-producer. Note that when categorizing the presumptive phenotypes, the interpretation of the ceftazidime result had to be disregarded, as no interpretative criteria were available.

In total, the categorization as ESBL-, pAmpC- or carbapenemase-producer was correct in almost all cases; i.e. out of 248 reported results, three were incorrect. The deviating results were presented by three different laboratories (#39, #56, and #60). For laboratory #39, the reported phenotypic results do not explain the selected presumptive phenotype for S-10.4 as 'unusual

**Table 1.** The number of AST performed and the percentage of correct results for each strain of *Salmonella* and *Campylobacter*.

EQAS 2015 – <i>Salmonella</i>			EQAS 2015 – <i>Campylobacter</i>		
Test strain	AST in total	% correct	Test strain	AST in total	% correct
S-10.1	400	99.3	C-10.1 ( <i>C. jejuni</i> )	185	97.3
S-10.2	400	99.8	C-10.2 ( <i>C. coli</i> )	186	98.4
S-10.3	587	98.8	C-10.3 ( <i>C. jejuni</i> )	186	97.3
S-10.4	552	99.6	C-10.4 ( <i>C. coli</i> )	155	100.0
S-10.5	398	98.5	C-10.5 ( <i>C. jejuni</i> )	186	98.4
S-10.6	399	99.5	C-10.6 ( <i>C. coli</i> )	185	98.4
S-10.7	585	99.1	C-10.7 ( <i>C. coli</i> )	186	97.8
S-10.8	584	99.7	C-10.8 ( <i>C. jejuni</i> )	186	100.0

**Table 2:** Percentage of correct antimicrobial susceptibility tests per antimicrobial by microorganism.

Antimicrobial	<i>Salmonella</i>	<i>Campylobacter</i>
Ampicillin	99.2	-
Cefotaxime	100.0	-
Cefoxitin	99.2	-
Ceftazidime	99.5	-
Chloramphenicol	99.6	-
Ciprofloxacin	98.8	99.6
Colistin	98.0	-
Ertapenem	100.0	-
Erythromycin	-	99.6
Gentamicin	99.6	97.6
Imipenem	100.0	-
Meropenem	100.0	-
Nalidixic acid	99.6	98.0
Streptomycin	-	97.2
Sulphonamides	98.4	-
Tetracycline	99.6	98.6
Tigecycline	97.8	-
Trimethoprim	99.2	-

phenotype'. As follow-up, the laboratory reported that the cause for the incorrect categorization was based on the laboratory SOP not being strictly followed (the SOP clearly states that resistance to meropenem leads to 'presumptive carbapenemase phenotype'). For laboratory #56 and #60, however, the deviating results were caused by an incorrect result as resistant for cefoxitin (#56; S-10.7), and laboratory #60 commented that the positive

synergy test and susceptibility to cefoxitin for strain S-10.3 indicated presumptive ESBL. The laboratory stated that the cefepime MIC value is a little lower than they usually see for this profile, and together with the fact that the strain is ceftazidime sensitive, they concluded that the strain was 'unusual phenotype'.

### 3.3.2 *Campylobacter* trial

For the *Campylobacter* strains, 98.4% of AST's were correctly tested. Table 1 presents that the variation in the obtained correct results ranged from 97.3-100% and Table 2 illustrates that the percentage of correct AST per antimicrobial was above 97.2% for the *Campylobacter* test strains with streptomycin exhibiting the lowest level.


The participants were requested to identify the *Campylobacter* species. All 31 laboratories delivered in total 248 results of which one identification was incorrect (laboratory #34). The incorrect species identification did not lead to incorrect AST results for this laboratory and strain.

## 3.4 Deviations by laboratory


Figure 3 and 4 illustrate the percentage of deviations for each participating laboratory. The laboratories are ranked according to their performance determined by the percentage of deviating results in the antimicrobial susceptibility tests.

**Table 3:** Overview of ESBL-, pAmpC- and carbapenemase-producing *Salmonella* test strains and proportion of laboratories that obtained the expected result; number and percentages of laboratories which correctly detected and confirmed the ESBL-, pAmpC- and carbapenemase-producing *Salmonella* strains. Fields shaded in grey with numbers in *italics* indicate an unexpected result.

		Strain S-10.3	Strain S-10.4	Strain S-10.7	Strain S-10.8
ESC-genes harboured in the test strain		<i>bla</i> <sub>TEM-1</sub> <i>bla</i> <sub>CTX-M-9</sub>	<i>bla</i> <sub>TEM-1</sub> <i>bla</i> <sub>OXA-48</sub>	<i>bla</i> <sub>TEM-1</sub> <i>bla</i> <sub>SHV-12</sub> <i>bla</i> <sub>CTX-M-15</sub>	<i>bla</i> <sub>TEM-52</sub>
ESBL-, pAmpC- and carbapenemase-producing strain – expected results		ESBL	carbapenemase	ESBL	ESBL
Obtained results	Confirmed ESBL-producer	30/31 (97%)	-	30/31 (97%)	31/31 (100%)
	Confirmed ESBL + pAmpC-producer	-	-	<i>1/31 (3%)</i>	-
	Confirmed pAmpC-producer	-	-	-	-
	Confirmed carbapenemase-producer	-	30/31 (97%)	-	-
	Confirmed unusual phenotype	<i>1/31 (3%)</i>	<i>1/31 (3%)</i>	-	-
	Not ESBL-, pAmpC- or carbapenemase-producing	-	-	-	-

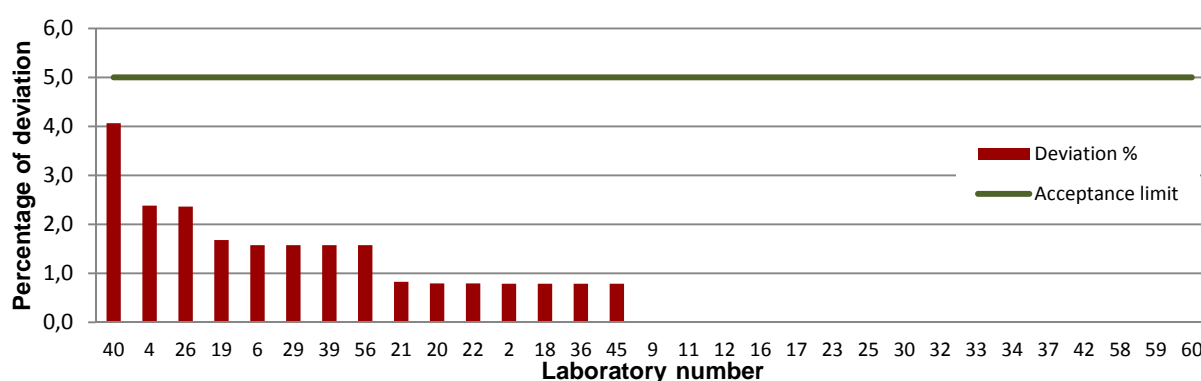

**Figure 3:** Individual participants' deviations in percent of their total number of *Salmonella* AST's.

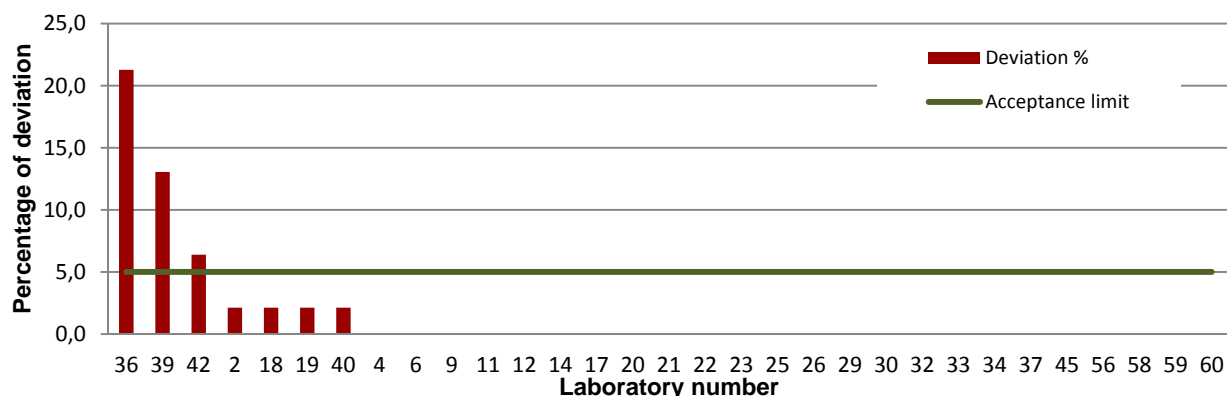

**Figure 4:** Individual participants' deviations in percent of their total number of *Campylobacter* AST's.





**Table 4** Obtained values for AST of *E. coli* ATCC 25922. AMP; ampicillin, FEP; cefepime FOT; cefotaxime, FOX; cefoxitin, TAZ; ceftazidime, CHL; chloramphenicol, CIP; ciprofloxacin, COL; colistin, ERT: ertapenem, GEN; gentamicin, IMI; imipenem, MER; meropenem, NAL; nalidixic acid, SMX; sulphonamides, TET; tetracycline, TGC; tigecycline, TMP; trimethoprim.

MIC determination <i>E. coli</i> ATCC 25922			
Antimicrobial	Proportion outside QC range	Obtained values in MIC steps (min/max)	
		Below lower QC limit	Above upper QC limit
Panel 1, AMP	0/29 (0%)	-	-
Panel 1, FOT	0/27 (0%)	-	-
Panel 1, TAZ	0/28 (0%)	-	-
Panel 1, CHL	0/29 (0%)	-	-
Panel 1, CIP	1/29 (3%)	-	1 step
Panel 1, COL	0/29 (0%)	-	-
Panel 1, GEN	0/29 (0%)	-	-
Panel 1, MER	0/28 (0%)	-	-
Panel 1, NAL	0/29 (0%)	-	-
Panel 1, SMX	1/28 (4%)	-	1 step
Panel 1, TET	0/29 (0%)	-	-
Panel 1, TGC	0/28 (0%)	-	-
Panel 1, TMP	3/29 (10%)	1 step	-
Panel 2, FEP	1/24 (4%)	-	2 steps
Panel 2, FOT	1/24 (4%)	-	1 steps
Panel 2, FOX	0/25 (0%)	-	-
Panel 2, TAZ	0/25 (0%)	-	-
Panel 2, ERT	0/25 (0%)	-	-
Panel 2, IMI	0/25 (0%)	-	-
Panel 2, MER	0/25 (0%)	-	-

**Table 5** Obtained values for AST of *C. jejuni* ATCC 33560. CIP; ciprofloxacin, ERY; erythromycin, GEN; gentamicin, NAL; nalidixic acid, TET; tetracycline.

MIC determination <i>C. jejuni</i> ATCC 33560			
Antimicrobial	Proportion outside QC range	Obtained values in MIC steps (min/max)	
		Below lower QC limit	Above upper QC limit
CIP	1/29 (3%)	-	1 step
ERY	0/29 (0%)	-	-
GEN	3/28 (11%)	2 steps	-
NAL	0/28 (0%)	-	-
TET	2/28 (7%)	1 step	2 steps

### 3.4.1 *Salmonella* trial

All 31 participating laboratories obtained a result within the acceptance limit at 5% deviations for the *Salmonella* strains. The maximum percentage of deviations was at 4.1%, presenting a very good result across the EURL-AR network.

### 3.4.2 *Campylobacter* trial

In the *Campylobacter* trial, most laboratories performed very well. Applying the 5% acceptance threshold, 28 of 31 participating laboratories performed acceptably, with 24 laboratories having no deviations (Figure 4).

Three laboratories present a deviation level above the 5% acceptance level (#36, #39, and #42). Of these, one laboratory with a deviation level at 21.3% (#36) was regarded as an outlier.


## 3.5 Deviations by reference strains

In the following section, deviations are defined as results of antimicrobial susceptibility tests on the reference strain that are outside the quality control (QC) acceptance intervals (App. 5).

Obtained values from the participants' testing of the QC strains are listed in Appendix 6a and 6b, and in Table 4 and 5. For both the *Salmonella* and *Campylobacter* trial, 29 laboratories uploaded data from QC-testing on the relevant reference strain.

Appendix 6a indicates that of the 29 laboratories submitting AST-results for the reference strain *E. coli* ATCC 25922, six laboratories produced in all seven values outside the QC-limit. Table 4 illustrates the obtained results which are fully presented in Appendix 6a.

Table 5 presents the proportion of laboratories with results for the *C. jejuni* reference strain ATCC 33560 below or above the QC interval. six deviations were seen from five different laboratories.

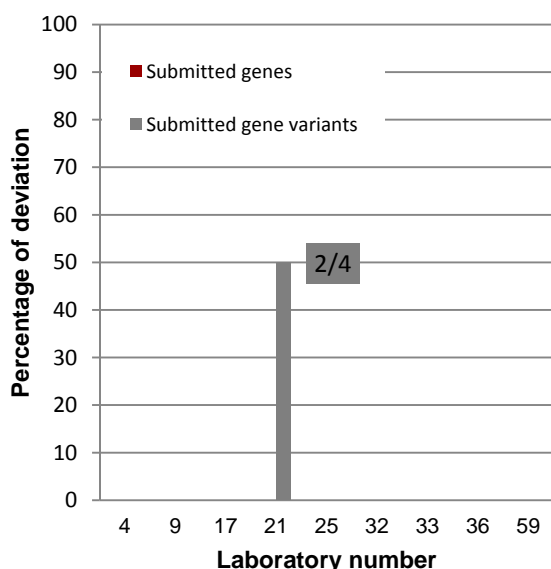

**Figure 5:** Individual participants' deviations in percent of their total number of results from the genotypic characterization.

### 3.6 Genotypic characterisation

For the optional genotypic characterisation of the ESC-producing *Salmonella* test strains, nine laboratories participated. In Appendix 9, information is collected on detected genes, genes which were tested but not detected, primers used, and references for the method used. One laboratory performed whole genome sequencing of the ESC-producing *Salmonella* whereas the remaining eight laboratories indicated the use of various types of conventional PCR to identify the relevant genes.

Table 6 indicate the obtained results, both on gene and variant level. Moreover, Figure 5 indicates that two discordant results related to the gene variant were submitted by laboratory #21. CTX-M-15 belongs to the CTX M-1-group, however, it is not clear if this was the background for the deviation.

**Table 6:** Results from the participation of nine laboratories in the optional genotypic characterisation component of the EQAS

Test strain	Expected gene	Proportion of correct results (gene level)	Proportion of correct results (variant level)	Additional genes/variants identified
S-10.3	TEM-1	7/7 (100%)	5/5 (100%)	CTX M-4
	CTX-M-9	9/9 (100%)	7/8 (88%)	
S-10.4	TEM-1	6/6 (100%)	4/4 (100%)	
	OXA-48	9/9 (100%)	9/9 (100%)	
S-10.7	TEM-1	8/8 (100%)	6/6 (100%)	CTX M-1
	SHV-12	8/8 (100%)	6/6 (100%)	
	CTX-M-15	9/9 (100%)	7/8 (88%)	
S-10.8	TEM-52	9/9 (100%)	8/8 (100%)	

## 4. Discussion

It is important to consider that the number of EQAS participants differs from year to year, which implies that comparisons among different EQAS iterations should be interpreted with caution.

As also specified in the EU regulation 2013/652/EU, all participants in the present

EQAS performed AST by dilution methods, primarily as microbroth determination.

This 2015 proficiency test is the second possibility of testing *Salmonella* and *Campylobacter* strains with the panels designed to follow the requirements of Decision 2013/652/EU. This allows for the possibility that

the experience obtained since the introduction of the legislation and the focus it has created on AST in the laboratories has had an impact on the generally satisfactory results obtained at the present EQAS.

#### 4.1 *Salmonella* trial

Overall, the percentage of correct antimicrobial susceptibility test results of *Salmonella* was 99.3%. All (n=31) participants obtained satisfactory results according to the level of acceptance (<5% deviation).

As indicated in Figure 2, the overall quality of the results in the 2015-EQAS would appear to have increased again to the level of the years 2012 and 2013, also, the measure when comparing results obtained from testing the internal control strain indicates a steady and very good quality of results.

As indicated by Figure 3, all laboratories exhibited very good results with deviation levels below 5%. Follow-up has therefore not been necessary based on these results, and none of the laboratories were defined as outliers.

For the *E. coli* reference strain, the obtained results were in general in agreement with the CLSI recommendations. Two laboratories (#4 and #59) did not submit values related to the quality control strain, but within the submitted results, trimethoprim appeared to have most results outside the QC-range, and for cefepime one result was two steps above the QC-range.

For the two laboratories #57 and #58 which had a deviation level for the AST results above the acceptance limit in EQAS 2014 with values of 7.9% and 9.4%, respectively, one (#58) has this year increased their performance and have no deviations this year, and the other (#57) did not participate in the 2015-iteration.

##### ESC-producing *Salmonella* test strains

The detection of ESC-producing microorganisms remains to be important and is a mandatory part of this EQAS.

Of the four *Salmonella* test strains relevant for this component of the EQAS (S-10.3, S-10.4, S-10.7, and S-10.8), three were ESBL-phenotypes and one was a carbapenemase phenotype. The testing and interpretation of results for these strains appeared not to cause major difficulties for any of the participating laboratories.

Of the 31 laboratories which tested *Salmonella*, three laboratories (#39, #56 and #60) each submitted one incorrect AmpC-, ESBL-, carbapenemase categorization (App. 8a). The deviations all had background in the laboratory handling of the strains; for one, the incorrect categorization was based on an incorrect result as resistant for cefoxitin (#56; S-10.7). Even if no acceptance limit has been defined for this component of the EQAS, the overall result that 98.8% of the obtained results were as expected appears satisfactory.

#### 4.2 *Campylobacter* trial

For the *Campylobacter* component of this year's EQAS, 31 laboratories submitted results leading to an overall percentage of correct AST results at 98.4%. The performance varied from no deviations up to 21.3% deviations, with 28 laboratories performing satisfactorily according to the established acceptance range.

It appears that there has been a decrease in the level of deviations for the overall AST result. Also, results obtained from testing the internal control strain indicate a steady and very good quality of results.

Three laboratories (#36, #39, and #42) obtained deviation levels above 5%, one of these was defined as an outlier (#36) with a deviation level at 21.3%. For none of these laboratories, the values obtained for the QC-strain did not indicate methodical issues to be the reason for the obtained deviations.

The EURL-AR have been in contact with all three laboratories presenting deviation levels above 5% to identify possible causes of this

unsatisfactory performance and to improve the quality of results. As for laboratory #36, the high level of deviations was caused by the switch of panels used and not switching to the corresponding result sheet for reading the plates. The laboratory transferred the obtained results to the correct result sheet and subsequently found that not all obtained results were as expected. As for laboratory #39, re-testing of strain C-10.3 that caused five deviations was performed, leading this time to one deviation, only. The laboratory will continue follow-up when fresh stocks of the recommended panel for *Campylobacter* are received. At laboratory #42, the obtained MICs were as expected, but the deviations were caused by typos when submitting the interpretation into the database.

All participating laboratories except two (#4 and #59) uploaded data from tests performed on the *C. jejuni* reference strain and the proportion of results within the QC intervals was 95.8%. Four of the six values outside the QC intervals were one step below or above the QC-limits, the remaining two were two dilution steps above or below the QC-limits. The laboratories obtaining these values should follow-up on these high/low values, and it is suggested that these values are monitored over time to ensure that the tests render a reliable result for the particular antimicrobial.

Laboratories #29, and #40 which were regarded

as outliers for the 2014 *Campylobacter* EQAS with deviation levels at 22.9% and 34.0%, respectively, both increased their performance extensively in the 2015-iteration and obtained deviation levels at 0.0%, and 2.1%, respectively.

### 4.3 Genotypic characterisation

The focus on genotypic characterization of microorganisms is increasing in the EU and worldwide. In EU, communication has been ongoing to improve laboratory detection and confirmation of ESBL- and AmpC-producing *Enterobacteriaceae*.

Furthermore, the agenda now is focusing at the implementation of detection of carbapenemase resistant organisms and the importance of determining the identity of the genes responsible for the carbapenemase production by molecular methods.

The optional genotypic characterisation offered as a supplementary part of this EQAS should therefore be seen as an important possibility for the NRL-AR's to introduce this method in the laboratory and thereby be at the forefront when the method proposals are adopted. This year, nine laboratories participated in this optional EQAS component and even if no acceptance limit has been defined, the 98.3% correct results (N=119) appears to be a satisfactory results.

## 5. Conclusions

The goal of the EURL-AR EQAS is to have all participating NRLs performing antimicrobial susceptibility testing of *Salmonella* and *Campylobacter* with a deviation level below 5%. This seems within reach for *Salmonella* as well as for *Campylobacter*.

Compared to the EQAS 2014, the performance of the NRL's in 2015 appears to have improved for *Salmonella* AST's (99.3% in 2015 and

97.6% in 2014) and is now again at the level from 2013 and before (Figure 2). Regarding *Campylobacter* AST's, the performance of the NRL's also appear to have improved from 2014 to 2015, with a change in deviation level from 4.0% (2014) to 1.6% (2015). For the *Campylobacter* AST, one laboratory (#36) was regarded as an outlier. Follow-up/re-testing internally at the NRL's with deviation levels above 5% has shown acceptable results for all

three laboratories (#36, #39, and #42).

The test covering the identification of the phenotype of *Salmonella* test strains producing beta-lactamases of the ESBL-, AmpC, and carbapenemase-type rendered acceptable results. This is a priority area within the EURL-AR activities, and it is encouraging to see acceptable results in identifying and categorizing these strains.

Nine NRLs participated in the EQAS

component consisting of genotypic testing of ESBL-, AmpC- and carbapenemase-producing *Enterobacteriaceae* presenting satisfactory results.

Finally, the EURL-AR is open to suggestions to improve future EQAS trials and invites the entire network to contribute with ideas for training courses and specific focus areas to expand the network's knowledge in antimicrobial resistance.

## 6. References

---

**EFSA**, Technical specifications on the harmonised monitoring and reporting of antimicrobial resistance in *Salmonella*, *Campylobacter* and indicator *Escherichia coli* and *Enterococcus* spp. bacteria transmitted through food. EFSA Journal 2012;10(6):2742 [64 pp.].

**European Commission**, 2013/652/EU: Commission Implementing Decision of 12

November 2013 on the monitoring and reporting of antimicrobial resistance in zoonotic and commensal bacteria

**Schwarz S**, Silley P, Simjee S, Woodford N, van DE, Johnson AP & Gaastra W. (2010) Editorial: assessing the antimicrobial susceptibility of bacteria obtained from animals. J Antimicrob Chemother 65: 601-604



G00-06-001/01.12.2014

## EURL-AR EQAS pre-notification

### **EQAS 2015 FOR *SALMONELLA*, *CAMPYLOBACTER* AND OPTIONAL GENOTYPIC CHARACTERISATION**

The EURL-AR announces the launch of another EQAS, thus providing the opportunity for proficiency testing which is considered an essential tool for the generation of reliable laboratory results of consistently good quality.

This EQAS consists of antimicrobial susceptibility testing of eight *Salmonella* isolates and eight *Campylobacter* isolates. Additionally, quality control (QC) strains *E. coli* ATCC 25922 (CCM 3954) and *C. jejuni* ATCC 33560 (CCM 6214) will be distributed to new participants.

This EQAS is specifically for NRL's on antimicrobial resistance. Laboratories designated to be NRL-AR do not need to sign up to participate but are automatically regarded as participants. You may contact the EQAS-Coordinator if you wish to inform of changes in relation to your level of participation in previous years. The EURL-AR will be able to cover the expenses for one parcel per EU Member State. Therefore, countries with more than one laboratory registered on the EURL-AR contact-list will be contacted directly to confirm which laboratory will be included for participation free of charge.

The invitation to participate in the proficiency test is extended to additional participants from official NRLs and participants from laboratories which are involved in the network but are not designated NRLs (cost for participation will be 100 EURO).

### **TO AVOID DELAY IN SHIPPING THE ISOLATES TO YOUR LABORATORY**

The content of the parcel is "UN3373, Biological Substance Category B": Eight *Salmonella* strains, eight *Campylobacter* and for new participants also the QC strains mentioned above. Please provide the EQAS coordinator with documents or other information that can simplify customs procedures (e.g. specific text that should be written on the proforma invoice). To avoid delays, we kindly ask you to send this information already at this stage.

### **TIMELINE FOR RESULTS TO BE RETURNED TO THE NATIONAL FOOD INSTITUTE**

Shipment of isolates and protocol: The isolates will be shipped in October 2015. The protocol for this proficiency test will be available for download from the website ([www.eurl-ar.eu](http://www.eurl-ar.eu)).

Submission of results: Results must be submitted to the National Food Institute **no later than December 4<sup>th</sup> 2015** via the password-protected website.



Upon reaching the deadline, each participating laboratory is kindly asked to enter the password-protected website once again to download an automatically generated evaluation report.

EQAS report: A report summarising and comparing results from all participants will be issued. In the report, laboratories will be presented coded, which ensures full anonymity. The EURL-AR and the EU Commission, only, will have access to un-coded results. The report will be publicly available.

Next EQAS: The next EURL-AR EQAS that we will send out to the EURL-AR network focuses on isolation of ESBL and ampC-producing *E.coli* from samples which is expected to be sent to participating laboratories around 1<sup>st</sup> November, 2015.

**Please contact me if you have comments or questions regarding the EQAS**

Sincerely,

Susanne Karlsmosen Pedersen (suska@food.dtu.dk)

**EQAS-Coordinator**

## Participant list

<i>Salmonella</i>	<i>Campylobacter</i>	Genotypic characterisation	Institute	Country
X	X	-	Austrian Agency for Health and Food Safety	Austria
X	X	X	Institute of Public Health	Belgium
X	X	-	National Diagnostic and Research Veterinary Institute	Bulgaria
X	X	-	Croatian Veterinary Institut	Croatia
X	X	-	Veterinary Services	Cyprus
X	X	X	State Veterinary Institute Praha	Czech Republic
X	X	-	Danish Veterinary and Food Administration, DVFA	Denmark
X	X	-	Estonian Veterinary and Food Laboratory	Estonia
X	X	-	Finnish Food Safety Authority EVIRA	Finland
X	-	-	Agence nationale de sécurité sanitaire ANSES - Fougères LERMVD	France
-	X	-	Agence nationale de sécurité sanitaire ANSES - Ploufragan - LERAP	France
X	X	X	Federal Institute for Risk Assessment	Germany
X	X	-	Veterinary Laboratory of Chalkis	Greece
X	X	-	Central Agricultural Office Veterinary Diagnostic Directorate	Hungary
X	X	-	University of Iceland	Iceland
X	X	-	Central Veterinary Research Laboratory	Ireland
X	X	X	Istituto Zooprofilattico Sperimentale delle Regioni Lazio e Toscana	Italy
X	X	-	Institute of Food Safety, Animal Health and Environment "BIOR"	Latvia
X	X	-	National Food and Veterinary Risk Assessment Institute	Lithuania
X	X	X	Laboratoire national de Santé	Luxembourg
X	X	-	Public Health Laboratory	Malta
X	X	X	Central Veterinary Institute of Wageningen UR	Netherlands
X	X	-	Food and Consumer Product Safety Authority (VWA)	Netherlands
X	X	X	Veterinærinstituttet	Norway
X	X	-	National Veterinary Research Institute	Poland
X	X	-	Laboratorio Nacional de Investigação Veterinária	Portugal
X	X	-	Institute for Diagnosis and Animal Health	Romania
X	X	-	State Veterinary and Food Institute (SVFI)	Slovakia
X	X	-	National Veterinary Institute	Slovenia
X	X	X	Laboratorio Central de Sanidad, Animal de Algete	Spain
X	X	-	VISAVET Health Surveillance Center, Complutense University	Spain
X	X	X	National Veterinary Institute, SVA	Sweden
X	X	-	Vetsuisse Faculty Bern, Institute of Veterinary Bacteriology	Switzerland
X	X	-	The Veterinary Laboratory Agency	United Kingdom

	Designated NRL-AR by the competent authority of the member state
	Non-NRL-AR enrolled by the EURL-AR
	Not a Member State of the EU

Reference values (MIC-value and interpretation) - *Salmonella*

	Ampicillin AMP		Azithromycin AZI		Cefepime FEP		Cefotaxime FOT		Cefotaxime/clav F/C		F:F/C ratio		Cefoxitin FOX		Ceftazidime TAZ		Ceftazidime/clav T/C		T:T/C ratio		Chloramphenicol CHL		Ciprofloxacin CIP		Colistin COL		Ertapenem	
EURL S-10.1	= 2	SUSC	= 8				<= 0.25	SUSC							<= 0.5	SUSC					<= 8	SUSC	= 0.06	SUSC	= 8	RESIST		
EURL S-10.2	> 64	RESIST	= 8				<= 0.25	SUSC							<= 0.5	SUSC					> 128	RESIST	= 0.03	SUSC	<= 1	SUSC		
EURL S-10.3	> 64	RESIST	= 8		= 2		= 8	RESIST	= 0.25/4	>=8	= 2	SUSC	= 1	SUSC	= 0.25/4	<8	<= 8	SUSC	= 0.5	RESIST	<= 1	SUSC	<= 0.015	SUSC				
EURL S-10.4	> 64	RESIST	= 8		= 1		= 2	RESIST	= 2/4	<8	= 4	SUSC	<= 0.5	SUSC	= 0.5/4	<8	<= 8	SUSC	= 0.03	SUSC	<= 1	SUSC	= 2	RESIST				
EURL S-10.5	> 64	RESIST	= 8				<= 0.25	SUSC					<= 0.5	SUSC			= 64	RESIST	= 0.03	SUSC	<= 1	SUSC						
EURL S-10.6	= 2	SUSC	= 8				<= 0.25	SUSC					<= 0.5	SUSC			<= 8	SUSC	= 0.03	SUSC	<= 1	SUSC						
EURL S-10.7	> 64	RESIST	> 64		> 32		> 64	RESIST	= 0.25/4	>=8	= 8	SUSC	> 128	RESIST	= 2/4	>=8	> 128	RESIST	= 0.25	RESIST	<= 1	SUSC	= 0.03	SUSC				
EURL S-10.8	> 64	RESIST	= 4		= 2		= 8	RESIST	= 0.06/4	>=8	= 2	SUSC	= 16	RESIST	= 0.25/4	>=8	<= 8	SUSC	= 0.03	SUSC	= 2	SUSC	= 0.015	SUSC				

	Gentamicin GEN		IMIPENEM IMI		MEROPENEM MER		Nalidixic acid NAL		Sulfamethoxazole SMX		TEMOCILLIN TRM		Tetracycline TETRA		TIGECYCLINE TGC		Trimethoprim TMP		ESBL-category		Relevant genes	
EURL S-10.1	<= 0.5	SUSC			= 0.06	SUSC	<= 4	SUSC	= 32	SUSC			<= 2	SUSC	= 0.5	SUSC	<= 0.25	SUSC	N/A		N/A	
EURL S-10.2	= 1	SUSC			= 0.06	SUSC	<= 4	SUSC	> 1024	RESIST			= 64	RESIST	= 0.5	SUSC	<= 0.25	SUSC	N/A		N/A	
EURL S-10.3	= 1	SUSC	= 0.5	SUSC	<= 0.03	SUSC	> 128	RESIST	= 32	SUSC	= 8		= 64	RESIST	= 0.5	SUSC	<= 0.25	SUSC	Presumptive ESBL-phenotype		TEM-1; CTX M-9	
EURL S-10.4	<= 0.5	SUSC	= 4	RESIST	= 2	RESIST	<= 4	SUSC	= 32	SUSC	> 128		<= 2	SUSC	= 0.5	SUSC	<= 0.25	SUSC	Presumptive carbapenemase phenotype		TEM-1; OXA-48	
EURL S-10.5	<= 0.5	SUSC			= 0.06	SUSC	<= 4	SUSC	> 1024	RESIST			> 64	RESIST	= 2	RESIST	> 32	RESIST	N/A		N/A	
EURL S-10.6	<= 0.5	SUSC			= 0.06	SUSC	<= 4	SUSC	= 64	SUSC			= 4	SUSC	= 0.5	SUSC	<= 0.25	SUSC	N/A		N/A	
EURL S-10.7	> 32	RESIST	= 0.5	SUSC	= 0.06	SUSC	= 8	SUSC	> 1024	RESIST	= 32		> 64	RESIST	= 1	SUSC	> 32	RESIST	Presumptive ESBL-phenotype		TEM-1; CTX M-15; SHV-12	
EURL S-10.8	= 1	SUSC	= 0.5	SUSC	= 0.06	SUSC	= 8	SUSC	= 32	SUSC	= 8		<= 2	SUSC	= 0.5	SUSC	= 0.5	SUSC	Presumptive ESBL-phenotype		TEM-52	

 Resistant

## Reference values (MIC-value and interpretation) - *Campylobacter*

Appendix 3b, page 1 of 1

Species	Code	Ciprofloxacin CIP		Erythromycin ERY		Gentamicin GEN		Nalidixic acid NAL		Streptomycin STR		Tetracycline TET	
<i>C. jejuni</i>	EURL C-10.1	= 16	RESIST	> 128	RESIST	= 0.25	SUSC	> 64	RESIST	= 1	SUSC	<= 0.5	SUSC
<i>C. coli</i>	EURL C-10.2	<= 0.12	SUSC	<= 1	SUSC	= 0.25	SUSC	= 4	SUSC	= 16	RESIST	<= 0.5	SUSC
<i>C. jejuni</i>	EURL C-10.3	<= 0.12	SUSC	<= 1	SUSC	= 0.25	SUSC	= 2	SUSC	= 1	SUSC	<= 0.5	SUSC
<i>C. coli</i>	EURL C-10.4	= 0.25	SUSC	> 128	RESIST	= 0.5	SUSC	= 8	SUSC	= 1	SUSC	= 2	SUSC
<i>C. jejuni</i>	EURL C-10.5	= 8	RESIST	<= 1	SUSC	= 0.25	SUSC	> 64	RESIST	= 1	SUSC	> 64	RESIST
<i>C. coli</i>	EURL C-10.6	> 16	RESIST	> 128	RESIST	= 0.5	SUSC	> 64	RESIST	= 2	SUSC	> 64	RESIST
<i>C. coli</i>	EURL C-10.7	= 16	RESIST	<= 1	SUSC	= 0.5	SUSC	> 64	RESIST	= 4	SUSC	<= 0.5	SUSC
<i>C. jejuni</i>	EURL C-10.8	= 16	RESIST	> 128	RESIST	> 16	RESIST	> 64	RESIST	> 16	RESIST	> 64	RESIST

 Resistant



G00-06-001/01.12.2014

**EURL-AR External Quality Assurance System 2015**  
- *Salmonella*, *Campylobacter* and optional genotypic characterisation

Id: «Lab\_no\_»

«Name»

«Institute\_\_»

«Country»

**Kgs. Lyngby, October 2015**

Dear «Name»,

Please find enclosed the bacterial strains for the EURL-AR EQAS 2015. Upon arrival to your laboratory, the strains should be stored dark and at 4°C for stabs, and dark and cool for lyophilized strains. Charcoal swabs must be subcultured straight away.

On the EURL-AR-website ([www.eurl-ar.eu](http://www.eurl-ar.eu)) the following documents relevant for the EURL-AR EQAS are available:

- Protocol for *Salmonella* and *Campylobacter* including test forms
- Instructions for Opening and Reviving Lyophilised Cultures
- Subculture and Maintenance of Quality Control Strains

We ask you to examine the eight *Salmonella* and the eight *Campylobacter* strains that we send to you by performing antimicrobial susceptibility testing. The ESBL-producing *Salmonella* strains should be characterised genotypically (optional) according to the description in the protocol. In the protocol you can find detailed description of the procedures to follow. Additionally, you can find a description of the procedure to enter your results into the interactive web database. For accessing the database, you need this username and password.

Your username: «Username»

Your password: «Password»

Please keep this document  
Your username and password will not appear in other documents

Results should be submitted to the database no later than **December 4<sup>th</sup> 2015**.

Please acknowledge receipt of this parcel immediately upon arrival (to [suska@food.dtu.dk](mailto:suska@food.dtu.dk)).  
Do not hesitate to contact us for further information.

Yours sincerely,

Susanne Karlsmosen Pedersen  
**EQAS-Coordinator**



# PROTOCOL

For antimicrobial susceptibility testing of *Salmonella*, *Campylobacter* and optional genotypic characterisation of AmpC-, ESBL- and carbapenemase-producing test strains

1	INTRODUCTION .....	1
2	OBJECTIVES .....	2
3	OUTLINE OF THE SALM/CAMP EQAS 2015 .....	2
3.1	Shipping, receipt and storage of strains .....	2
3.2	QC reference strains .....	2
3.3	Antimicrobial susceptibility testing .....	3
3.4	Optional genotypic characterisation.....	6
4	REPORTING OF RESULTS AND EVALUATION .....	6
5	HOW TO ENTER RESULTS IN THE INTERACTIVE DATABASE .....	7

## 1 INTRODUCTION

The organisation and implementation of an External Quality Assurance System (EQAS) on antimicrobial susceptibility testing (AST) of *Salmonella* and *Campylobacter* is among the tasks of the EU Reference Laboratory for Antimicrobial Resistance (EURL-AR). The Salm/Camp EQAS 2015 will include AST of eight *Salmonella* and *Campylobacter* strains and AST of reference strains *E. coli* ATCC 25922 (CCM 3954) and *C. jejuni* ATCC 33560 (CCM 6214).

The above-mentioned reference strains are included in the parcel only for new participants of the EQAS who did not receive them previously. The reference strains are original CERTIFIED cultures provided free of charge, and should be used for future internal quality control for antimicrobial susceptibility testing in your laboratory. The reference strains will not be included in the years to come. Therefore, please take proper care of these strains. Handle and maintain them as suggested in the manual 'Subculture and Maintenance of QC Strains' available on the EURL-AR website (see [www.eurl-ar.eu](http://www.eurl-ar.eu)).

## EU Reference Laboratory for Antimicrobial Resistance External Quality Assurance System (EQAS) 2015



Various aspects of the proficiency test scheme may from time to time be subcontracted. When subcontracting occurs it is placed with a competent subcontractor and the National Food Institute is responsible to the scheme participants for the subcontractor's work.

## 2 OBJECTIVES

This EQAS aims to support laboratories to assess and, if necessary, to improve the quality of results obtained by AST of pathogens of food- and animal-origin, with special regard to *Salmonella* and *Campylobacter*. Further objectives are to evaluate and improve the comparability of surveillance data on antimicrobial susceptibility of *Salmonella* and *Campylobacter* reported to EFSA by different laboratories.

## 3 OUTLINE OF THE SALM/CAMP EQAS 2015

### 3.1 Shipping, receipt and storage of strains

In October 2015, the National Reference Laboratories for Antimicrobial Resistance (NRL-AR) will receive a parcel containing eight *Salmonella* and *Campylobacter* strains from the National Food Institute. This parcel will also contain reference strains, but only for participants who did not receive them previously. All strains belong to UN3373, Biological substance, category B. Extended spectrum beta-lactamase (ESBL)-producing strains as well as carbapenemase producing strains are included in the selected material and are part of the optional EQAS-item, consisting of characterization of genes conferring ESBL- or carbapenemase production.

The reference strains are shipped lyophilised, the *Campylobacter* test strains are shipped as a charcoal swabs and the *Salmonella* test strains are stab cultures. On arrival, the stab cultures and the charcoal swabs must be subcultured, and all cultures should be adequately stored until testing. A suggested procedure for reconstitution of the lyophilised reference strains is presented below.

### 3.2 QC reference strains

For a suggested procedure for reconstitution of the lyophilised, please refer to the document 'Instructions for opening and reviving lyophilised cultures' on the EURL-AR-website (see [www.eurl-ar.eu](http://www.eurl-ar.eu)).

Note that, for the testing of the *E. coli* ATCC25922 reference strain, the two compounds, sulfamethoxazole and sulfisoxazole, are regarded as comparable, i.e. the obtained MIC-value from the testing of sulfamethoxazole will be evaluated against the acceptance range listed in CLSI M100 for sulfisoxazole.



### 3.4 Antimicrobial susceptibility testing

The strains should be tested for susceptibility to the antimicrobials listed in Tables 1, 2 and 3, using the method implemented in your laboratory for performing monitoring for EFSA and applying the interpretative criteria listed below.

Participants should perform minimum inhibitory concentration (MIC) determination using the methods stated in the EC regulation EC 652/2013. For interpretation of the results, use the cut-off values listed in Tables 1, 2 and 3 (except where indicated) represent the current epidemiological cut-off values developed by EUCAST ([www.eucast.org](http://www.eucast.org)), and allow categorisation of bacterial isolates into two categories; resistant or susceptible. A categorisation as intermediate is not accepted.

As the current regulation and recommendations focus on MIC testing only, results obtained by disk diffusion cannot be submitted.

#### 3.4.1 *Salmonella*

The interpretative criteria that should be applied for categorizing the *Salmonella* test strain as resistant or susceptible are those listed in Tables 1 and 2.

Table 1: Antimicrobials recommended for AST of *Salmonella* spp. and interpretative criteria according to table 1 in EC regulation 652/2013

Antimicrobial	MIC (µg/mL) (R>)
Ampicillin (AMP)	8
Azithromycin (AZI)	Not available*
Cefotaxime (FOT)	0.5
Ceftazidime (TAZ)	2
Chloramphenicol (CHL)	16
Ciprofloxacin (CIP)	0.06
Colistin (COL)	2
Gentamicin (GEN)	2
Meropenem (MERO)	0.125
Nalidixic acid (NAL)	16
Sulfonamides (SMX)	256**
Tetracycline (TET)	8
Tigecycline (TGC)	1***
Trimethoprim (TMP)	2

\* Participants are requested to upload the MIC value obtained without selecting an interpretation.

\*\* CLSI M100 Table 2A

\*\*\* Data from EUCAST is available for *S. Enteritidis*, *S. Typhimurium*, *S. Typhi* and *S. Paratyphi* (for the purpose of this proficiency test, the ECOFF at 1 is applied)

## EU Reference Laboratory for Antimicrobial Resistance External Quality Assurance System (EQAS) 2015



Table 2: Antimicrobials recommended for additional AST of *Salmonella* spp. resistant to cefotaxime, ceftazidime or meropenem and interpretative criteria according to table 4 in EC regulation 652/2013

Antimicrobial	MIC (µg/mL) (R>)
Cefepime, FEP	Not available*
Cefotaxime, FOT	0.5
Cefotaxime + clavulanic acid (F/C)	Not applicable
Cefoxitin, FOX	8
Ceftazidime, TAZ	2
Ceftazidime+ clavulanic acid (T/C)	Not applicable
Ertapenem, ETP	0.06
Imipenem, IMI	1
Meropenem, MERO	0.125
Temocillin, TRM	Not available*

\* Participants are requested to upload the MIC value obtained without selecting an interpretation

### Plasmid-mediated quinolone resistance

When performing antimicrobial susceptibility testing of the *Salmonella* test strains, the interpretative criteria listed in Table 1 should be able to detect plasmid mediated quinolone resistant test strains.

### Extended-beta-lactam- and carbapenem resistance

Confirmatory tests for AmpC-, ESBL- and carbapenemase production are **mandatory** on all strains resistant to cefotaxime (FOT), ceftazidime (TAZ) or meropenem and should be performed by testing the second panel of antimicrobials (Table 2 in this document corresponding to Table 4 in EC regulation 652/2013).

Confirmatory tests for AmpC-, ESBL- and carbapenemase production require the use of both cefotaxime (FOT) and ceftazidime (TAZ) alone and in combination with a  $\beta$ -lactamase inhibitor (clavulanic acid). Synergy is defined either as a  $\geq 3$  twofold concentration decrease in an MIC for either antimicrobial agent tested in combination with clavulanic acid vs. its MIC when tested alone (MIC FOT : FOT/Cl or TAZ : TAZ/Cl ratio  $\geq 8$ ) (CLSI M100 Table 2A; Enterobacteriaceae). The presence of synergy indicates ESBL production. Resistance to cefepime gives further indication of ESBL production, but is not essential. Confirmatory test for carbapenemase production requires the testing of meropenem (MERO).

Detection of AmpC-type beta-lactamases can be performed by testing the bacterium for susceptibility to cefoxitin (FOX). Resistance to FOX could indicate the presence of an AmpC-type beta-lactamase that may be verified by PCR and sequencing.



The classification of the phenotypic results should be based on the most recent EFSA recommendations (EFSA 2012<sup>1</sup>) indicating the strains as:

- Presumptive ESBL: strains with positive synergy test, susceptible to cefoxitin and resistant to cefepime
- Presumptive ESBL+pAmpC: strains with positive or negative synergy test, resistant to cefoxitin and resistant to cefepime
- Presumptive pAmpC phenotype: strains with negative synergy test, resistant to cefoxitin and susceptible to cefepime
- Presumptive carbapenemase phenotype: strain resistant to meropenem
- Unusual phenotype: any other combinations

*We recommend, however, that strains showing synergy with clavulanic acid for at least one of the third generation cephalosporins (cefotaxime or ceftazidime) should be considered ESBL-producing, independently of the cefepime result.*

### 3.4.2 *Campylobacter*

For AST of *Campylobacter*, MIC methods should be applied, i.e. broth or agar dilution methods using incubation at 36-37°C for 48 hours or 42°C for 24 hours.

Table 3: Antimicrobials recommended for AST of *Campylobacter jejuni* and *C. coli* and interpretative criteria according to table 1 in EC regulation 652/2013

Antimicrobial	<i>C. jejuni</i>	<i>C. coli</i>
	MIC (µg/mL) (R>)	MIC (µg/mL) (R>)
Ciprofloxacin (CIP)	0.5	0.5
Erythromycin (ERY)	4	8
Gentamicin (GEN)	2	2
Nalidixic acid (NAL)	16	16
Streptomycin (STR)	4	4
Tetracycline (TET)	1	2

### Identification of *Campylobacter* species

Species identification of the *Campylobacter* test strains must be performed by the NRLs using in-house methods or adopting the protocol available on the EURL-AR website under: <http://eurl-ar.eu/233-protocols.htm>.

<sup>1</sup> European Food Safety Authority; Technical specifications on the harmonised monitoring and reporting of antimicrobial resistance in *Salmonella*, *Campylobacter* and indicator *Escherichia coli* and *Enterococcus* spp. bacteria transmitted through food. EFSA Journal 2012; 10(6):2742. [64 pp.] doi:10.2903/j.efsa.2012.2742. Available online: [www.efsa.europa.eu/efsajournal](http://www.efsa.europa.eu/efsajournal)



### 3.5 Optional genotypic characterisation

For the optional genotypic characterisation of the AmpC-, ESBL- or carbapenemase producing *Salmonella* test strains, the requested results are the genes conferring AmpC-, ESBL- or carbapenemase -production harboured in the test strains. The genes included in the test are the following: ACC, ACT, CMY, CTX, DHA, FOX, GES, IMP, KPC, MOX, NDM, OXA, PER, SHV, TEM, VEB, and VIM. The database lists the relevant variants of the genes.

When uploading the results in the database, the identified genes will be evaluated against the expected results. The results will be evaluated on the detected gene (ACC-, ACT-, CMY-, etc.) as well as the variant identified.

The method used for the genotypic characterisation should be your laboratory's routine method. The expected results listed in the database are those obtained by the EURL-AR.

## 4 REPORTING OF RESULTS AND EVALUATION

Test forms are available for recording your results before you enter them into the interactive web database.

We recommend reading carefully the description reported in paragraph 5 before entering your results in the web database. **Results must be submitted no later than December 4<sup>th</sup> 2015.** After the deadline when all participants have uploaded results, you will be able to login to the database once again, and to view and print an automatically generated report evaluating your results. Results in agreement with the expected interpretation are categorised as 'correct', while results deviating from the expected interpretation are categorised as 'incorrect'.

If you experience difficulties in entering your results, please contact us directly.

All results will be summarized in a report which will be publicly available. The data in the report will be presented with laboratory codes. A laboratory code is known to the individual laboratory, whereas the complete list of laboratories and their codes is confidential and known only to the EURL-AR and the EU Commission. All conclusions will be public.

If you have questions, please do not hesitate to contact the EQAS Coordinator:

Susanne Karlsmosen Pedersen  
National Food Institute,  
Technical University of Denmark  
Søtofts Plads, Building 221, DK-2800 Lyngby  
Denmark  
Tel: +45 3588 6601  
Fax: +45 3588 6341  
E-mail: [suska@food.dtu.dk](mailto:suska@food.dtu.dk)



## EU Reference Laboratory for Antimicrobial Resistance External Quality Assurance System (EQAS) 2015



### 5 HOW TO ENTER RESULTS IN THE INTERACTIVE DATABASE

Please read carefully this paragraph before entering the web page.

Remember that you need by your side the completed test forms.

Enter the EURL-AR EQAS 2015 start web page (<http://eurl-ar.food.dtu.dk>), write your username and password (lower-case) and press enter. Your username and password are indicated in the letter following your strains. Do not hesitate to contact us if you experience problems with the login.

You can browse back and forth by using the Home or back keys, but please remember to save your inputs before.

Click on either “*Salmonella* test results” or “*Campylobacter* test results” for input of test results.

Click on "Start of Data Entry - Methods"

In the next page, you navigate among fields with the Tab-key and the mouse.

Complete the fields related to the method used for antimicrobial susceptibility testing and the brand of MIC trays, etc.

When submitting *Campylobacter* results, fill in the incubation conditions applied for susceptibility testing of *Campylobacter* – 36°C/48h or 42°C/24h.

Click on "save and go to next page"

In the data entry pages, you enter the species (for *Campylobacter* only), the obtained MIC-value and the interpretation (R, resistant or S, susceptible) for each *Salmonella* and *Campylobacter* strain.

For *Salmonella*, remember to also report the results for the ESBL detection tests.

If you did not test for susceptibility to a given antimicrobial, please leave the field empty.

Click on "save and go to next page"

When uploading data on the reference strains, please enter MIC values in µg/ml. Remember to use the operator keys to show symbols like “equal to”, etc.

Click on “save“.

Review the input pages by browsing through them and make corrections if necessary. Remember to save a page if you make corrections. If you press home a page without saving changes, you will see an error screen. In this case, click on “save“ to save your results, browse back to the page and then continue.

Please complete the evaluation form.

Before approving your input, please be sure that you have filled in all the relevant fields as **YOU CAN ONLY APPROVE ONCE!** The approval blocks your data entry in the interactive database.

If you have performed the optional genotypic characterisation:

## EU Reference Laboratory for Antimicrobial Resistance External Quality Assurance System (EQAS) 2015



Click on “Gene test” and follow the description in the database for upload of the results of the optional genotypic characterization. Approve your input. Be sure that you have filled in all the results before approval. The approval blocks your data entry in the interactive database, but allows you to see the submitted results.

## *Salmonella, Campylobacter* and genetic characterisation

# TEST FORMS

Name:

Name of laboratory:

Name of institute:

City:

Country:

E-mail:

Fax:

Comments:



## TEST FORM

Does your laboratory have an accreditation for performing *Salmonella* AST? ☐ Yes ☐ No

Which method did you use for antimicrobial susceptibility testing of *Salmonella* in this EQAS:

☐ Broth microdilution

☐ Agar dilution

Brand of microbroth plates/agar:

Incubation conditions: °C/ h

How many *Salmonella* isolates does your laboratory annually isolate:

How many *Salmonella* isolates does your laboratory annually test for antimicrobial susceptibility by a MIC method:

Which method was followed for the preparation of the inoculum (please describe)

- Which standard was followed (TREK, CLSI...)
- Which solvent was used for the preparation of the 0.5 McFarland solution (water, saline)
- Please describe in detail how you prepared the dilution of the inoculum (including the volume in final MH-dilution and intended dilution level; e.g. diluted 1:1000 by adding 10µl of 0.5 McFarland solution in 10ml MH broth, for an expected inoculum of  $1 \times 10^5$  CFU/ml)

Comments or additional information:



## TEST FORM

Does your laboratory have an accreditation for *Campylobacter* AST? ☐ Yes ☐ No

Incubation conditions: ☐ 36-37°C / 48h ☐ 42°C / 24h

Method used for antimicrobial susceptibility testing of *Campylobacter* in this EQAS::

- ☐ Broth microdilution  
☐ Agardilution

Brand of microbroth plates/agar:

How many *Campylobacter* isolates does your laboratory annually isolate:

How many *Campylobacter* isolates does your laboratory annually susceptibility test:

Which method was followed for the preparation of the inoculum (please describe)

- Which standard was followed (TREK, CLSI...)
- Which solvent was used for the preparation of the 0.5 McFarland solution (water, saline)
- Please describe in detail how you prepared the dilution of the inoculum (including the volume in final MH-dilution and intended dilution level; e.g. diluted 1:1000 by adding 10µl of 0.5 McFarland solution in 10ml MH broth, for an expected inoculum of  $1 \times 10^5$  CFU/ml)

Comments or additional information:

## TEST FORM

Strain	Antimicrobial	Results and interpretation		
		≤ >	MIC-value (µg/ml)	S / R
<i>Salmonella</i> EURL S. 10.X	Ampicillin, AMP			
	Azithromycin, AZI			
	Cefotaxime, FOT			
	Ceftazidime, TAZ			
	Chloramphenicol, CHL			
	Ciprofloxacin CIP			
	Colistin, COL			
	Gentamicin, GEN			
	Meropenem, MERO			
	Nalidixic acid, NAL			
	Sulfamethoxazole, SMX			
	Tetracycline, TET			
	Tigecycline, TGC			
	Trimethoprim, TMP			

All strains resistant to cefotaxime (FOT), ceftazidime (TAZ) or meropenem (MERO) must be included for testing in the second panel as part of confirmatory tests for ESBL-, AmpC or carbapenemase production. See further description in the protocol section '3.3.1 *Salmonella*'.

Strain	Antimicrobial	Results and interpretation		
		≤ >	MIC-value (µg/ml)	S / R
<i>Salmonella</i> EURL S. 10.X	Cefepime, FEP			
	Cefotaxime, FOT			
	Cefotaxime + clavulanic acid (F/C)			
	Cefoxitin, FOX			
	Ceftazidime, TAZ			
	Ceftazidime+ clavulanic acid (T/C)			
	Ertapenem, ETP			
	Imipenem, IMI			
	Meropenem, MERO			
	Temocillin, TRM			

- |  |  |  |
|--|--|--|
| <input type="checkbox"/> Presumptive ESBL        | <input type="checkbox"/> Presumptive pAmpC         | <input type="checkbox"/> Unusual phenotype               |
| <input type="checkbox"/> Presumptive ESBL+ pAmpC | <input type="checkbox"/> Presumptive carbapenemase | <input type="checkbox"/> No ESBL, AmpC- or carbapenemase |

Comments (include optional genotype or other results):

## TEST FORM

Antimicrobial susceptibility testing of reference strain *E. coli* ATCC 25922

	Antimicrobial	MIC-value (µg/ml)
1 <sup>st</sup> panel	Ampicillin, AMP	
	Azithromycin, AZI	
	Cefotaxime, FOT	
	Ceftazidime, TAZ	
	Chloramphenicol, CHL	
	Ciprofloxacin, CIP	
	Colistin, COL	
	Gentamicin, GEN	
	Meropenem, MERO	
	Nalidixic acid, NAL	
	Sulfamethoxazole, SMX*	
	Tetracycline, TET	
	Tigecycline, TGC	
	Trimethoprim, TMP	
2 <sup>nd</sup> panel	Cefepime, FEP	
	Cefotaxime, FOT	
	Cefotaxime + clavulanic acid (F/C)	
	Cefoxitin, FOX	
	Ceftazidime, TAZ	
	Ceftazidime+ clavulanic acid (T/C)	
	Ertapenem, ETP	
	Imipenem, IMI	
	Meropenem, MERO	
	Temocillin, TRM	

\* for the testing of the *E. coli* ATCC25922 reference strain, sulfamethoxazole and sulfisoxazole, are regarded as comparable, i.e. the obtained MIC-value from the testing of sulfamethoxazole will be evaluated against the acceptance range listed in CLSI M100 for sulfisoxazole (CLSI M100, Table 3).





## TEST FORM

Strain	Antimicrobial	Interpretation	
		MIC-value (µg/ml)	S / R
<b><i>Campylobacter</i></b> <b>EURL C-10.X</b>  <input type="checkbox"/> <i>C. jejuni</i>  <input type="checkbox"/> <i>C. coli</i>	Ciprofloxacin		
	Erythromycin		
	Gentamicin		
	Nalidixic acid		
	Streptomycin		
	Tetracycline		
<b><i>Campylobacter</i></b> <b>EURL C-10.X</b>  <input type="checkbox"/> <i>C. jejuni</i>  <input type="checkbox"/> <i>C. coli</i>	Ciprofloxacin		
	Erythromycin		
	Gentamicin		
	Nalidixic acid		
	Streptomycin		
	Tetracycline		
<b><i>Campylobacter</i></b> <b>EURL C-10.X</b>  <input type="checkbox"/> <i>C. jejuni</i>  <input type="checkbox"/> <i>C. coli</i>	Ciprofloxacin		
	Erythromycin		
	Gentamicin		
	Nalidixic acid		
	Streptomycin		
	Tetracycline		
<b><i>Campylobacter</i></b> <b>EURL C-10.X</b>  <input type="checkbox"/> <i>C. jejuni</i>  <input type="checkbox"/> <i>C. coli</i>	Ciprofloxacin		
	Erythromycin		
	Gentamicin		
	Nalidixic acid		
	Streptomycin		
	Tetracycline		

## TEST FORM

Susceptibility testing of *Campylobacter jejuni* reference strain ATCC 33560

Strain	Antimicrobial	MIC-value (µg/ml)	
		36 °C/48 hours	42 °C/24 hours
<i>C. jejuni</i> ATCC 33560	Ciprofloxacin		
	Erythromycin		
	Nalidixic acid		
	Tetracycline		

### For Agar dilution:

Susceptibility testing of *Campylobacter jejuni* reference strain ATCC 33560

Strain	Antimicrobial	MIC-value (µg/ml)
<i>C. jejuni</i> ATCC 33560	Ciprofloxacin	
	Erythromycin	
	Gentamicin	
	Nalidixic acid	
	Tetracycline	

## TEST FORM – genotypic characterisation

Genotypic characterisation of the test strains

Strain code:	Method used: If PCR-methods, additional information should be given below			
Gene:  <input type="checkbox"/> Found <input type="checkbox"/> Tested, not found	<input type="checkbox"/> Published method , reference:			
	<input type="checkbox"/> In-house method			
	Primer used 5'→3':			
	Primer used 3'→5':			
Gene:  <input type="checkbox"/> Found <input type="checkbox"/> Tested, not found	<input type="checkbox"/> Published method , reference:			
	<input type="checkbox"/> In-house method			
	Primer used 5'→3':			
	Primer used 3'→5':			
Gene:  <input type="checkbox"/> Found <input type="checkbox"/> Tested, not found	<input type="checkbox"/> Published method , reference:			
	<input type="checkbox"/> In-house method			
	Primer used 5'→3':			
	Primer used 3'→5':			
Gene:  <input type="checkbox"/> Found <input type="checkbox"/> Tested, not found	<input type="checkbox"/> Published method , reference:			
	<input type="checkbox"/> In-house method			
	Primer used 5'→3':			
	Primer used 3'→5':			
Gene:  <input type="checkbox"/> Found <input type="checkbox"/> Tested, not found	<input type="checkbox"/> Published method , reference:			
	<input type="checkbox"/> In-house method			
	Primer used 5'→3':			
	Primer used 3'→5':			

Comments:



# INSTRUCTIONS FOR OPENING AND REVIVING LYOPHILISED CULTURES

*Instructions adjusted from Czech Collection of Microorganisms (CCM) document 'Instructions for Opening and Reviving of Freeze-Dried Bacteria and Fungi' available on <http://www.sci.muni.cz>.*

Lyophilised cultures are supplied in vacuum-sealed ampoules. Care should be taken in opening the ampoule. All instructions given below should be followed closely to ensure the safety of the person who opens the ampoule and to prevent contamination of the culture.

- a. Check the number of the culture on the label inside the ampoule
- b. Make a file cut on the ampoule near the middle of the plug (see Figure 1)
- c. Disinfect the ampoule with alcohol-dampened gauze or alcohol-dampened cotton wool from just below the plug to the pointed end
- d. Apply a red-hot glass rod to the file cut to crack the glass and allow air to enter slowly into the ampoule
- e. Remove the pointed end of the ampoule into disinfectant
- f. Add about 0.3 ml appropriate broth to the dried suspension using a sterile Pasteur pipette and mix carefully to avoid creating aerosols. Transfer the contents to one or more suitable solid and /or liquid media
- g. Incubate the inoculated medium at appropriate conditions for several days
- h. Autoclave or disinfect effectively the used Pasteur pipette, the plug and all the remains of the original ampoule before discarding

## Notes:

- Cultures should be grown on media and under conditions as recommended in the CCM catalogue (see <http://www.sci.muni.cz>)
- Cultures may need at least one subculturing before they can be optimally used in experiments
- Unopened ampoules should be kept in a dark and cool place!

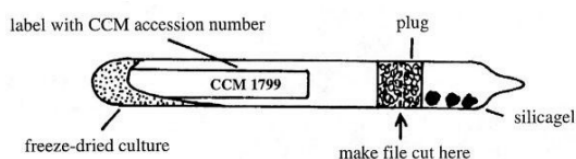


Figure 1: from CCM document 'Instructions for Opening and Reviving of Freeze-Dried Bacteria and Fungi' available on <http://www.sci.muni.cz>



# SUBCULTURE AND MAINTENANCE OF QUALITY CONTROL STRAINS

---

## 1.1 Purpose

Improper storage and repeated subculturing of bacteria can produce alterations in antimicrobial susceptibility test results. The Clinical and Laboratory Standards Institute (CLSI, formerly NCCLS) has published a guideline for Quality Control (QC) stock culture maintenance to ensure consistent antimicrobial susceptibility test results.

## 1.2 References

M100-S24, January 2014 (Performance Standards for Antimicrobial Susceptibility Testing)

M7-A9, January 2012 (Methods for Dilution Antimicrobial Susceptibility Test for Bacteria That Grow Aerobically; Approved Standard)

## 1.3 Definition of Terms

Reference Culture: A reference culture is a microorganism preparation that is acquired from a culture type collection.

Reference Stock Culture: A reference stock culture is a microorganism preparation that is derived from a reference culture. Guidelines and standards outline how reference stock cultures must be processed and stored.

Working Stock Cultures: A working stock culture is growth derived from a reference stock culture. Guidelines and standards outline how working stock cultures must be processed and how often they can be subcultured.

Subcultures (Passages): A subculture is simply the transfer of established microorganism growth on media to fresh media. The subsequent growth on the fresh media constitutes a subculture or passage. Growing a reference culture or reference stock culture from its preserved status (frozen or lyophilized) is not a subculture. The preserved microorganism is not in a stage of established growth until it is thawed or hydrated and grown for the first time

## 1.4 Important Considerations

- Do not use disc diffusion strains for MIC determination.
- Obtain QC strains from a reliable source such as ATCC
- CLSI requires that QC be performed either on the same day or weekly (only after 30 day QC validation)
- Any changes in materials or procedure must be validated with QC before implemented
- For example: Agar and broth methods may give different QC ranges for drugs such as glycopeptides, aminoglycosides and macrolides
- Periodically perform colony counts to check the inoculum preparation procedure

## EU Reference Laboratory for Antimicrobial Resistance External Quality Assurance System (EQAS)



- Ideally, test values should be in the middle of the acceptable range
- Graphing QC data points over time can help identify changes in data helpful for troubleshooting problems

### 1.5 Storage of Reference Strains

#### Preparation of stock cultures

- Use a suitable stabilizer such as 50% fetal calf serum in broth, 10-15% glycerol in tryptic soy broth, defibrinated sheep blood or skim milk to prepare multiple aliquots.
- Store at -20°C, -70°C or liquid nitrogen. (Alternatively, freeze dry.)
- Before using rejuvenated strains for QC, subculture to check for purity and viability.

#### Working cultures

- Set up on agar slants with appropriate medium, store at 4-8°C and subculture weekly.
- Replace the working strain with a stock culture at least monthly.
- If a change in the organisms inherent susceptibility occurs, obtain a fresh stock culture or a new strain from a reference culture collection e.g. ATCC.

### 1.6 Frequency of Testing

#### Weekly vs. daily testing

Weekly testing is possible if the lab can demonstrate satisfactory performance with daily testing as follows:

- Documentation showing reference strain results from 30 consecutive test days were within the acceptable range.
- For each antimicrobial/organism combination, no more than 3 out of 30 MIC values may be outside the acceptable range.

When the above are fulfilled, each quality control strain may be tested once a week and whenever any reagent component is changed.

#### Corrective Actions

If an MIC is outside the range in weekly testing, corrective action is required as follows:

- Repeat the test if there is an obvious error e.g. wrong strain or incubation conditions used
- If there is no obvious error, return to daily control testing

The problem is considered resolved only after the reference strain is tested for 5 consecutive days and each drug/organism result is within specification on each day.

If the problem cannot be resolved, continue daily testing until the errors are identified.

Repeat the 30 days validation before resuming weekly testing.




**EU Reference Laboratory for Antimicrobial Resistance  
External Quality Assurance System (EQAS)**



**DAILY MIC QC CHART**

**Appendix A. Quality Control Protocol Flow Charts**

**Quality Control (QC) Protocol: Daily Testing**

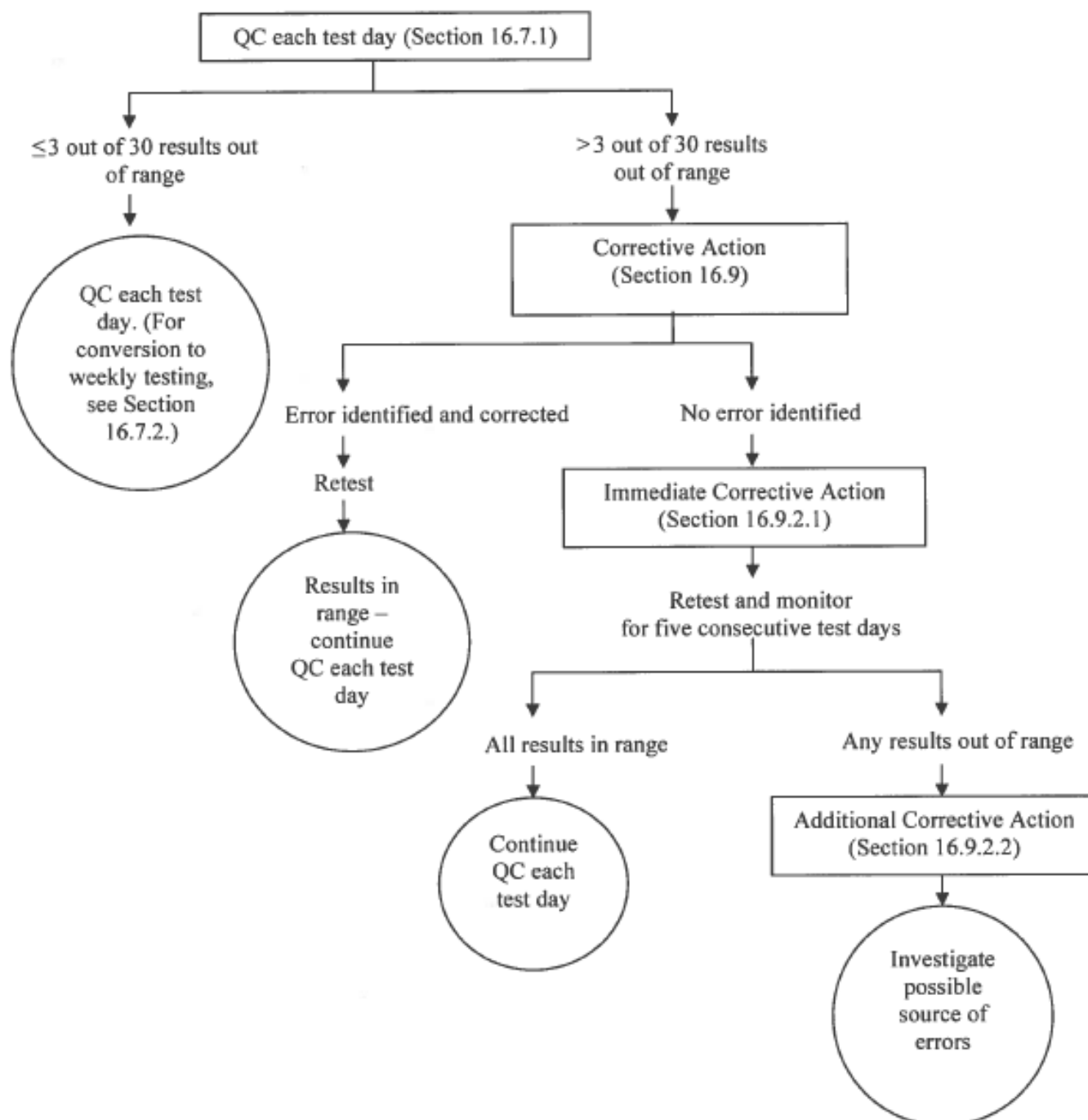


Reference: CLSI M7-A9, page 46



## Appendix A. (Continued)

### QC Protocol: Weekly Testing

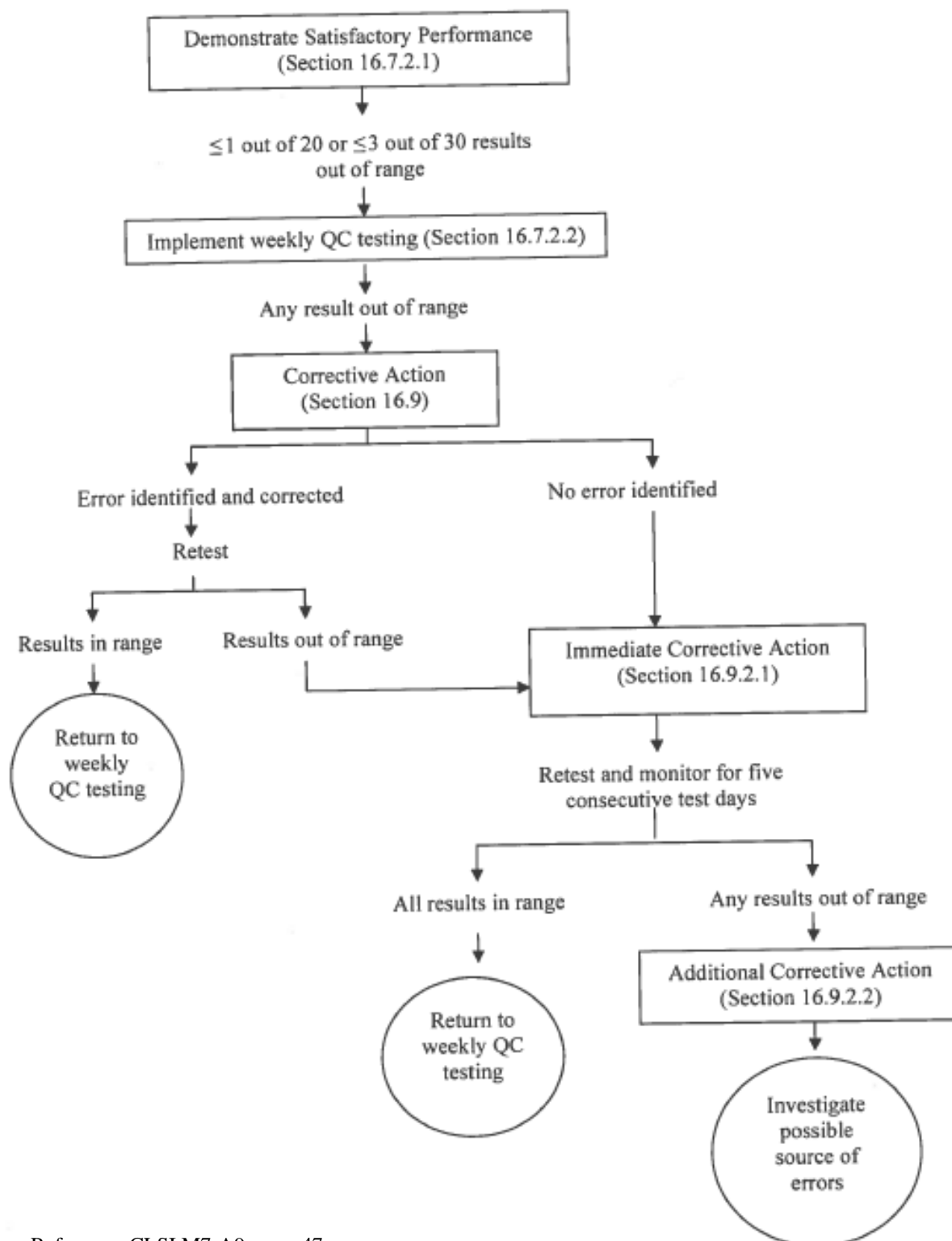

Reference: CLSI M7-A9, page 47

## Quality Control ranges for ATCC reference strains

<i>E. coli</i> ATCC 25922	
Antimicrobial	MIC
Ampicillin, AMP	2-8
Azithromycin, AZT	none
Cefepime, FEP	0.015-0.12
Cefotaxime, FOT	0.03-0.12
Cefotaxime + clavulanic acid, F/C	none
Cefoxitin, FOX	2-8
Ceftazidime, CAZ	0.06-0.5
Ceftazidime + clavulanic acid, T/C	none
Chloramphenicol, CHL	2-8
Ciprofloxacin, CIP	0.004-0.016
Colistin, COL	0.25-2
Ertapenem, ETP	0.004-0.016
Gentamicin, GEN	0.25-1
Imipenem, IMI	0.06-0.25
Meropenem, MERO	0.008-0.06
Nalidixic acid, NAL	1-4
Sulfisoxazole, FIS	8-32
Temocillin, TRM	none
Tetracycline, TET	0.5-2
Tigecycline, TGC	0.03-0.25
Trimethoprim, TMP	0.5-2

MIC ranges (µg/mL) are according to CLSI M100 S25 (range for ciprofloxacin and ertapenem extended to include 0.016).

<i>Campylobacter jejuni</i> ATCC 33560				
Antimicrobial	Microbroth (36-37°C/48h)	Microbroth (42°C/24h)	Agar dilution (36-37°C/48h)	Agar dilution (42°C/24h)
Ciprofloxacin, CIP	0.06-0.25	0.03-0.12	0.12-1	0.06-0.5
Erythromycin, ERY	0.5-2	0.25-2	1-8	1-4
Gentamicin, GEN	0.5-2	0.25-2	0.5-2	0.5-4
Nalidixic acid, NAL	4-16	4-16	None	None
Tetracycline, TET	0.25-2	0.25-1	None	None

MIC ranges (µg/mL) are according to CLSI (VET01-S2)

Test results from the reference strain *E. coli* ATCC 25922

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
2	1	Ampicillin	=	4	2	8	1	MIC
2	1	Azithromycin	=	4	-	-	-	MIC
2	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
2	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
2	1	Chloramphenicol	<=	8	2	8	1	MIC
2	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
2	1	Colistin	<=	1	0.25	2	1	MIC
2	1	Gentamicin	<=	0.5	0.25	1	1	MIC
2	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
2	1	Nalidixic acid	<=	4	1	4	1	MIC
2	1	Sulfamethoxazole	=	16	8	32	1	MIC
2	1	Tetracycline	<=	2	0.5	2	1	MIC
2	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
2	1	Trimethoprim	=	0.5	0.5	2	1	MIC
6	1	Ampicillin	=	4	2	8	1	MIC
6	1	Azithromycin	=	4	-	-	-	MIC
6	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
6	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
6	1	Chloramphenicol	<=	8	2	8	1	MIC
6	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
6	1	Colistin	<=	1	0.25	2	1	MIC
6	1	Gentamicin	<=	0.5	0.25	1	1	MIC
6	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
6	1	Nalidixic acid	<=	4	1	4	1	MIC
6	1	Sulfamethoxazole	=	32	8	32	1	MIC
6	1	Tetracycline	<=	2	0.5	2	1	MIC
6	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
6	1	Trimethoprim	=	1	0.5	2	1	MIC
6	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
6	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
6	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
6	2	Cefoxitin	=	4	2	8	1	MIC
6	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
6	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
6	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
6	2	Imipenem	=	0.25	0.06	0.25	1	MIC
6	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
6	2	Temocillin	=	4	-	-	-	MIC
9	1	Ampicillin	=	4	2	8	1	MIC
9	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
9	1	Chloramphenicol	=	8	2	8	1	MIC
9	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
9	1	Colistin	<=	1	0.25	2	1	MIC
9	1	Gentamicin	<=	0.5	0.25	1	1	MIC
9	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
9	1	Nalidixic acid	<=	4	1	4	1	MIC
9	1	Sulfamethoxazole	=	16	8	32	1	MIC
9	1	Tetracycline	<=	2	0.5	2	1	MIC
9	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
9	1	Trimethoprim	=	1	0.5	2	1	MIC
9	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
9	2	Cefoxitin	=	4	2	8	1	MIC
9	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
9	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
9	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
9	2	Meropenem	<=	0.03	0.008	0.06	1	MIC

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
11	1	Ampicillin	=	4	2	8	1	MIC
11	1	Azithromycin	=	8	-	-	-	MIC
11	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
11	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
11	1	Chloramphenicol	<=	8	2	8	1	MIC
11	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
11	1	Colistin	<=	1	0.25	2	1	MIC
11	1	Gentamicin	=	1	0.25	1	1	MIC
11	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
11	1	Nalidixic acid	<=	4	1	4	1	MIC
11	1	Sulfamethoxazole	=	32	8	32	1	MIC
11	1	Tetracycline	<=	2	0.5	2	1	MIC
11	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
11	1	Trimethoprim	=	0.5	0.5	2	1	MIC
11	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
11	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
11	2	Cefotaxime/clavulanic acid	<=	0.12	-	-	-	MIC
11	2	Cefoxitin	=	4	2	8	1	MIC
11	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
11	2	Ceftazidime/clavulanic acid	=	0.25	-	-	-	MIC
11	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
11	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
11	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
11	2	Temocillin	=	32	-	-	-	MIC
12	1	Ampicillin	=	8	2	8	1	MIC
12	1	Azithromycin	=	4	-	-	-	MIC
12	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
12	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
12	1	Chloramphenicol	<=	8	2	8	1	MIC
12	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
12	1	Colistin	<=	1	0.25	2	1	MIC
12	1	Gentamicin	=	1	0.25	1	1	MIC
12	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
12	1	Nalidixic acid	<=	4	1	4	1	MIC
12	1	Sulfamethoxazole	=	16	8	32	1	MIC
12	1	Tetracycline	<=	2	0.5	2	1	MIC
12	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
12	1	Trimethoprim	=	0.5	0.5	2	1	MIC
12	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
12	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
12	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
12	2	Cefoxitin	=	4	2	8	1	MIC
12	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
12	2	Ceftazidime/clavulanic acid	=	0.25	-	-	-	MIC
12	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
12	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
12	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
12	2	Temocillin	=	16	-	-	-	MIC
16	1	Ampicillin	=	4	2	8	1	MIC
16	1	Azithromycin	=	4	-	-	-	MIC
16	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
16	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
16	1	Chloramphenicol	<=	8	2	8	1	MIC
16	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
16	1	Colistin	<=	1	0.25	2	1	MIC
16	1	Gentamicin	<=	0.5	0.25	1	1	MIC
16	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
16	1	Nalidixic acid	<=	4	1	4	1	MIC
16	1	Sulfamethoxazole	=	32	8	32	1	MIC
16	1	Tetracycline	<=	2	0.5	2	1	MIC
16	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
16	1	Trimethoprim	=	0.5	0.5	2	1	MIC
16	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
16	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
16	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
16	2	Cefoxitin	=	4	2	8	1	MIC
16	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
16	2	Ceftazidime/clavulanic acid	=	0.25	-	-	-	MIC
16	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
16	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
16	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
16	2	Temocillin	=	16	-	-	-	MIC

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
17	1	Ampicillin	=	8	2	8	1	MIC
17	1	Azithromycin	=	4	-	-	-	MIC
17	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
17	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
17	1	Chloramphenicol	<=	8	2	8	1	MIC
17	1	Ciprofloxacin	<=	0.01	0.004	0.016	1	MIC
17	1	Colistin	<=	1	0.25	2	1	MIC
17	1	Gentamicin	=	1	0.25	1	1	MIC
17	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
17	1	Nalidixic acid	<=	4	1	4	1	MIC
17	1	Sulfamethoxazole	=	32	8	32	1	MIC
17	1	Tetracycline	<=	2	0.5	2	1	MIC
17	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
17	1	Trimethoprim	=	0.5	0.5	2	1	MIC
17	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
17	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
17	2	Cefotaxime/clavulanic acid	=	0.12	-	-	-	MIC
17	2	Cefoxitin	=	4	2	8	1	MIC
17	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
17	2	Ceftazidime/clavulanic acid	=	0.25	-	-	-	MIC
17	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
17	2	Imipenem	=	0.25	0.06	0.25	1	MIC
17	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
17	2	Temocillin	=	16	-	-	-	MIC
18	1	Ampicillin	=	2	2	8	1	MIC
18	1	Azithromycin	<=	2	-	-	-	MIC
18	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
18	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
18	1	Chloramphenicol	<=	8	2	8	1	MIC
18	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
18	1	Colistin	<=	1	0.25	2	1	MIC
18	1	Gentamicin	<=	0.5	0.25	1	1	MIC
18	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
18	1	Nalidixic acid	<=	4	1	4	1	MIC
18	1	Sulfamethoxazole	=	32	8	32	1	MIC
18	1	Tetracycline	<=	2	0.5	2	1	MIC
18	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
18	1	Trimethoprim	=	1	0.5	2	1	MIC
18	2	Cefepime	=	0.12	0.015	0.125	1	MIC
18	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
18	2	Cefoxitin	=	4	2	8	1	MIC
18	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
18	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
18	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
18	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
18	2	Temocillin	=	8	-	-	-	MIC
19	1	Ampicillin	=	4	2	8	1	MIC
19	1	Azithromycin	=	4	-	-	-	MIC
19	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
19	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
19	1	Chloramphenicol	<=	8	2	8	1	MIC
19	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
19	1	Colistin	<=	1	0.25	2	1	MIC
19	1	Gentamicin	<=	0.5	0.25	1	1	MIC
19	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
19	1	Nalidixic acid	<=	4	1	4	1	MIC
19	1	Sulfamethoxazole	=	32	8	32	1	MIC
19	1	Tetracycline	<=	2	0.5	2	1	MIC
19	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
19	1	Trimethoprim	=	5	0.5	2	0	MIC
19	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
19	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
19	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
19	2	Cefoxitin	=	4	2	8	1	MIC
19	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
19	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
19	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
19	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
19	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
19	2	Temocillin	=	8	-	-	-	MIC



Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
20	1	Ampicillin	=	4	2	8	1	MIC
20	1	Azithromycin	=	8	-	-	-	MIC
20	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
20	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
20	1	Chloramphenicol	<=	8	2	8	1	MIC
20	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
20	1	Colistin	<=	1	0.25	2	1	MIC
20	1	Gentamicin	<=	0.5	0.25	1	1	MIC
20	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
20	1	Nalidixic acid	<=	4	1	4	1	MIC
20	1	Sulfamethoxazole	=	32	8	32	1	MIC
20	1	Tetracycline	<=	2	0.5	2	1	MIC
20	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
20	1	Trimethoprim	=	0.5	0.5	2	1	MIC
20	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
20	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
20	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
20	2	Cefoxitin	=	2	2	8	1	MIC
20	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
20	2	Ceftazidime/clavulanic acid	=	0.25	-	-	-	MIC
20	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
20	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
20	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
20	2	Temocillin	=	8	-	-	-	MIC
21	1	Ampicillin	=	4	2	8	1	MIC
21	1	Azithromycin	=	4	-	-	-	MIC
21	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
21	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
21	1	Chloramphenicol	<=	8	2	8	1	MIC
21	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
21	1	Colistin	<=	1	0.25	2	1	MIC
21	1	Gentamicin	<=	0.5	0.25	1	1	MIC
21	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
21	1	Nalidixic acid	<=	4	1	4	1	MIC
21	1	Sulfamethoxazole	=	32	8	32	1	MIC
21	1	Tetracycline	<=	2	0.5	2	1	MIC
21	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
21	1	Trimethoprim	=	0.5	0.5	2	1	MIC
21	2	Cefepime	=	0.5	0.015	0.125	0	MIC
21	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
21	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
21	2	Cefoxitin	=	8	2	8	1	MIC
21	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
21	2	Ceftazidime/clavulanic acid	=	0.25	-	-	-	MIC
21	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
21	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
21	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
21	2	Temocillin	=	16	-	-	-	MIC
22	1	Ampicillin	=	2	2	8	1	MIC
22	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
22	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
22	1	Chloramphenicol	<=	8	2	8	1	MIC
22	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
22	1	Colistin	<=	1	0.25	2	1	MIC
22	1	Gentamicin	<=	0.5	0.25	1	1	MIC
22	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
22	1	Nalidixic acid	<=	4	1	4	1	MIC
22	1	Tetracycline	<=	2	0.5	2	1	MIC
22	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
22	1	Trimethoprim	=	0.5	0.5	2	1	MIC
22	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
22	2	Cefoxitin	=	4	2	8	1	MIC
22	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
22	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
22	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
22	2	Meropenem	<=	0.03	0.008	0.06	1	MIC

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
23	1	Ampicillin	=	4	2	8	1	MIC
23	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
23	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
23	1	Chloramphenicol	<=	8	2	8	1	MIC
23	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
23	1	Colistin	<=	1	0.25	2	1	MIC
23	1	Gentamicin	=	1	0.25	1	1	MIC
23	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
23	1	Nalidixic acid	<=	4	1	4	1	MIC
23	1	Sulfamethoxazole	=	16	8	32	1	MIC
23	1	Tetracycline	<=	2	0.5	2	1	MIC
23	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
23	1	Trimethoprim	=	0.5	0.5	2	1	MIC
23	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
23	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
23	2	Cefoxitin	=	2	2	8	1	MIC
23	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
23	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
23	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
23	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
25	1	Ampicillin	=	4	2	8	1	MIC
25	1	Azithromycin	=	4	-	-	-	MIC
25	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
25	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
25	1	Chloramphenicol	<=	8	2	8	1	MIC
25	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
25	1	Colistin	<=	1	0.25	2	1	MIC
25	1	Gentamicin	<=	0.5	0.25	1	1	MIC
25	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
25	1	Nalidixic acid	<=	4	1	4	1	MIC
25	1	Sulfamethoxazole	<=	8	8	32	1	MIC
25	1	Tetracycline	<=	2	0.5	2	1	MIC
25	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
25	1	Trimethoprim	=	0.5	0.5	2	1	MIC
26	1	Ampicillin	=	4	2	8	1	MIC
26	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
26	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
26	1	Chloramphenicol	<=	8	2	8	1	MIC
26	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
26	1	Colistin	<=	1	0.25	2	1	MIC
26	1	Gentamicin	=	1	0.25	1	1	MIC
26	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
26	1	Nalidixic acid	<=	4	1	4	1	MIC
26	1	Sulfamethoxazole	=	16	8	32	1	MIC
26	1	Tetracycline	<=	2	0.5	2	1	MIC
26	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
26	1	Trimethoprim	=	0.5	0.5	2	1	MIC
29	1	Ampicillin	=	8	2	8	1	MIC
29	1	Chloramphenicol	=	8	2	8	1	MIC
29	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
29	1	Colistin	=	1	0.25	2	1	MIC
29	1	Gentamicin	=	0.5	0.25	1	1	MIC
29	1	Nalidixic acid	<=	4	1	4	1	MIC
29	1	Sulfamethoxazole	=	16	8	32	1	MIC
29	1	Tetracycline	=	2	0.5	2	1	MIC
29	1	Trimethoprim	=	1	0.5	2	1	MIC
29	2	Cefepime	=	0.12	0.015	0.125	1	MIC
29	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
29	2	Ceftazidime	=	0.25	0.06	0.5	1	MIC
29	2	Ertapenem	=	0.015	0.004	0.016	1	MIC
29	2	Imipenem	=	0.12	0.06	0.25	1	MIC
29	2	Meropenem	=	0.03	0.008	0.06	1	MIC

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
30	1	Ampicillin	=	4	2	8	1	MIC
30	1	Azithromycin	=	4	-	-	-	MIC
30	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
30	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
30	1	Chloramphenicol	<=	8	2	8	1	MIC
30	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
30	1	Colistin	<=	1	0.25	2	1	MIC
30	1	Gentamicin	<=	0.5	0.25	1	1	MIC
30	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
30	1	Nalidixic acid	<=	4	1	4	1	MIC
30	1	Sulfamethoxazole	=	16	8	32	1	MIC
30	1	Tetracycline	<=	2	0.5	2	1	MIC
30	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
30	1	Trimethoprim	=	0.5	0.5	2	1	MIC
30	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
30	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
30	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
30	2	Cefoxitin	=	2	2	8	1	MIC
30	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
30	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
30	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
30	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
30	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
30	2	Temocillin	=	8	-	-	-	MIC
32	1	Ampicillin	=	4	2	8	1	MIC
32	1	Azithromycin	=	4	-	-	-	MIC
32	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
32	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
32	1	Chloramphenicol	<=	8	2	8	1	MIC
32	1	Ciprofloxacin	<=	0.15	0.004	0.016	1	MIC
32	1	Colistin	<=	1	0.25	2	1	MIC
32	1	Gentamicin	<=	0.5	0.25	1	1	MIC
32	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
32	1	Nalidixic acid	<=	4	1	4	1	MIC
32	1	Sulfamethoxazole	=	16	8	32	1	MIC
32	1	Tetracycline	<=	2	0.5	2	1	MIC
32	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
32	1	Trimethoprim	=	0.5	0.5	2	1	MIC
32	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
32	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
32	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
32	2	Cefoxitin	=	4	2	8	1	MIC
32	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
32	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
32	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
32	2	Imipenem	=	0.25	0.06	0.25	1	MIC
32	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
32	2	Temocillin	=	8	-	-	-	MIC
33	1	Ampicillin	=	4	2	8	1	MIC
33	1	Azithromycin	=	4	-	-	-	MIC
33	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
33	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
33	1	Chloramphenicol	<=	8	2	8	1	MIC
33	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
33	1	Colistin	<=	1	0.25	2	1	MIC
33	1	Gentamicin	<=	0.5	0.25	1	1	MIC
33	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
33	1	Nalidixic acid	<=	4	1	4	1	MIC
33	1	Sulfamethoxazole	=	32	8	32	1	MIC
33	1	Tetracycline	<=	2	0.5	2	1	MIC
33	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
33	1	Trimethoprim	=	0.5	0.5	2	1	MIC
33	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
33	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
33	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
33	2	Cefoxitin	=	4	2	8	1	MIC
33	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
33	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
33	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
33	2	Imipenem	=	0.25	0.06	0.25	1	MIC
33	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
33	2	Temocillin	=	8	-	-	-	MIC

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
34	1	Ampicillin	=	4	2	8	1	MIC
34	1	Azithromycin	=	4	-	-	-	MIC
34	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
34	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
34	1	Chloramphenicol	<=	8	2	8	1	MIC
34	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
34	1	Colistin	<=	1	0.25	2	1	MIC
34	1	Gentamicin	<=	0.5	0.25	1	1	MIC
34	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
34	1	Nalidixic acid	<=	4	1	4	1	MIC
34	1	Sulfamethoxazole	=	16	8	32	1	MIC
34	1	Tetracycline	<=	2	0.5	2	1	MIC
34	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
34	1	Trimethoprim	=	0.5	0.5	2	1	MIC
34	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
34	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
34	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
34	2	Cefoxitin	=	4	2	8	1	MIC
34	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
34	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
34	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
34	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
34	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
34	2	Temocillin	=	16	-	-	-	MIC
36	1	Ampicillin	=	4	2	8	1	MIC
36	1	Azithromycin	<=	2	-	-	-	MIC
36	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
36	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
36	1	Chloramphenicol	<=	8	2	8	1	MIC
36	1	Ciprofloxacin	=	0.03	0.004	0.016	0	MIC
36	1	Colistin	<=	1	0.25	2	1	MIC
36	1	Gentamicin	=	1	0.25	1	1	MIC
36	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
36	1	Nalidixic acid	<=	4	1	4	1	MIC
36	1	Sulfamethoxazole	<=	8	8	32	1	MIC
36	1	Tetracycline	<=	2	0.5	2	1	MIC
36	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
36	1	Trimethoprim	<=	0.25	0.5	2	0	MIC
36	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
36	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
36	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
36	2	Cefoxitin	=	4	2	8	1	MIC
36	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
36	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
36	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
36	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
36	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
36	2	Temocillin	=	16	-	-	-	MIC
37	1	Ampicillin	=	4	2	8	1	AGA
37	1	Azithromycin	=	4	-	-	-	AGA
37	1	Cefotaxime	<=	0.25	0.03	0.125	1	AGA
37	1	Ceftazidime	<=	0.5	0.06	0.5	1	AGA
37	1	Chloramphenicol	<=	8	2	8	1	AGA
37	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	AGA
37	1	Colistin	<=	1	0.25	2	1	AGA
37	1	Gentamicin	<=	0.5	0.25	1	1	AGA
37	1	Meropenem	<=	0.03	0.008	0.06	1	AGA
37	1	Nalidixic acid	<=	4	1	4	1	AGA
37	1	Sulfamethoxazole	=	32	8	32	1	AGA
37	1	Tetracycline	<=	2	0.5	2	1	AGA
37	1	Tigecycline	<=	0.25	0.03	0.25	1	AGA
37	1	Trimethoprim	=	1	0.5	2	1	AGA
37	2	Cefepime	=	0.03	0.015	0.125	1	AGA
37	2	Cefotaxime	<=	0.25	0.03	0.125	1	AGA
37	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	AGA
37	2	Cefoxitin	=	4	2	8	1	AGA
37	2	Ceftazidime	<=	0.25	0.06	0.5	1	AGA
37	2	Ceftazidime/clavulanic acid	<=	0.125	-	-	-	AGA
37	2	Ertapenem	<=	0.015	0.004	0.016	1	AGA
37	2	Imipenem	<=	0.125	0.06	0.25	1	AGA
37	2	Meropenem	<=	0.03	0.008	0.06	1	AGA
37	2	Temocillin	=	2	-	-	-	AGA

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
39	1	Ampicillin	=	4	2	8	1	MIC
39	1	Azithromycin	=	4	-	-	-	MIC
39	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
39	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
39	1	Chloramphenicol	<=	8	2	8	1	MIC
39	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
39	1	Colistin	<=	1	0.25	2	1	MIC
39	1	Gentamicin	<=	0.5	0.25	1	1	MIC
39	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
39	1	Nalidixic acid	<=	4	1	4	1	MIC
39	1	Sulfamethoxazole	=	16	8	32	1	MIC
39	1	Tetracycline	<=	2	0.5	2	1	MIC
39	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
39	1	Trimethoprim	=	0.5	0.5	2	1	MIC
39	2	Cefepime	=	0.06	0.015	0.125	1	MIC
39	2	Cefotaxime	=	0.25	0.03	0.125	0	MIC
39	2	Cefoxitin	=	4	2	8	1	MIC
39	2	Ceftazidime	=	0.25	0.06	0.5	1	MIC
39	2	Ertapenem	=	0.015	0.004	0.016	1	MIC
39	2	Imipenem	=	0.25	0.06	0.25	1	MIC
39	2	Meropenem	=	0.03	0.008	0.06	1	MIC
39	2	Temocillin	=	16	-	-	-	MIC
40	1	Ampicillin	=	2	2	8	1	MIC
40	1	Cefotaxime	=	0.12	0.03	0.125	1	MIC
40	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
40	1	Chloramphenicol	<=	8	2	8	1	MIC
40	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
40	1	Colistin	<=	1	0.25	2	1	MIC
40	1	Gentamicin	<=	0.5	0.25	1	1	MIC
40	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
40	1	Nalidixic acid	<=	4	1	4	1	MIC
40	1	Sulfamethoxazole	=	16	8	32	1	MIC
40	1	Tetracycline	<=	2	0.5	2	1	MIC
40	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
40	1	Trimethoprim	=	0.5	0.5	2	1	MIC
40	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
40	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
40	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
40	2	Cefoxitin	=	2	2	8	1	MIC
40	2	Ceftazidime	=	0.5	0.06	0.5	1	MIC
40	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
40	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
40	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
40	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
42	1	Ampicillin	=	4	2	8	1	MIC
42	1	Azithromycin	=	8	-	-	-	MIC
42	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
42	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
42	1	Chloramphenicol	<=	8	2	8	1	MIC
42	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
42	1	Colistin	<=	1	0.25	2	1	MIC
42	1	Gentamicin	<=	0.5	0.25	1	1	MIC
42	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
42	1	Nalidixic acid	<=	4	1	4	1	MIC
42	1	Sulfamethoxazole	=	64	8	32	0	MIC
42	1	Tetracycline	<=	2	0.5	2	1	MIC
42	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
42	1	Trimethoprim	=	1	0.5	2	1	MIC
42	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
42	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
42	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
42	2	Cefoxitin	=	4	2	8	1	MIC
42	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
42	2	Ceftazidime/clavulanic acid	<=	0.125	-	-	-	MIC
42	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
42	2	Imipenem	=	0.25	0.06	0.25	1	MIC
42	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
42	2	Temocillin	=	8	-	-	-	MIC

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
45	1	Ampicillin	=	4	2	8	1	MIC
45	1	Azithromycin	=	4	-	-	-	MIC
45	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
45	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
45	1	Chloramphenicol	<=	8	2	8	1	MIC
45	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
45	1	Colistin	<=	1	0.25	2	1	MIC
45	1	Gentamicin	<=	0.5	0.25	1	1	MIC
45	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
45	1	Nalidixic acid	<=	4	1	4	1	MIC
45	1	Sulfamethoxazole	=	32	8	32	1	MIC
45	1	Tetracycline	<=	2	0.5	2	1	MIC
45	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
45	1	Trimethoprim	=	1	0.5	2	1	MIC
45	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
45	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
45	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
45	2	Cefoxitin	=	4	2	8	1	MIC
45	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
45	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
45	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
45	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
45	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
45	2	Temocillin	=	16	-	-	-	MIC
56	1	Ampicillin	=	4	2	8	1	MIC
56	1	Azithromycin	=	4	-	-	-	MIC
56	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
56	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
56	1	Chloramphenicol	<=	8	2	8	1	MIC
56	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
56	1	Colistin	<=	1	0.25	2	1	MIC
56	1	Gentamicin	<=	0.5	0.25	1	1	MIC
56	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
56	1	Nalidixic acid	<=	4	1	4	1	MIC
56	1	Sulfamethoxazole	=	16	8	32	1	MIC
56	1	Tetracycline	<=	2	0.5	2	1	MIC
56	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
56	1	Trimethoprim	=	1	0.5	2	1	MIC
56	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
56	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
56	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
56	2	Cefoxitin	=	2	2	8	1	MIC
56	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
56	2	Ceftazidime/clavulanic acid	<=	0.12	-	-	-	MIC
56	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
56	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
56	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
56	2	Temocillin	=	4	-	-	-	MIC
58	1	Ampicillin	=	8	2	8	1	MIC
58	1	Azithromycin	=	8	-	-	-	MIC
58	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
58	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
58	1	Chloramphenicol	<=	8	2	8	1	MIC
58	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
58	1	Colistin	<=	1	0.25	2	1	MIC
58	1	Gentamicin	<=	0.5	0.25	1	1	MIC
58	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
58	1	Nalidixic acid	<=	4	1	4	1	MIC
58	1	Sulfamethoxazole	=	32	8	32	1	MIC
58	1	Tetracycline	<=	2	0.5	2	1	MIC
58	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
58	1	Trimethoprim	=	1	0.5	2	1	MIC
58	2	Cefepime	<=	0.06	0.015	0.125	1	MIC
58	2	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
58	2	Cefotaxime/clavulanic acid	<=	0.06	-	-	-	MIC
58	2	Cefoxitin	=	4	2	8	1	MIC
58	2	Ceftazidime	<=	0.25	0.06	0.5	1	MIC
58	2	Ceftazidime/clavulanic acid	=	0.25	-	-	-	MIC
58	2	Ertapenem	<=	0.015	0.004	0.016	1	MIC
58	2	Imipenem	<=	0.12	0.06	0.25	1	MIC
58	2	Meropenem	<=	0.03	0.008	0.06	1	MIC
58	2	Temocillin	=	16	-	-	-	MIC

Lab no.	Panel	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method
60	1	Ampicillin	=	4	2	8	1	MIC
60	1	Azithromycin	=	4	-	-	-	MIC
60	1	Cefotaxime	<=	0.25	0.03	0.125	1	MIC
60	1	Ceftazidime	<=	0.5	0.06	0.5	1	MIC
60	1	Chloramphenicol	<=	8	2	8	1	MIC
60	1	Ciprofloxacin	<=	0.015	0.004	0.016	1	MIC
60	1	Colistin	<=	1	0.25	2	1	MIC
60	1	Gentamicin	<=	0.5	0.25	1	1	MIC
60	1	Meropenem	<=	0.03	0.008	0.06	1	MIC
60	1	Nalidixic acid	<=	4	1	4	1	MIC
60	1	Sulfamethoxazole	<=	8	8	32	1	MIC
60	1	Tetracycline	<=	2	0.5	2	1	MIC
60	1	Tigecycline	<=	0.25	0.03	0.25	1	MIC
60	1	Trimethoprim	<=	0.25	0.5	2	<b>0</b>	MIC

MIC: Microbroth dilution

AGA: Agar dilution



Test results from the reference strain *C. jejuni* ATCC 33560

Lab no.	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method	36-37°C/48h	42°C/24h
2	Ciprofloxacin	=	0.12	0.06	0.25	1	MIC	X	
2	Erythromycin	=	1	0.5	2	1	MIC	X	
2	Gentamicin	=	0.25	0.5	2	0	MIC	X	
2	Nalidixic acid	=	8	4	16	1	MIC	X	
2	Streptomycin	=	1	-	-	-	MIC	X	
2	Tetracycline	=	2	0.25	2	1	MIC	X	
6	Ciprofloxacin	<=	0.12	0.03	0.125	1	MIC		X
6	Erythromycin	<=	1	0.25	2	1	MIC		X
6	Gentamicin	=	1	0.25	2	1	MIC		X
6	Nalidixic acid	=	8	4	16	1	MIC		X
6	Streptomycin	=	4	-	-	-	MIC		X
6	Tetracycline	<=	0.5	0.25	1	1	MIC		X
9	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
9	Erythromycin	=	2	0.5	2	1	MIC	X	
9	Gentamicin	=	1	0.5	2	1	MIC	X	
9	Nalidixic acid	=	4	4	16	1	MIC	X	
9	Streptomycin	=	1	-	-	-	MIC	X	
9	Tetracycline	=	1	0.25	2	1	MIC	X	
11	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
11	Erythromycin	<=	1	0.5	2	1	MIC	X	
11	Gentamicin	=	1	0.5	2	1	MIC	X	
11	Nalidixic acid	=	8	4	16	1	MIC	X	
11	Streptomycin	=	4	-	-	-	MIC	X	
11	Tetracycline	=	1	0.25	2	1	MIC	X	
12	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
12	Erythromycin	<=	1	0.5	2	1	MIC	X	
12	Gentamicin	=	2	0.5	2	1	MIC	X	
12	Nalidixic acid	=	8	4	16	1	MIC	X	
12	Streptomycin	=	4	-	-	-	MIC	X	
12	Tetracycline	<=	0.5	0.25	2	1	MIC	X	
14	Ciprofloxacin	<=	0.125	0.03	0.125	1	MIC		X
14	Erythromycin	<=	1	0.25	2	1	MIC		X
14	Gentamicin	=	1	0.25	2	1	MIC		X
14	Nalidixic acid	=	8	4	16	1	MIC		X
14	Streptomycin	=	4	-	-	-	MIC		X
14	Tetracycline	=	1	0.25	1	1	MIC		X
17	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
17	Erythromycin	<=	1	0.5	2	1	MIC	X	
17	Gentamicin	=	2	0.5	2	1	MIC	X	
17	Nalidixic acid	=	8	4	16	1	MIC	X	
17	Streptomycin	=	8	-	-	-	MIC	X	
17	Tetracycline	<=	0.5	0.25	2	1	MIC	X	
18	Ciprofloxacin	<=	0.12	0.03	0.125	1	MIC		X
18	Erythromycin	<=	1	0.25	2	1	MIC		X
18	Gentamicin	=	1	0.25	2	1	MIC		X
18	Nalidixic acid	=	8	4	16	1	MIC		X
18	Tetracycline	=	1	0.25	1	1	MIC		X
19	Ciprofloxacin	<=	0.12	0.06	0.25	1	MIC	X	
19	Erythromycin	<=	1	0.5	2	1	MIC	X	
19	Gentamicin	=	1	0.5	2	1	MIC	X	
19	Nalidixic acid	=	8	4	16	1	MIC	X	
19	Streptomycin	=	4	-	-	-	MIC	X	
19	Tetracycline	=	2	0.25	2	1	MIC	X	

Lab no.	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method	36-37°C/48h	42°C/24h
20	Ciprofloxacin	=	0.25	0.03	0.125	0	MIC		X
20	Erythromycin	<=	1	0.25	2	1	MIC		X
20	Gentamicin	=	1	0.25	2	1	MIC		X
20	Nalidixic acid	=	8	4	16	1	MIC		X
20	Streptomycin	=	2	-	-	-	MIC		X
20	Tetracycline	=	4	0.25	1	0	MIC		X
21	Ciprofloxacin	<=	0.12	0.03	0.125	1	MIC		X
21	Erythromycin	<=	1	0.25	2	1	MIC		X
21	Gentamicin	=	0.5	0.25	2	1	MIC		X
21	Nalidixic acid	=	4	4	16	1	MIC		X
21	Streptomycin	=	2	-	-	-	MIC		X
21	Tetracycline	<=	0.5	0.25	1	1	MIC		X
22	Ciprofloxacin	<=	0.06	0.03	0.125	1	MIC		X
22	Erythromycin	<=	1	0.25	2	1	MIC		X
22	Nalidixic acid	=	4	4	16	1	MIC		X
22	Tetracycline	=	1	0.25	1	1	MIC		X
23	Ciprofloxacin	=	0.12	0.03	0.125	1	MIC		X
23	Erythromycin	<=	0.5	0.25	2	1	MIC		X
23	Gentamicin	=	1	0.25	2	1	MIC		X
23	Nalidixic acid	=	4	4	16	1	MIC		X
23	Tetracycline	=	0.5	0.25	1	1	MIC		X
25	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
25	Erythromycin	<=	1	0.5	2	1	MIC	X	
25	Gentamicin	=	0.25	0.5	2	0	MIC	X	
25	Nalidixic acid	=	8	4	16	1	MIC	X	
25	Streptomycin	=	1	-	-	-	MIC	X	
25	Tetracycline	=	2	0.25	2	1	MIC	X	
26	Ciprofloxacin	<=	0.12	0.06	0.25	1	MIC	X	
26	Erythromycin	<=	1	0.5	2	1	MIC	X	
26	Gentamicin	=	0.5	0.5	2	1	MIC	X	
26	Nalidixic acid	=	4	4	16	1	MIC	X	
26	Streptomycin	=	2	-	-	-	MIC	X	
26	Tetracycline	<=	0.5	0.25	2	1	MIC	X	
29	Ciprofloxacin	=	0.06	0.03	0.125	1	MIC		X
29	Erythromycin	=	1	0.25	2	1	MIC		X
29	Gentamicin	=	2	0.25	2	1	MIC		X
29	Nalidixic acid	=	4	4	16	1	MIC		X
29	Tetracycline	=	0.25	0.25	1	1	MIC		X
30	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
30	Erythromycin	<=	1	0.5	2	1	MIC	X	
30	Gentamicin	=	1	0.5	2	1	MIC	X	
30	Nalidixic acid	=	8	4	16	1	MIC	X	
30	Streptomycin	=	2	-	-	-	MIC	X	
30	Tetracycline	=	2	0.25	2	1	MIC	X	
32	Ciprofloxacin	<=	0.12	0.06	0.25	1	MIC	X	
32	Erythromycin	<=	1	0.5	2	1	MIC	X	
32	Gentamicin	=	0.5	0.5	2	1	MIC	X	
32	Nalidixic acid	=	4	4	16	1	MIC	X	
32	Streptomycin	=	4	-	-	-	MIC	X	
32	Tetracycline	=	1	0.25	2	1	MIC	X	
33	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
33	Erythromycin	<=	1	0.5	2	1	MIC	X	
33	Gentamicin	=	0.5	0.5	2	1	MIC	X	
33	Nalidixic acid	=	8	4	16	1	MIC	X	
33	Streptomycin	=	2	-	-	-	MIC	X	
33	Tetracycline	<=	0.5	0.25	2	1	MIC	X	

Lab no.	Antimicrobial	Operator	Value	Low limit	High limit	Mark	Method	36-37°C/48h	42°C/24h
34	Ciprofloxacin	<=	0.12	0.06	0.25	1	MIC	X	
34	Erythromycin	<=	1	0.5	2	1	MIC	X	
34	Gentamicin	=	0.5	0.5	2	1	MIC	X	
34	Nalidixic acid	=	8	4	16	1	MIC	X	
34	Streptomycin	=	0.5	-	-	-	MIC	X	
34	Tetracycline	=	2	0.25	2	1	MIC	X	
36	Ciprofloxacin	<=	0.06	0.03	0.125	1	MIC		X
36	Erythromycin	<=	0.5	0.25	2	1	MIC		X
36	Gentamicin	=	2	0.25	2	1	MIC		X
36	Nalidixic acid	=	16	4	16	1	MIC		X
36	Streptomycin	=	2	-	-	-	MIC		X
36	Tetracycline	<=	0.12	0.25	1	0	MIC		X
37	Ciprofloxacin	=	0.25	0.12	1	1	AGA	X	
37	Erythromycin	=	2	1	8	1	AGA	X	
37	Gentamicin	=	2	0.5	2	1	AGA	X	
37	Nalidixic acid	=	8	-	-	-	AGA	X	
37	Streptomycin	=	4	-	-	-	AGA	X	
37	Tetracycline	=	1	-	-	-	AGA	X	
39	Ciprofloxacin	=	0.12	0.06	0.25	1	MIC	X	
39	Erythromycin	=	1	0.5	2	1	MIC	X	
39	Gentamicin	=	0.5	0.5	2	1	MIC	X	
39	Nalidixic acid	=	8	4	16	1	MIC	X	
39	Tetracycline	=	1	0.25	2	1	MIC	X	
40	Ciprofloxacin	<=	0.12	0.03	0.125	1	MIC		X
40	Erythromycin	<=	1	0.25	2	1	MIC		X
40	Gentamicin	=	0.25	0.25	2	1	MIC		X
40	Nalidixic acid	=	4	4	16	1	MIC		X
40	Tetracycline	<=	0.5	0.25	1	1	MIC		X
42	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
42	Erythromycin	<=	1	0.5	2	1	MIC	X	
42	Gentamicin	<=	0.12	0.5	2	0	MIC	X	
42	Nalidixic acid	=	8	4	16	1	MIC	X	
42	Streptomycin	=	1	-	-	-	MIC	X	
42	Tetracycline	=	2	0.25	2	1	MIC	X	
45	Ciprofloxacin	<=	0.12	0.06	0.25	1	MIC	X	
45	Erythromycin	<=	1	0.5	2	1	MIC	X	
45	Gentamicin	=	1	0.5	2	1	MIC	X	
45	Nalidixic acid	=	8	4	16	1	MIC	X	
45	Streptomycin	=	4	-	-	-	MIC	X	
45	Tetracycline	=	2	0.25	2	1	MIC	X	
56	Ciprofloxacin	<=	0.12	0.06	0.25	1	MIC	X	
56	Erythromycin	<=	1	0.5	2	1	MIC	X	
56	Gentamicin	=	1	0.5	2	1	MIC	X	
56	Nalidixic acid	=	8	4	16	1	MIC	X	
56	Streptomycin	=	4	-	-	-	MIC	X	
56	Tetracycline	=	1	0.25	2	1	MIC	X	
58	Ciprofloxacin	<=	0.12	0.06	0.25	1	MIC	X	
58	Erythromycin	<=	1	0.5	2	1	MIC	X	
58	Gentamicin	=	0.5	0.5	2	1	MIC	X	
58	Nalidixic acid	=	8	4	16	1	MIC	X	
58	Streptomycin	=	2	-	-	-	MIC	X	
58	Tetracycline	<=	0.5	0.25	2	1	MIC	X	
60	Ciprofloxacin	=	0.25	0.06	0.25	1	MIC	X	
60	Erythromycin	<=	1	0.5	2	1	MIC	X	
60	Gentamicin	=	1	0.5	2	1	MIC	X	
60	Nalidixic acid	=	8	4	16	1	MIC	X	
60	Streptomycin	=	4	-	-	-	MIC	X	
60	Tetracycline	=	1	0.25	2	1	MIC	X	

*Salmonella* - expected and obtained interpretation

Antimicrobial	Strain	Panel	Expected	% R	% S	No. correct	No. incorrect
Ampicillin AMP	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	R	100	0	31	0
	EURL S-10.3	Panel 1	R	100	0	31	0
	EURL S-10.4	Panel 1	R	97	3	31	1
	EURL S-10.5	Panel 1	R	100	0	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	R	97	3	31	1
	EURL S-10.8	Panel 1	R	100	0	31	0
Cefotaxime FOT	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	S	0	100	31	0
	EURL S-10.3	Panel 1	R	100	0	31	0
	EURL S-10.3	Panel 2	R	100	0	31	0
	EURL S-10.4	Panel 1	R	100	0	31	0
	EURL S-10.4	Panel 2	R	100	0	31	0
	EURL S-10.5	Panel 1	S	0	100	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	R	100	0	31	0
	EURL S-10.7	Panel 2	R	100	0	31	0
	EURL S-10.8	Panel 1	R	100	0	31	0
	EURL S-10.8	Panel 2	R	100	0	31	0
Cefoxitin FOX	EURL S-10.3	Panel 2	S	0	100	31	0
	EURL S-10.4	Panel 2	S	0	100	31	0
	EURL S-10.7	Panel 2	S	3	97	31	1
	EURL S-10.8	Panel 2	S	0	100	31	0
Ceftazidime TAZ	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	S	0	100	31	0
	EURL S-10.3	Panel 1	S	3	97	31	1
	EURL S-10.3	Panel 2	S	3	97	31	1
	EURL S-10.4	Panel 1	S	0	100	30	0
	EURL S-10.4	Panel 2	S	0	100	31	0
	EURL S-10.5	Panel 1	S	0	100	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	R	100	0	31	0
	EURL S-10.7	Panel 2	R	100	0	31	0
	EURL S-10.8	Panel 1	R	100	0	31	0
	EURL S-10.8	Panel 2	R	100	0	31	0
Chloramphenicol CHL	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	R	100	0	31	0
	EURL S-10.3	Panel 1	S	0	100	31	0
	EURL S-10.4	Panel 1	S	0	100	31	0
	EURL S-10.5	Panel 1	R	97	3	31	1
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	R	100	0	31	0
	EURL S-10.8	Panel 1	S	0	100	31	0
Ciprofloxacin CIP	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	S	0	100	31	0
	EURL S-10.3	Panel 1	R	94	6	31	2
	EURL S-10.4	Panel 1	S	0	100	31	0
	EURL S-10.5	Panel 1	S	0	100	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	R	97	3	31	1
	EURL S-10.8	Panel 1	S	0	100	31	0

Antimicrobial	Strain	Panel	Expected	% R	% S	No. correct	No. incorrect
Colistin COL	EURL S-10.1	Panel 1	R	90	10	31	3
	EURL S-10.2	Panel 1	S	0	100	30	0
	EURL S-10.3	Panel 1	S	3	97	31	1
	EURL S-10.4	Panel 1	S	0	100	31	0
	EURL S-10.5	Panel 1	S	0	100	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	S	3	97	31	1
	EURL S-10.8	Panel 1	S	0	100	31	0
Ertapenem ETP	EURL S-10.3	Panel 2	S	0	100	31	0
	EURL S-10.4	Panel 2	R	100	0	31	0
	EURL S-10.7	Panel 2	S	0	100	31	0
	EURL S-10.8	Panel 2	S	0	100	31	0
Gentamicin GEN	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	S	0	100	31	0
	EURL S-10.3	Panel 1	S	3	97	31	1
	EURL S-10.4	Panel 1	S	0	100	31	0
	EURL S-10.5	Panel 1	S	0	100	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	R	100	0	31	0
	EURL S-10.8	Panel 1	S	0	100	31	0
Imipenem IMI	EURL S-10.3	Panel 2	S	0	100	31	0
	EURL S-10.7	Panel 2	S	0	100	31	0
	EURL S-10.8	Panel 2	S	0	100	31	0
	<i>EURL S-10.4*</i>	<i>Panel 2</i>	<i>R</i>	<i>70</i>	<i>30</i>	<i>21</i>	<i>9</i>
Meropenem MER	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	S	0	100	31	0
	EURL S-10.3	Panel 1	S	0	100	31	0
	EURL S-10.3	Panel 2	S	0	100	31	0
	EURL S-10.4	Panel 1	R	100	0	31	0
	EURL S-10.4	Panel 2	R	100	0	31	0
	EURL S-10.5	Panel 1	S	0	100	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 2	S	0	100	31	0
	EURL S-10.8	Panel 1	S	0	100	31	0
	EURL S-10.8	Panel 2	S	0	100	31	0
Nalidixic acid NAL	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	S	0	100	31	0
	EURL S-10.3	Panel 1	R	100	0	31	0
	EURL S-10.4	Panel 1	S	3	97	31	1
	EURL S-10.5	Panel 1	S	0	100	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	S	0	100	30	0
	EURL S-10.8	Panel 1	S	0	100	31	0
Sulfamethoxazole SMX	EURL S-10.1	Panel 1	S	0	100	30	0
	EURL S-10.2	Panel 1	R	100	0	31	0
	EURL S-10.3	Panel 1	S	3	97	31	1
	EURL S-10.4	Panel 1	S	0	100	29	0
	EURL S-10.5	Panel 1	R	100	0	31	0
	EURL S-10.6	Panel 1	S	7	93	30	2
	EURL S-10.7	Panel 1	R	100	0	31	0
	EURL S-10.8	Panel 1	S	3	97	30	1
Tetracycline TET	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	R	97	3	31	1
	EURL S-10.3	Panel 1	R	100	0	31	0
	EURL S-10.4	Panel 1	S	0	100	30	0
	EURL S-10.5	Panel 1	R	100	0	31	0
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	R	100	0	31	0
	EURL S-10.8	Panel 1	S	0	100	31	0

Antimicrobial	Strain	Panel	Expected	% R	% S	No. correct	No. incorrect
Tigecycline TGC	EURL S-10.1	Panel 1	S	0	100	29	0
	EURL S-10.2	Panel 1	S	0	100	29	0
	EURL S-10.3	Panel 1	S	0	100	29	0
	EURL S-10.4	Panel 1	S	0	100	29	0
	EURL S-10.5	Panel 1	R	85	15	26	4
	EURL S-10.6	Panel 1	S	0	100	28	0
	EURL S-10.7	Panel 1	S	4	96	28	1
	EURL S-10.8	Panel 1	S	0	100	28	0
Trimethoprim TMP	EURL S-10.1	Panel 1	S	0	100	31	0
	EURL S-10.2	Panel 1	S	0	100	31	0
	EURL S-10.3	Panel 1	S	0	100	31	0
	EURL S-10.4	Panel 1	S	0	100	31	0
	EURL S-10.5	Panel 1	R	97	3	31	1
	EURL S-10.6	Panel 1	S	0	100	31	0
	EURL S-10.7	Panel 1	R	100	0	31	0
	EURL S-10.8	Panel 1	S	3	97	30	1

\*Strain/antimicrobial-combination excluded from the evaluation

*Campylobacter* - expected and obtained interpretation

Antimicrobial	Strain	Expected	% R	% S	No. correct	No. incorrect
Ciprofloxacin, CIP	EURL C-10.1	R	100	0	31	0
	EURL C-10.2	S	0	100	31	0
	EURL C-10.3	S	3	97	30	1
	EURL C-10.4	S	0	100	31	0
	EURL C-10.5	R	100	0	31	0
	EURL C-10.6	R	100	0	31	0
	EURL C-10.7	R	100	0	31	0
	EURL C-10.8	R	100	0	31	0
Erythromycin, ERY	EURL C-10.1	R	100	0	31	0
	EURL C-10.2	S	0	100	31	0
	EURL C-10.3	S	3	97	30	1
	EURL C-10.4	R	100	0	31	0
	EURL C-10.5	S	0	100	31	0
	EURL C-10.6	R	100	0	31	0
	EURL C-10.7	S	0	100	31	0
	EURL C-10.8	R	100	0	31	0
Gentamicin, GEN	EURL C-10.1	S	6	94	29	2
	EURL C-10.2	S	0	100	31	0
	EURL C-10.3	S	3	97	30	1
	EURL C-10.4	S	0	100	31	0
	EURL C-10.5	S	3	97	30	1
	EURL C-10.6	S	3	97	30	1
	EURL C-10.7	S	3	97	30	1
	EURL C-10.8	R	100	0	31	0
Nalidixic acid, NAL	EURL C-10.1	R	97	3	29	1
	EURL C-10.2	S	3	97	30	1
	EURL C-10.3	S	0	100	31	0
	EURL C-10.4	S	0	100	31	0
	EURL C-10.5	R	97	3	30	1
	EURL C-10.6	R	97	3	30	1
	EURL C-10.7	R	97	3	30	1
	EURL C-10.8	R	100	0	31	0
Streptomycin, STR	EURL C-10.1	S	3	97	30	1
	EURL C-10.2	R	94	6	29	2
	EURL C-10.3	S	3	97	30	1
	EURL C-10.4	S	0	100	31	0
	EURL C-10.5	S	0	100	31	0
	EURL C-10.6	S	3	97	29	1
	EURL C-10.7	S	6	94	29	2
	EURL C-10.8	R	100	0	31	0
Tetracycline, TET	EURL C-10.1	S	3	97	30	1
	EURL C-10.2	S	0	100	31	0
	EURL C-10.3	S	3	97	30	1
	EURL C-10.4*	S	26	74	31	8
	EURL C-10.5	R	97	3	30	1
	EURL C-10.6	R	100	0	31	0
	EURL C-10.7	S	0	100	31	0
	EURL C-10.8	R	100	0	31	0

\*Strain/antimicrobial-combination excluded from the evaluation

Deviations - *Salmonella*

Lab no.	Strain	Panel	Antimicrobial	Obtained MIC value	Obtained interpretation	Expected MIC-value	Expected interpretation
2	EURL S-10.1	1	Colistin COL	8	S	8	R
4	EURL S-10.3	1	Ciprofloxacin CIP	= 0.5	S	= 0.5	R
4	EURL S-10.5	1	Tigecycline TGC	2	S	2	R
4	EURL S-10.5	1	Trimethoprim TMP	> 32	S	> 32	R
6	EURL S-10.4	1	Ampicillin AMP	> 64	S	> 64	R
6	EURL S-10.7	1	Ampicillin AMP	> 64	S	> 64	R
18	EURL S-10.6	1	Sulfamethoxazole SMX	> 1024	R	64	S
19	EURL S-10.8	1	Sulfamethoxazole SMX	1024	R	32	S
19	EURL S-10.8	1	Trimethoprim TMP	> 32	R	= 0.5	S
20	EURL S-10.3	1	Gentamicin GEN	<= 0.5	R	1	S
21	EURL S-10.4	1	Nalidixic acid NAL	128	R	<= 4	S
22	EURL S-10.3	1	Colistin COL	8	R	<= 1	S
26	EURL S-10.5	1	Tigecycline TGC	1	S	2	R
26	EURL S-10.7	1	Colistin COL	8	R	<= 1	S
26	EURL S-10.7	1	Tigecycline TGC	<= 0.25	R	1	S
29	EURL S-10.3	1	Ceftazidime TAZ	1	R	1	S
29	EURL S-10.3	2	Ceftazidime TAZ	1	R	1	S
36	EURL S-10.5	1	Tigecycline TGC	<= 0.25	S	2	R
39	EURL S-10.3	1	Ciprofloxacin CIP	= 0.5	S	= 0.5	R
39	EURL S-10.4		ESBL test conclusion	Unusual phenotype		Presumptive carbapenemase	
39	EURL S-10.5	1	Chloramphenicol CHL	64	S	64	R
40	EURL S-10.1	1	Colistin COL	2	S	8	R
40	EURL S-10.2	1	Tetracycline TET	<= 2	S	64	R
40	EURL S-10.3	1	Sulfamethoxazole SMX	> 1024	R	32	S
40	EURL S-10.5	1	Tigecycline TGC	1	S	2	R
40	EURL S-10.7	1	Ciprofloxacin CIP	= 0.06	S	= 0.25	R
45	EURL S-10.6	1	Sulfamethoxazole SMX	> 1024	R	64	S
56	EURL S-10.1	1	Colistin COL	2	S	8	R
56	EURL S-10.7	2	Cefoxitin FOX	16	R	8	S
56	EURL S-10.7		ESBL test conclusion	Presumptive ESBL + pAmpC		Presumptive ESBL	
60	EURL S-10.3		ESBL test conclusion	Unusual phenotype		Presumptive ESBL	



Deviations - *Campylobacter*

Lab no.	Strain	Antimicrobial	Obtained MIC value	Obtained interpretation	Expected MIC-value	Expected interpretation
2	EURL C-10.6	Streptomycin STR	16	R	2	S
18	EURL C-10.7	Streptomycin STR	8	R	4	S
19	EURL C-10.7	Streptomycin STR	8	R	4	S
36	EURL C-10.1	Gentamicin GEN	16	R	= 0.25	S
36	EURL C-10.1	Nalidixic acid NAL	4	S	> 64	R
36	EURL C-10.2	Nalidixic acid NAL	> 64	R	4	S
36	EURL C-10.2	Streptomycin STR	2	S	16	R
36	EURL C-10.5	Gentamicin GEN	> 16	R	= 0.25	S
36	EURL C-10.5	Nalidixic acid NAL	8	S	> 64	R
36	EURL C-10.6	Gentamicin GEN	> 16	R	= 0.5	S
36	EURL C-10.6	Nalidixic acid NAL	8	S	> 64	R
36	EURL C-10.7	Gentamicin GEN	16	R	= 0.5	S
36	EURL C-10.7	Nalidixic acid NAL	8	S	> 64	R
39	EURL C-10.2	Streptomycin STR	<= 0.5	S	16	R
39	EURL C-10.3	Ciprofloxacin CIP	> 8	R	<= 0.12	S
39	EURL C-10.3	Erythromycin ERY	> 64	R	<= 1	S
39	EURL C-10.3	Gentamicin GEN	> 64	R	= 0.25	S
39	EURL C-10.3	Streptomycin STR	8	R	1	S
39	EURL C-10.3	Tetracycline TET	> 16	R	<= 0.5	S
40	EURL C-10.5	Tetracycline TET	1	S	> 64	R
42	EURL C-10.1	Gentamicin GEN	= 0.5	R	= 0.25	S
42	EURL C-10.1	Streptomycin STR	1	R	1	S
42	EURL C-10.1	Tetracycline TET	<= 0.5	R	<= 0.5	S

## Genotypic characterization (optional); obtained results

Labno	Strain	Genotype	Gene number	Not detected	Method	Reference	Primer 5 3	Primer 3 5
4	EURL-S10.3	CTX	M-9		PCR (published)	Hasman 2005. JAC 56:115-121	-	-
4	EURL-S10.3	OXA	-	X	PCR (published)	Hasman 2005. JAC 56:115-121	-	-
4	EURL-S10.3	SHV	-	X	PCR (published)	Arlet1997.FEMS ML. 152:163-7	-	-
4	EURL-S10.3	TEM	-1		PCR (published)	Olesen 2004. MDR. 10:334-340	-	-
4	EURL-S10.4	CTX	-	X	PCR (published)	Hasman 2005. JAC 56:115-121	-	-
4	EURL-S10.4	IMP	-	X	PCR (published)	Poirel. 2011 Diag Micro Infect Dis.70(1):119-23.	-	-
4	EURL-S10.4	KPC	-	X	PCR (published)	Poirel. 2011 Diag Micro Infect Dis.70(1):119-23.	-	-
4	EURL-S10.4	OXA	-48		PCR (published)	Poirel. 2011 Diag Micro Infect Dis.70(1):119-23.	-	-
4	EURL-S10.4	SHV	-	X	PCR (published)	Arlet1997.FEMS ML. 152:163-7	-	-
4	EURL-S10.4	TEM	-1		PCR (published)	Olesen 2004. MDR. 10:334-340	-	-
4	EURL-S10.4	VIM	-	X	PCR (published)	Poirel. 2011 Diag Micro Infect Dis.70(1):119-23.	-	-
4	EURL-S10.7	CTX	M-15		PCR (published)	Hasman 2005. JAC 56:115-121	-	-
4	EURL-S10.7	OXA	-	X	PCR (published)	Hasman 2005. JAC 56:115-121	-	-
4	EURL-S10.7	SHV	-12		PCR (published)	Arlet1997.FEMS ML. 152:163-7	-	-
4	EURL-S10.7	TEM	-1		PCR (published)	Olesen 2004. MDR. 10:334-340	-	-
4	EURL-S10.8	CTX	-	X	-	Hasman 2005. JAC 56:115-121	-	-
4	EURL-S10.8	OXA	-	X	-	Hasman 2005. JAC 56:115-121	-	-
4	EURL-S10.8	SHV	-	X	-	Arlet1997.FEMS ML. 152:163-7	-	-
4	EURL-S10.8	TEM	-52		-	Olesen 2004. MDR. 10:334-340	-	-
9	EURL-S10.3	CTX	M-9		PCR (published)	JAC 2010;65;490-495	-	-
9	EURL-S10.3	TEM	-1		PCR (published)	JAC 2010;65;490-495	-	-
9	EURL-S10.4	OXA	-48		PCR (published)	JAC 2010;65;490-495	-	-
9	EURL-S10.7	CTX	M-15		PCR (published)	JAC 2010;65;490-495	-	-
9	EURL-S10.7	SHV	-12		PCR (published)	JAC 2010;65;490-495	-	-
9	EURL-S10.7	TEM	-1		PCR (published)	JAC 2010;65;490-495	-	-
9	EURL-S10.8	TEM	-52		PCR (published)	JAC 2010;65;490-495	-	-
17	EURL-S10.3	CTX	M-9		PCR (published)	Paauw et al. (2006)	TGGTGACAAAGAGAGTGCAACG	ATGCCGAAGGGCTGTGA
17	EURL-S10.3	TEM	-1		PCR (published)	Guerra et al. (2001)	TTGGGTGCACGAGTGGGT	GCTTCCCGCAACAATTA
17	EURL-S10.4	OXA	-48		PCR (published)	Guerra et al. (2000)	AGCAGCGCCAGTGCATCA	GGAAACTTGGGTCTGAAT
17	EURL-S10.4	TEM	-1		PCR (published)	Guerra et al. (2001)	TTGGGTGCACGAGTGGGT	GCTTCCCGCAACAATTA
17	EURL-S10.7	CTX	M-15		PCR (published)	Carrattoli et al 2008	CCCATGGTTAAAAATCACTGC	CTTAGACGGCAAAAGCGCTG
17	EURL-S10.7	SHV	-12		PCR (published)	Weill et al., (2004)	TTATCTCCCTGTTAGCCACC	CCGAGCGAAATCAGCAAATC
17	EURL-S10.7	TEM	-1		PCR (published)	Guerra et al. (2001)	TTGGGTGCACGAGTGGGT	GCTTCCCGCAACAATTA
17	EURL-S10.8	TEM	-52		PCR (published)	Guerra et al. (2001)	TTGGGTGCACGAGTGGGT	GCTTCCCGCAACAATTA
21	EURL-S10.3	CTX	M-4		PCR (published)	doi: 10.1093/jac/dki412	attgaaagcgttcatcacc	caaagagagtgaacggatg
21	EURL-S10.4	OXA	-48		PCR (published)	doi: 10.3201/eid1710.110655	CATCAAGTTCAACCAACCG	GCGTGGTTAAGGATGAACAC
21	EURL-S10.7	CTX	M-1		PCR (published)	doi:10.1128/JCM.42.12.5715-5721.2004	agccgacgacgtaataca	gacgatgtactggctgagc
21	EURL-S10.7	SHV	-		PCR (published)	doi:10.1016/j.diagmicrobio.2006.04.016	tcctccgatgccgccagatca	gccgggttattctattgt
21	EURL-S10.8	TEM	-52		PCR (in-house)	-	ttaccaatgcttaataca	atgagtattcaacatttccg
25	EURL-S10.3	CTX	M-9		PCR (in-house)	-	-	-
25	EURL-S10.3	TEM	-		PCR (in-house)	-	-	-
25	EURL-S10.4	OXA	-48		PCR (in-house)	-	-	-
25	EURL-S10.4	TEM	-		PCR (in-house)	-	-	-
25	EURL-S10.7	CTX	M-15		PCR (in-house)	-	-	-
25	EURL-S10.7	TEM	-		PCR (in-house)	-	-	-
25	EURL-S10.8	TEM	-52		PCR (in-house)	-	-	-
32	EURL-S10.3	ACC	-	X	PCR (published)	(Hasman et al. 2005)	-	-
32	EURL-S10.3	ACT	-	X	PCR (published)	Voets et al. 2011	-	-
32	EURL-S10.3	CMY	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.3	CTX	M-9		PCR (published)	PediatrInfectDisJ28:814-818	-	-

Labno	Strain	Genotype	Gene number	Not detected	Method	Reference	Primer 5 3	Primer 3 5
32	EURL-S10.3	DHA	-	X	PCR (published)	(Gonzalez-Sanz et al.2009)	-	-
32	EURL-S10.3	FOX	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.3	IMP	-	X	PCR (published)	L. Poirel et al 2011	-	-
32	EURL-S10.3	KPC	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.3	MOX	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.3	NDM	-	X	PCR (published)	L. Poirel et al 2011	-	-
32	EURL-S10.3	OXA	-48	X	PCR (published)	Voets et al 2011	-	-
32	EURL-S10.3	OXA	-10	X	PCR (published)	Voets et al. 2011	-	-
32	EURL-S10.3	OXA	-	X	PCR (published)	J. Antimic.Chemothe(2009) 64	-	-
32	EURL-S10.3	SHV	-	X	PCR (published)	FEMSMicrobiolLett1997152:1637	-	-
32	EURL-S10.3	TEM*	-	X	PCR (published)	AntimicrAgentsChemotherap2009	-	-
32	EURL-S10.3	VEB	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.3	VIM	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.4	ACC	-	X	PCR (published)	Hasman et al. 2005)	-	-
32	EURL-S10.4	ACT	-	X	PCR (published)	Voets et al. 2011	-	-
32	EURL-S10.4	CMY	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.4	CTX	-	X	PCR (published)	PediatrInfectDisJ28:814-818	-	-
32	EURL-S10.4	DHA	-	X	PCR (published)	(Gonzalez-Sanz et al.2009)	-	-
32	EURL-S10.4	FOX	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.4	IMP	-	X	PCR (published)	L. Poirel et al 2011	-	-
32	EURL-S10.4	KPC	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.4	MOX	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.4	NDM	-	X	PCR (published)	L. Poirel et al 2011	-	-
32	EURL-S10.4	OXA	-48		PCR (published)	Voets et al. 2011	-	-
32	EURL-S10.4	OXA	-10	X	PCR (published)	J. Antimic.Chemothe(2009) 64	-	-
32	EURL-S10.4	OXA	-10	X	PCR (published)	Voets et al. 2011	-	-
32	EURL-S10.4	SHV	-	X	PCR (published)	FEMSMicrobiolLett1997152:1637	-	-
32	EURL-S10.4	TEM	-	X	PCR (published)	AntimicrAgentsChemotherap2009	-	-
32	EURL-S10.4	VEB	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.4	VIM	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.7	ACC	-	X	PCR (published)	(Hasman et al. 2005)	-	-
32	EURL-S10.7	ACT	-	X	PCR (published)	Voets et al. 2011	-	-
32	EURL-S10.7	CMY	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.7	CTX	M-15		PCR (published)	PediatrInfectDisJ28:814-818	-	-
32	EURL-S10.7	DHA	-	X	PCR (published)	(Gonzalez-Sanz et al.2009)	-	-
32	EURL-S10.7	FOX	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.7	IMP	-	X	PCR (published)	L. Poirel et al 2011	-	-
32	EURL-S10.7	KPC	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.7	MOX	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.7	NDM	-	X	PCR (published)	L. Poirel et al 2011	-	-
32	EURL-S10.7	OXA	-48	X	PCR (published)	Voets et al 2011	-	-
32	EURL-S10.7	OXA	-10	X	PCR (published)	Voets et al 2011	-	-
32	EURL-S10.7	OXA	-	X	PCR (published)	J. Antimic.Chemothe(2009) 64	-	-
32	EURL-S10.7	SHV	-12		PCR (published)	FEMSMicrobiolLett1997152:1637	-	-
32	EURL-S10.7	TEM	-1		PCR (published)	AntimicrAgentsChemotherap2009	-	-
32	EURL-S10.7	VEB	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.7	VIM	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.8	ACC	-	X	PCR (published)	(Hasman et al. 2005)	-	-
32	EURL-S10.8	ACT	-	X	PCR (published)	Voets et al. 2011	-	-
32	EURL-S10.8	CMY	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.8	CTX	-	X	PCR (published)	PediatrInfectDisJ28:814-818	-	-
32	EURL-S10.8	DHA	-	X	PCR (published)	(Gonzalez-Sanz et al.2009)	-	-

Labno	Strain	Genotype	Gene number	Not detected	Method	Reference	Primer 5 3	Primer 3 5
32	EURL-S10.8	FOX	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.8	IMP	-	X	PCR (published)	L. Poirel et al 2011	-	-
32	EURL-S10.8	KPC	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.8	MOX	-	X	PCR (published)	A.AGE.CHEM2006.ClinMicr1999-37	-	-
32	EURL-S10.8	NDM	-	X	PCR (published)	L. Poirel et al 2011	-	-
32	EURL-S10.8	OXA	-48	X	PCR (published)	Voets et al 2011	-	-
32	EURL-S10.8	OXA	-10	X	PCR (published)	Voets et al 2011	-	-
32	EURL-S10.8	OXA	-	X	PCR (published)	J. Antimic.Chemothe(2009) 64	-	-
32	EURL-S10.8	SHV	-	X	PCR (published)	FEMSMicrobiolLett1997152:1637	-	-
32	EURL-S10.8	TEM	-52		PCR (published)	AntimicrAgentsChemotherap2009	-	-
32	EURL-S10.8	VEB	-	X	PCR (published)	Dallanne et al 2010	-	-
32	EURL-S10.8	VIM	-	X	PCR (published)	Dallanne et al 2010	-	-
33	EURL-S10.3	CTX	-		PCR (published)	Woodford,et al. (2006)	CAAAGAGARTGCAACGGATG	ATTGGAAGCGTTCATCACC
33	EURL-S10.3	TEM	-		PCR (published)	Fang, et al. (2008)	CGCCGCATACACTATTCTCAGAATGA	ACGCTCACC GGCTCCAGATTAT
33	EURL-S10.4	OXA	-48		PCR (published)	Poirel et al.(2011)	GCGTGGTTAAGGATGAACAC	CATCAAGTTCAACCAACCG
33	EURL-S10.4	TEM	-		PCR (published)	Fang,et al. (2008).	CGCCGCATACACTATTCTCAGAATGA	ACGCTCACC GGCTCCAGATTAT
33	EURL-S10.7	CTX	-		PCR (published)	Woodford,et al. (2006)	AAAAATCACTGCGYCA GTTC	AGCTTATTCATCGCCACGTT
33	EURL-S10.7	SHV	-		PCR (published)	Fang,et al. (2008).	CTTTATCGGCCCTCACTCAA	AGGTGCTCATCATGGGAAAG
33	EURL-S10.7	TEM	-		PCR (published)	Fang,et al. (2008)	CGCCGCATACACTATTCTCAGAATGA	ACGCTCACC GGCTCCAGATTAT
33	EURL-S10.8	TEM	-		PCR (published)	Fang,et al. (2008)	CGCCGCATACACTATTCTCAGAATGA	ACGCTCACC GGCTCCAGATTAT
36	EURL-S10.3	CTX	M-9		PCR (published)	Hasman et al. JAC. 2005 Jul;56(1):115-21	ATGTGCAGYACCAAGTAARGTKATGGC	TGGGTRAARTARGTSACCAGAAYCAGCGG
36	EURL-S10.3	SHV	-	X	PCR (published)	Briñas et al. AAC. 2002 Oct;46(10):3156-63	CACTCAAGGATGTATTGTG	TTAGCGTTGCCAGTGCTCG
36	EURL-S10.3	TEM	-1		PCR (published)	Briñas et al. AAC. 2002 Oct;46(10):3156-63	TTCTGAAGACGAAAGGGC	ACGCTCAGTGGAAACGAAAC
36	EURL-S10.4	CTX	-	X	PCR (published)	Hasman et al. JAC 2005 Jul;56(1):115-21	ATGTGCAGYACCAAGTAARGTKATGGC	TGGGTRAARTARGTSACCAGAAYCAGCGG
36	EURL-S10.4	IMP	-	X	PCR (published)	Ellington JAC (2007)59(2):321-322	GGAATAGAGTGGCTTAAYTCTC	CCAAACYACTASGTTATCT
36	EURL-S10.4	NDM	-	X	PCR (published)	Mushtaq et. al. JAC (2011)66(9):2002-2005	GGGCAGTCGCTTCCAACGGT	GTAGTGCTCAGTGTCCGCAT
36	EURL-S10.4	OXA	-48		PCR (published)	Poirel et al. DiaMicInfDis. 2011 May;70(1):119-23	GCGTGGTTAAGGATGAACAC	CATCAAGTTCAACCAACCG
36	EURL-S10.4	SHV	-	X	PCR (published)	Briñas et al. AAC 2002 Oct;46(10):3156-63	CACTCAAGGATGTATTGTG	TTAGCGTTGCCAGTGCTCG
36	EURL-S10.4	TEM	-1		PCR (published)	Briñas et al. AAC 2002 Oct;46(10):3156-63	TTCTGAAGACGAAAGGGC	ACGCTCAGTGGAAACGAAAC
36	EURL-S10.4	VIM	-	X	PCR (published)	Ellington JAC (2007)59(2):321-322	GATGGTGTGGTGCATATA	CGAATGCGCAGCACCAG
36	EURL-S10.7	CTX	M-15		PCR (published)	Hasman et al. JAC. 2005 Jul;56(1):115-21.	ATGTGCAGYACCAAGTAARGTKATGGC	TGGGTRAARTARGTSACCAGAAYCAGCGG
36	EURL-S10.7	SHV	-12		PCR (published)	Briñas et al. AAC. 2002 Oct;46(10):3156-63.	CACTCAAGGATGTATTGTG	TTAGCGTTGCCAGTGCTCG
36	EURL-S10.7	TEM	-1		PCR (published)	Briñas et al. AAC. 2002 Oct;46(10):3156-63.	TTCTGAAGACGAAAGGGC	ACGCTCAGTGGAAACGAAAC
36	EURL-S10.8	CTX	-	X	PCR (published)	Hasman et al. JAC. 2005 Jul;56(1):115-21.	ATGTGCAGYACCAAGTAARGTKATGGC	TGGGTRAARTARGTSACCAGAAYCAGCGG
36	EURL-S10.8	SHV	-	X	PCR (published)	Briñas et al. AAC. 2002 Oct;46(10):3156-63.	CACTCAAGGATGTATTGTG	TTAGCGTTGCCAGTGCTCG
36	EURL-S10.8	TEM	-52		PCR (published)	Briñas et al. AAC. 2002 Oct;46(10):3156-63.	TTCTTGAAGACGAAAGGGC	ACGCTCAGTGGAAACGAAAC
59	EURL-S10.3	CTX	M-9		Whole genome sequenced	-	-	-
59	EURL-S10.3	TEM	-1		Whole genome sequenced	-	-	-
59	EURL-S10.4	OXA	-48		Whole genome sequenced	-	-	-
59	EURL-S10.4	TEM	-1		Whole genome sequenced	-	-	-
59	EURL-S10.7	CTX	M-15		Whole genome sequenced	-	-	-
59	EURL-S10.7	SHV	-12		Whole genome sequenced	-	-	-
59	EURL-S10.7	TEM	-1		Whole genome sequenced	-	-	-
59	EURL-S10.8	TEM	-52		Whole genome sequenced	-	-	-

## Legend:

Fields shaded grey indicate that the gene was expected

Genes in bold and white font, were detected but not expected

\*TEM-1 does not confer ESBL-production and is as such not included as an expected result. TEM-1 was, however, present in S-10.3, S-10.4 and S-10.7

## Genotypic characterization (optional); comments by participants

Labno	Strain	Comment
4	S-9.5	TEM-1 detected
4	S-9.5	TEM-1 detected
17	S-9.3	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. CTX-M9 primer is group specific, validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.3	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. CTX-M9 primer is group specific, validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.3	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. CTX-M9 primer is group specific, validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.4	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. Validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.4	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. Validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.4	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. Validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.7	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. All used primers are specific for a respective group of genes, validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.7	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. All used primers are specific for a respective group of genes, validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.7	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. All used primers are specific for a respective group of genes, validation was done via sequencing TEM-1 was additionally found, it's just added for additional information, not as ESBL-gene
17	S-9.8	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. TEM primer is group specific, validation was done via sequencing
17	S-9.8	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. TEM primer is group specific, validation was done via sequencing
17	S-9.8	All primers were 5'->3'. The Primer in 3'->5' box is the reverse complement to the original primer. TEM primer is group specific, validation was done via sequencing
21	S-9.7	qnrB positive.
21	S-9.7	qnrB positive.
32	S-9.3	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.3	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.3	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.4	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.4	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.4	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.7	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.7	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.7	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.8	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.8	SPM Gene tested for/not detected (L. Poirel et al 2011)
32	S-9.8	SPM Gene tested for/not detected (L. Poirel et al 2011)
59	S-9.3	genenumber for genotype TEM: TEM-1B
59	S-9.4	genenumber for genotype TEM: TEM-1D
59	S-9.7	genenumber for genotype TEM: TEM-1B
59	S-9.8	genenumber for genotype TEM: TEM-52B

National Food Institute  
Technical University of Denmark  
Mørkhøj Bygade 19  
2860 Søborg

Tel. 35 88 70 00  
Fax 35 88 70 01

[www.food.dtu.dk](http://www.food.dtu.dk)

ISBN: 978-87-93109-86-5