



## Installation Qualification, Performance Qualification and Routine Control

**Miller, Arne**

*Publication date:*  
2011

*Document Version*  
Publisher's PDF, also known as Version of record

[Link back to DTU Orbit](#)

*Citation (APA):*  
Miller, A. (Author). (2011). Installation Qualification, Performance Qualification and Routine Control. Sound/Visual production (digital)

---

### General rights

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from the public portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain
- You may freely distribute the URL identifying the publication in the public portal

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.



*Cameron Quality Consultants*

## **Validation and routine control of sterilization by irradiation**

RISØ

Arne Miller  
Risø High Dose Reference Laboratory  
Risø DTU  
Technical University of Denmark  
DK-4000 Roskilde  
Denmark





*Cameron Quality Consultants*

## **Content of presentation:**

- **Installation Qualification - IQ**
- **Operational Qualification - OQ**
- **Performance Qualification - PQ**
- **Process Control and Monitoring**

RISO





*Cameron Quality Consultants*

## **Installation qualification**

*undertaken to demonstrate that the sterilization equipment has been supplied and installed in accordance with specification.*

Whether or not data are “in accordance with their specification” depends on agreement

Measurements often the same as for Operational Qualification

Responsibility: Facility Operator

ISO 11137-1, sect 9.1

RISQ





*Cameron Quality Consultants*

## **Operational qualification**

- *- carried out by irradiating appropriate test material to demonstrate the capability of the equipment to deliver the sterilization process that has been defined.*
- *- provides baseline data to show consistent operation of the facility*

**Responsibility: Facility Operator**

**ISO 11137-1, sect 9.2**

RISO





*Cameron Quality Consultants*

## Operational qualification Gamma

- dose distribution in reference product(s)
- dose as function of dwell time
- process interruption
- mix of products

RISO






Cameron Quality Consultants

## OQ gamma cont..

### Gamma Facility Layout

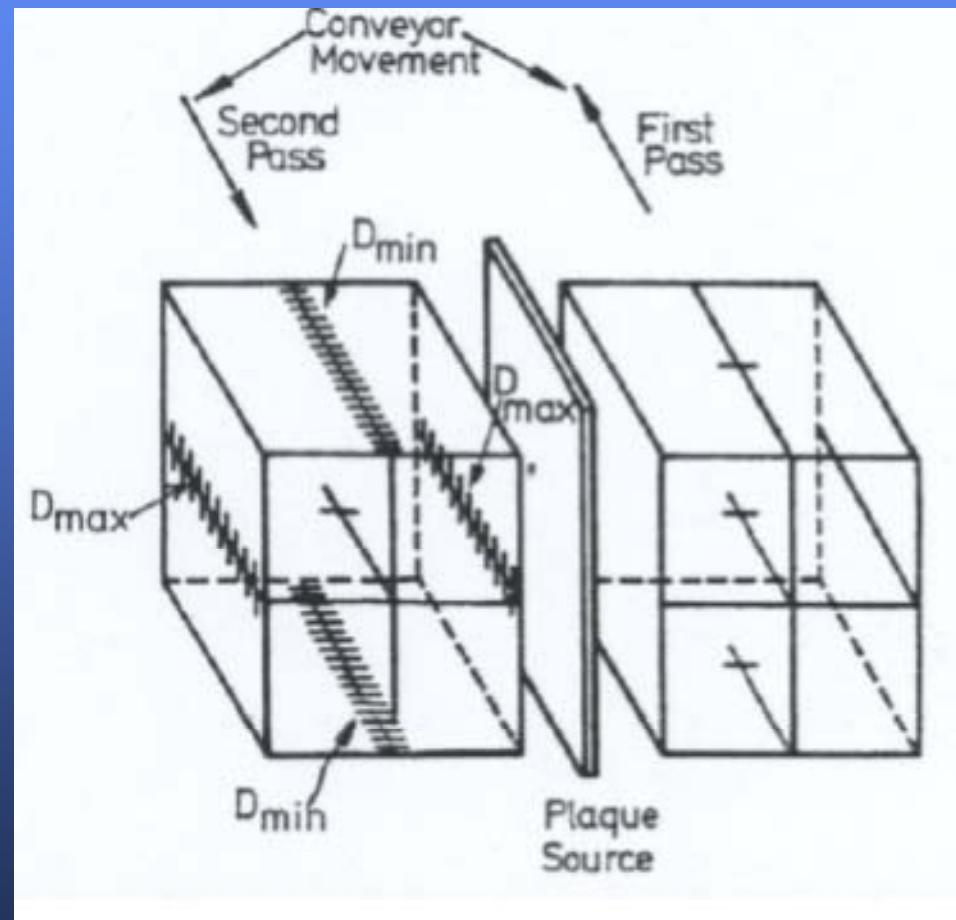

RISO










Cameron Quality Consultants

## OO Gamma cont..

Dose distribution  
in reference product

Isodose curves

RISO

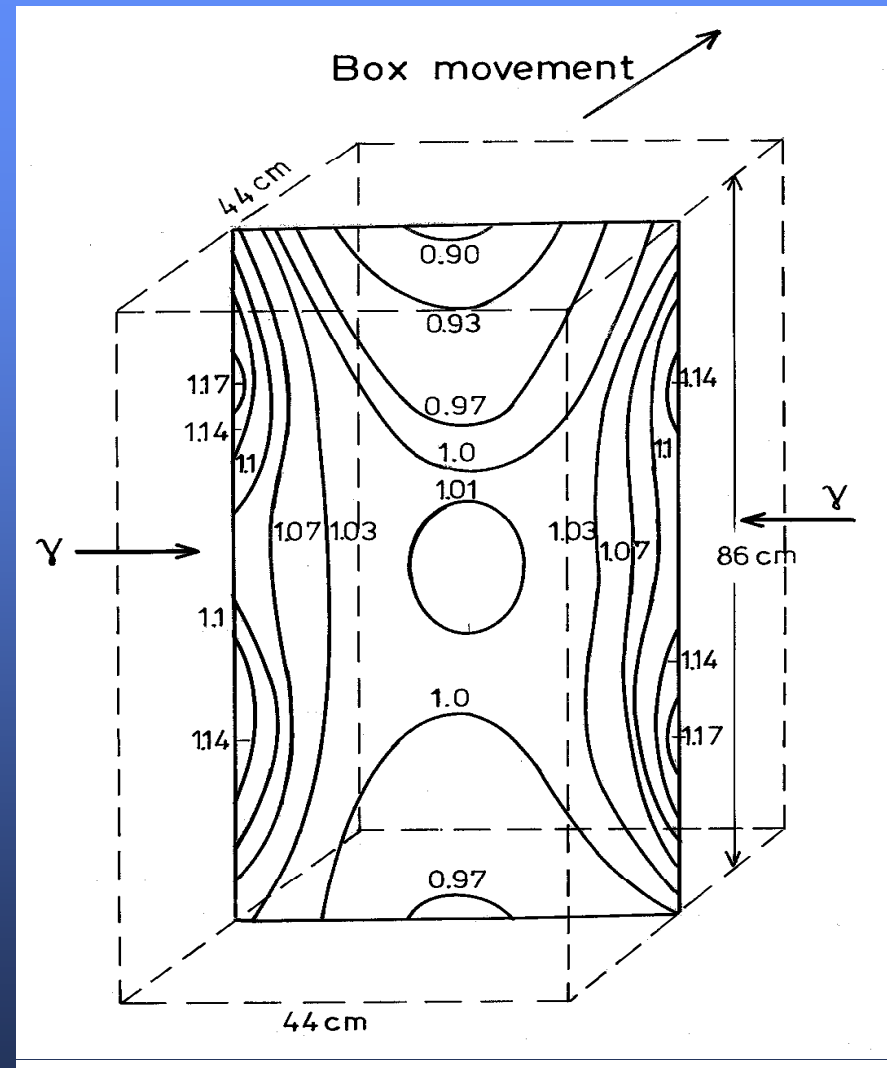



*Cameron Quality Consultants*

## **OQ gamma cont..**

### **Dose mapping**

- **Dose maps must be made with fully loaded irradiation chamber**
- **Range of densities of reference product must cover extremes of products to be irradiated**
- **Additional measurement:**
  - **effect of process interruption**

RISO





*Cameron Quality Consultants*

## **OQ gamma cont..**

### **Dose mapping**

- Dose mapping shall be carried out on a sufficient number of irradiation containers to allow determination of the distribution and variability of dose between containers
- Dose mapping must be repeated when significant changes of the irradiation plant are made, *e.g.* source change

RISO





*Cameron Quality Consultants*

## Operational qualification Electron beam

- Characteristic parameters to be measured
  - dose distribution reference product
  - scan width
  - energy
  - dose as function of speed, current, beam width
  - beam spot
  - process interruption

RISO






Cameron Quality Consultants

## OQ E-beam cont..

### Placement of dosimeters in reference product

RISO

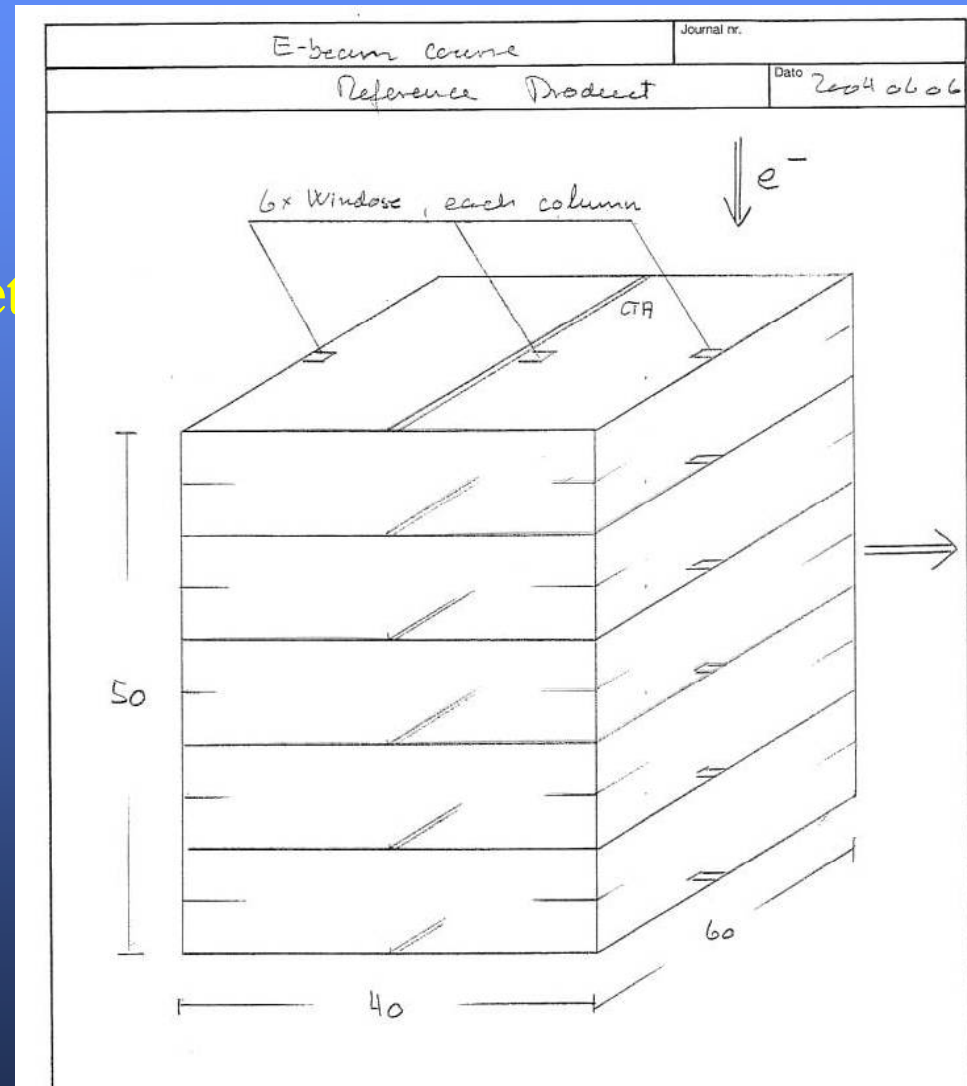




*Cameron Quality Consultants*

## OQ E-beam cont..

### Dose distribution in reference product

RISO

Limits for  
acceptable  
variations can be  
defined.

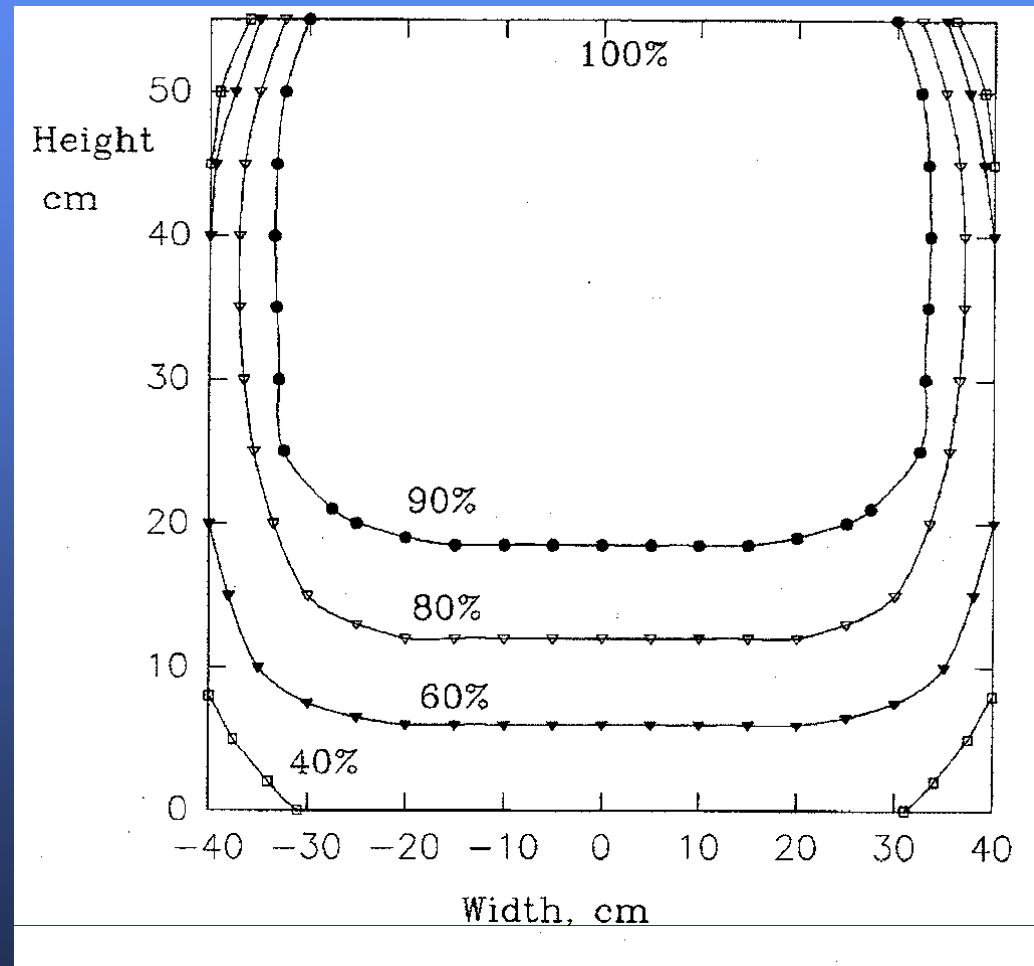




*Cameron Quality Consultants*

## OQ E-beam cont..

### Scan width (=extended beam width)

RISO

Limits for acceptable variations can be defined.

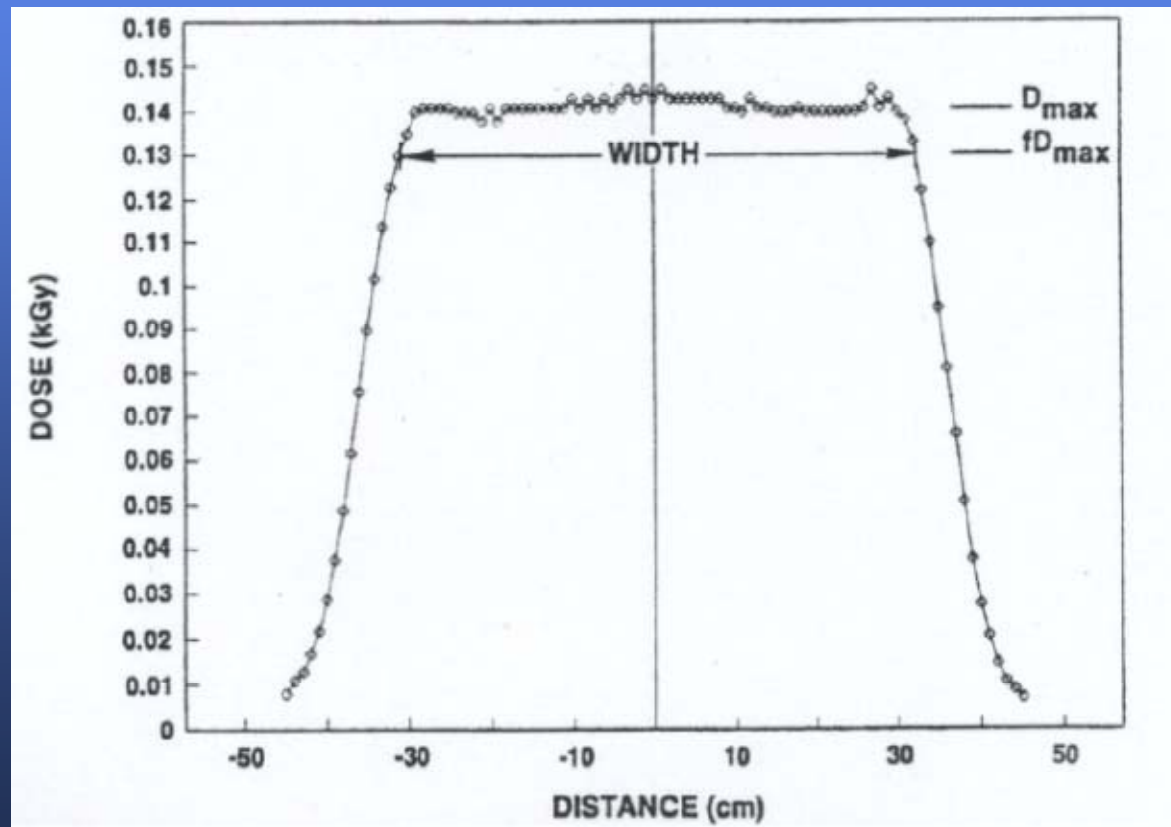




*Cameron Quality Consultants*

## OO E-beam cont..

### Energy

### Wedge and stack for energy measurement

RISO

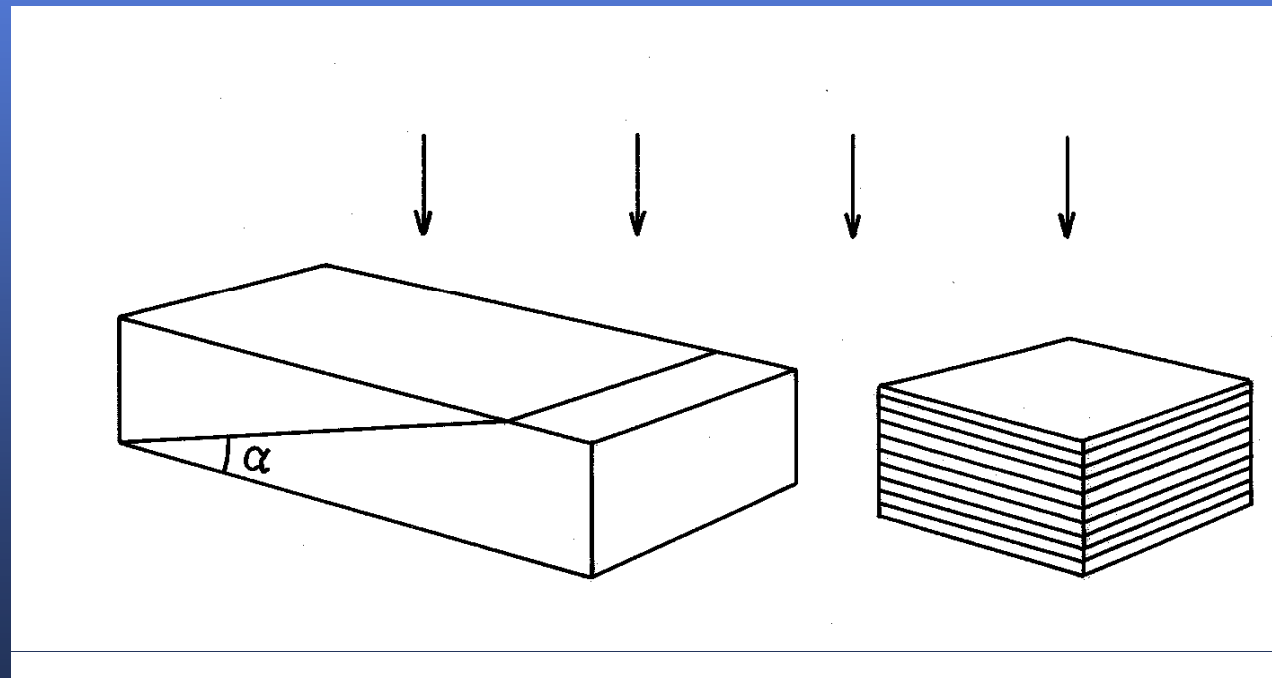






*Cameron Quality Consultants*

## OO E-beam cont..

Typical depth-dose curve

Energy – range relationships

RISO

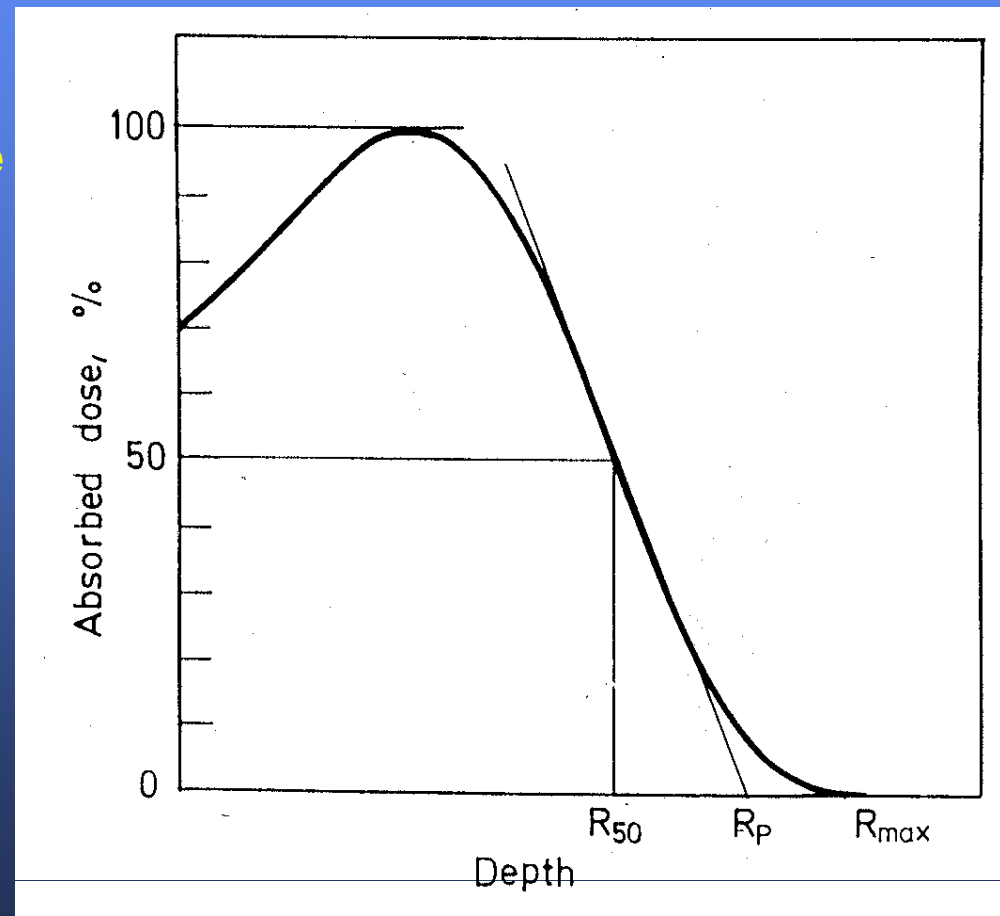




*Cameron Quality Consultants*

## OO E-beam cont..

Dose as a function of

- beam current (I)
- conv speed (V)
- beam width (Sw)

RISO

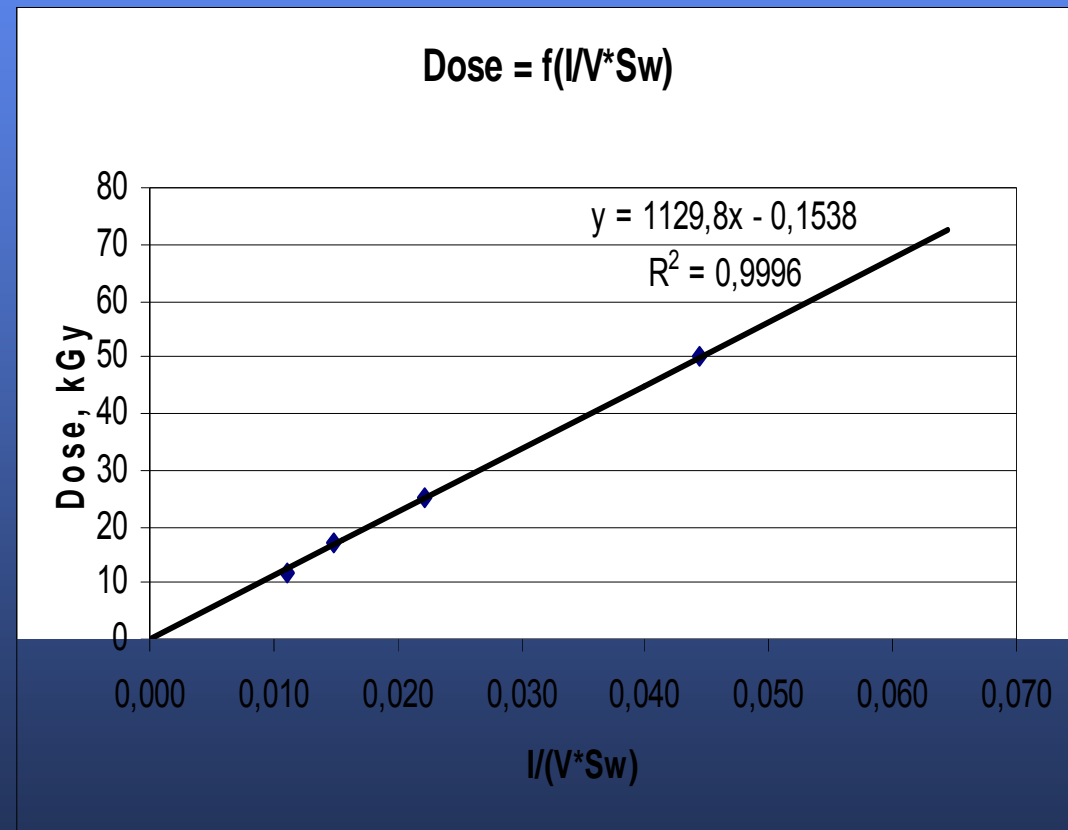




*Cameron Quality Consultants*

**OQ E-beam cont..**

**Beam spot pattern**

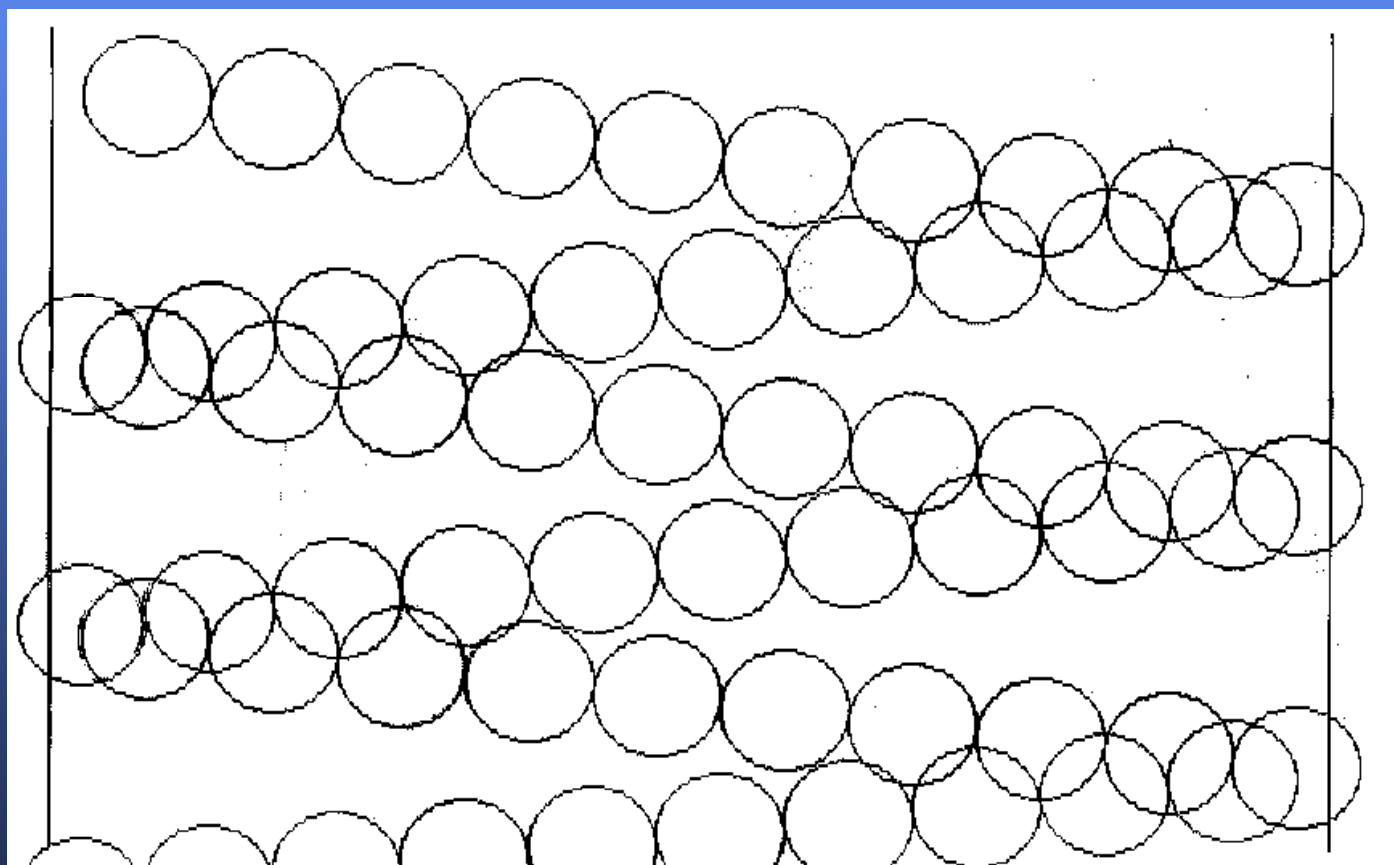

RISO






*Cameron Quality Consultants*

**OQ E-beam cont..**

**Process interruption**

RISO

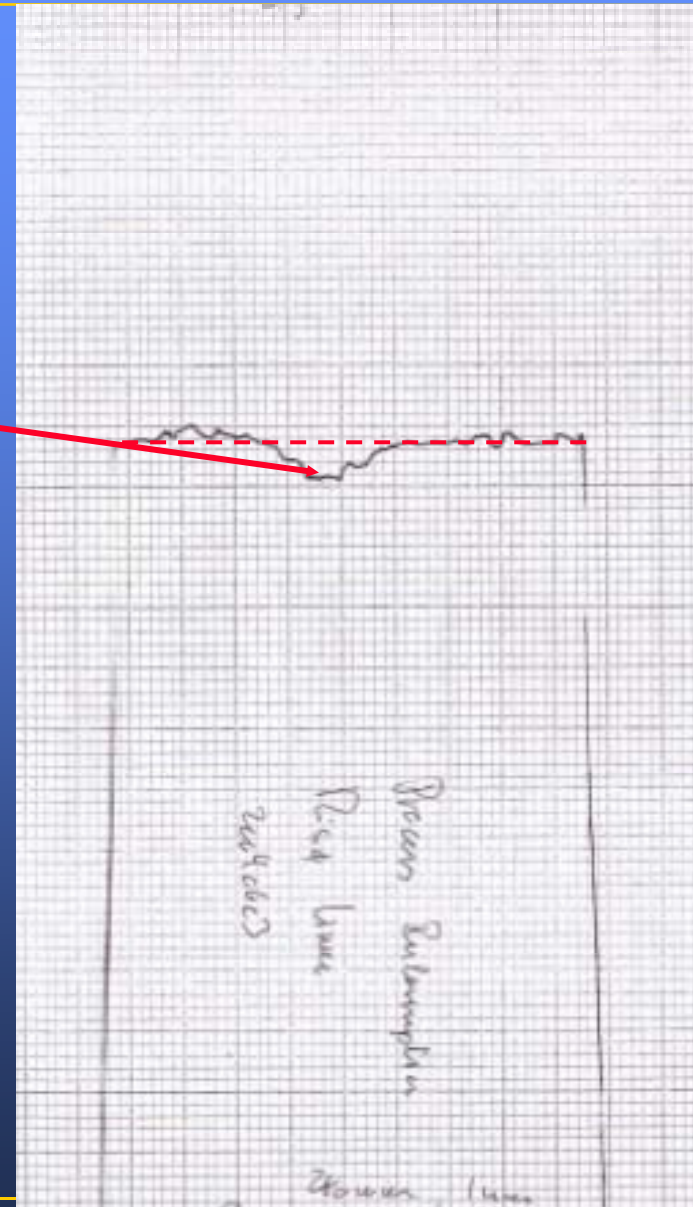



*Cameron Quality Consultants*

## **OQ E-beam cont..**

- E-beam OQ should be repeated when significant changes of the irradiation plant are made that may influence dose or dose distribution.
- It is recommended to repeat OQ at regular intervals, e.g. annually in order to demonstrate consistent operation.

RISO





*Cameron Quality Consultants*

## **Performance qualification**

Concerns dose mapping of real product

to identify the location and magnitude of minimum and maximum doses

and

to determine the relationship between the min and max doses and the routine dose

**Responsibility: Primary Manufacturer**

ISO 11137-1, sect 9.3

RISO





*Cameron Quality Consultants*

## **Performance qualification ..cont..**

Purpose of dose mapping

- to demonstrate that:

- $D(\text{min}) > D(\text{Steril})$  and  $D(\text{max}) < D(\text{acceptable})$
- $D(\text{min})$  determined by sterilization requirements
- $D(\text{max})$  determined by radiation-induced changes in product

RISO





*Cameron Quality Consultants*

## **Performance qualification ..cont..**

### **Strategies for dose mapping**

**- based on:**

- OQ measurements
- Inhomogeneous product distribution, orientation, voids, interfaces.
- Experience
  
- Processing categories

### **Limitations of dose mapping**

- Ability to measure “the true” value of D(min) or D(max)
- Location of dosimeter
- Resolution of dosimeter system

RISO







*Cameron Quality Consultants*

## Performance qualification ..cont..

Examples from real life .....

RISO





*Cameron Quality Consultants*

## Dose map examples



RISO





*Cameron Quality Consultants*

## **Performance qualification ..cont..**

### Measurement Uncertainty

*Dose mapping shall be carried out using representative irradiation containers sufficient in number to determine the variability of dose between containers*

### Practical approach:

Either dose map in detail 3-10 containers

or

- dose map in detail one container
- repeat measurement of D (min) and D (max) in several (e.g. 10) containers
- evaluate measurement uncertainty taking into account all uncertainty components.

RISO





*Cameron Quality Consultants*

## Performance qualification ..cont..

### Setting Process Limits

#### Select process parameters

- gamma: e.g. timer setting
- e-beam: e.g. conveyor speed
  
- so that  $D(\min)$  will exceed  $D(\text{sterile})$  and
- so that  $D(\max)$  will not exceed  $D(\text{acceptable})$

Select process limits based on the measured uncertainty.

#### Practical approach

- select process parameters so that  $D(\min)$  is two standard deviations higher than  $D(\text{steril})$

RISO





*Cameron Quality Consultants*


## Performance qualification ..cont..

### Selection of uncertainty limits for process control

RISØ



Limit based on  
2 standard deviations

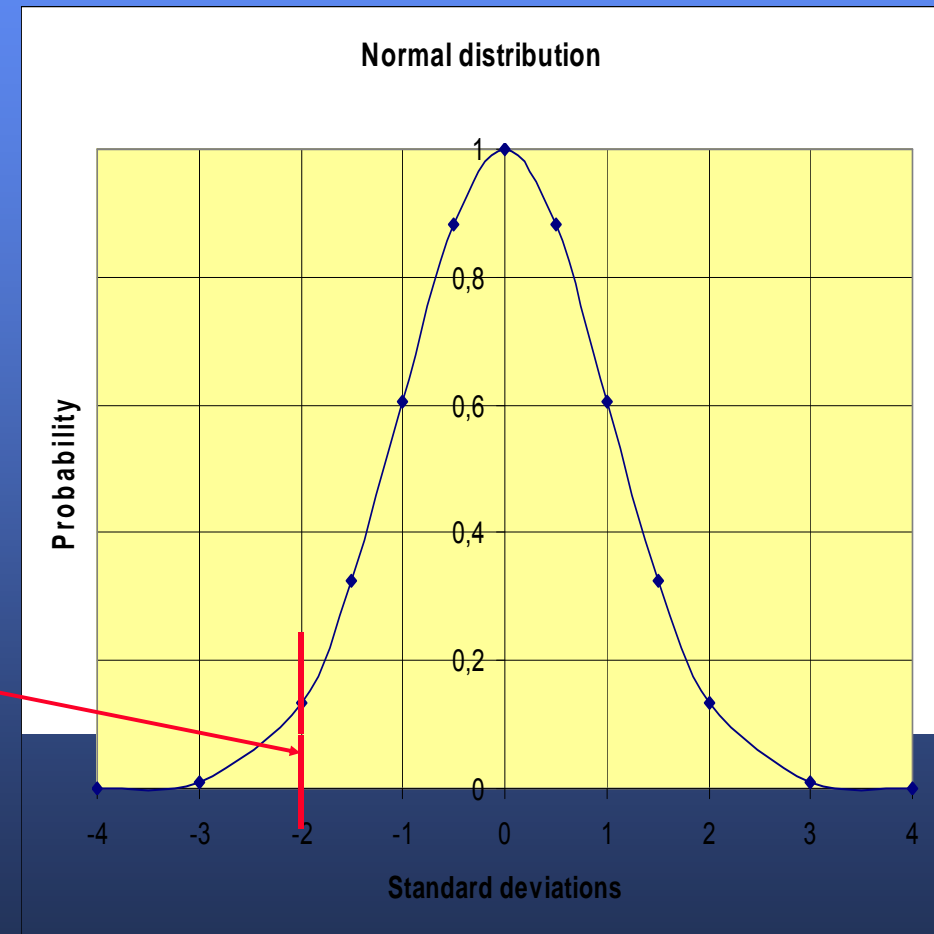



*Cameron Quality Consultants*

## **Performance qualification ..cont..**

### Repetition of Performance Qualification Dose Map

- Only needed if product is changed  
or if facility is changed
- review of documentation at defined intervals is recommended,  
e.g. every 3 years

RISO

