



## Direct, label-free detection of Ochratoxin A using Ag capped silicon nanopillar SERS substrates

Rostami, Simindokht; Zhai, Demi Shuang; Zor, Kinga; Schmidt, Michael Stenbæk; Rindzevicius, Tomas; Smedsgaard, Jørn; Boisen, Anja

*Publication date:*  
2018

*Document Version*  
Peer reviewed version

[Link back to DTU Orbit](#)

*Citation (APA):*  
Rostami, S., Zhai, D. S., Zor, K., Schmidt, M. S., Rindzevicius, T., Smedsgaard, J., & Boisen, A. (2018). *Direct, label-free detection of Ochratoxin A using Ag capped silicon nanopillar SERS substrates*. Abstract from Biosensors 2018, Miami, Florida, United States.

---

### General rights

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from the public portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain
- You may freely distribute the URL identifying the publication in the public portal

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.

**Title:****Direct, label-free detection of Ochratoxin A using Ag capped silicon nanopillar SERS substrates****Authors & affiliations:**

*Simindokht Rostami<sup>1,2,\*</sup>, Demi Shuang Zhai<sup>3</sup>, Kinga Zór<sup>1</sup>, Michael S. Schmidt<sup>1</sup>, Tomas Rindzevicius<sup>1</sup>, Jørn Smedsgaard<sup>3</sup>, Anja Boisen<sup>1</sup>*

<sup>1</sup>*Center for Intelligent Drug Delivery and Sensing Using Microcontainers and Nanomechanics (IDUN), Department of Micro- and Nanotechnology, Technical University of Denmark, Kgs. Lyngby, Denmark*

<sup>2</sup>*Department of Chemistry, Faculty of Science, K. N. Toosi University of Technology, Tehran, Iran*

<sup>3</sup>*National Food Institute Technical University of Denmark, Kgs. Lyngby*  
*\*simros@nanotech.dtu.dk*

**Abstract:**

Ochratoxin A (OTA), one of the mycotoxins secreted by fungi species, is classified as one of the most significant contaminants in various foods like coffee, maize, cereals and cereal products, milk, wine and beer. It is known that OTA is responsible for several diseases, hence different analytical methods including HPLC, LC-MS, immunoassays, spectroscopic methods and biosensors have been developed to quantify OTA in samples [1].

Various, nanostructure based biosensors have been evaluated to create new analytical probes for detection of OTA [2] in combination with OTA aptamers has been widely used to improve sensitivity and accuracy [3]. However, in order to simplify the detection procedure and to avoid the need for complex surface modification procedures, there is still the need to develop simple and fast methods for direct detection of OTA.

Surface-enhanced Raman Spectroscopy (SERS) is a powerful technique, capable of direct detection of molecular fingerprints of analytes, with high sensitivity and fast response time [4], therefore it would be a suitable approach for the detection of OTA. Based on our knowledge, the SERS sensors for OTA used today, are mostly based on modifying the surfaces with recognition elements such as antibodies or aptamers in order to increase sensitivity.

In this work, we present the development of a SERS sensor for direct, label-free detection of OTA using silver capped silicon nanopillars, as a low cost, fast and easy-to-use approach (Fig.1). To evaluate the dynamic range of the sensor, OTA was measured in ethanol (EtOH). The calibration curve of OTA in EtOH showed that OTA can be detected in a linear range between 0.02 to 1ppm OTA (Fig.2).

The here presented results shows a promising path for quick and portable detection of OTA in remote rural areas. In addition, the sensor sensitivity can be further improved by using the electrochemical-SERS concept.

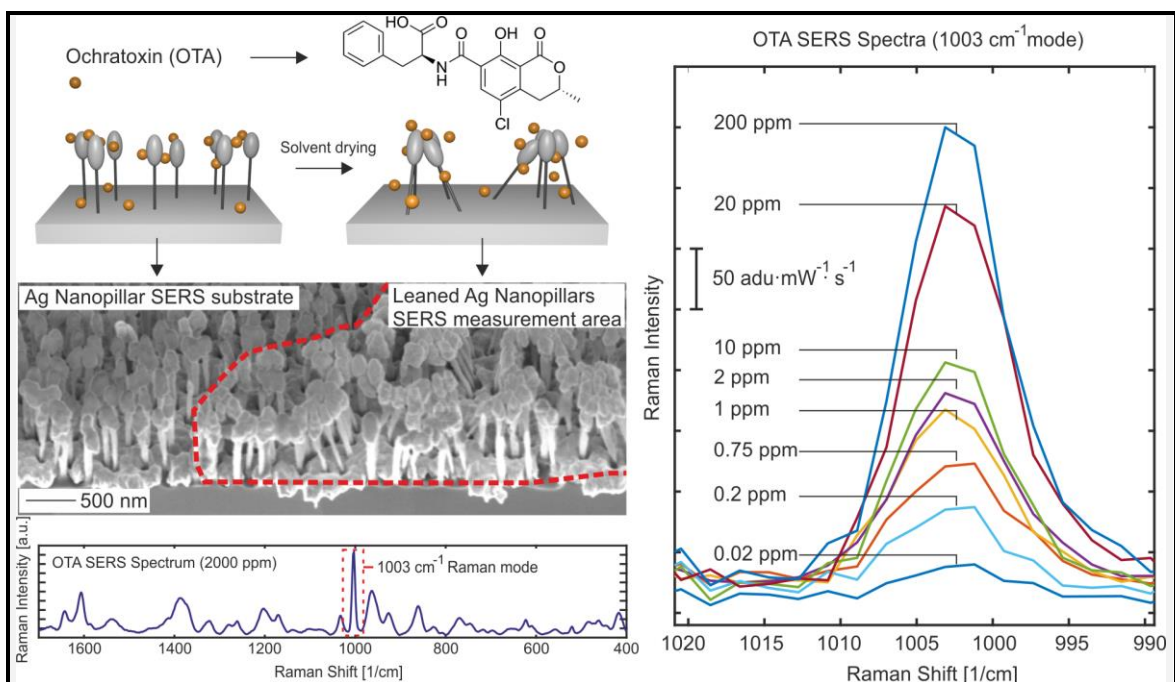


Figure 1. Left: schematic illustration of direct ochratoxin A (OTA) detection using leaning silver-capped silicon nanopillar SERS substrates. A droplet of analyte is deposited on the SERS substrate and left for drying. As a solvent evaporates, capillary forces pull the Ag nanopillars together forming microsized clusters that trap OTA in the vicinity of the so-called electromagnetic “hot spots” [5]. Right: OTA SERS spectra for concentrations varying between 0.02 – 200 ppm (1003  $\text{cm}^{-1}$  mode).

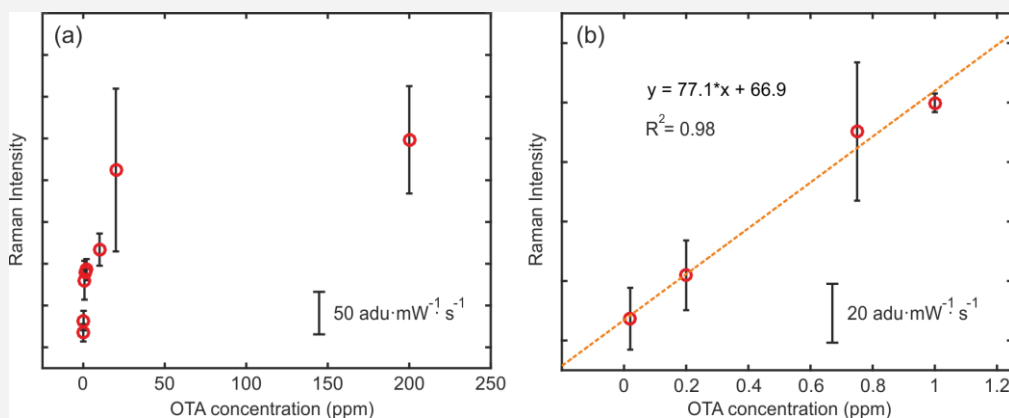


Figure 2. A) Plot of Raman intensity as a function of OTA concentration. The error bar represents the standard deviation of triplicate measurement; B) Calibration curve of OTA in concentration range between 0.02 to 1 ppm; error bar represents the standard deviations of three independent measurements.

#### References:

- [1] T.H Ha, *Toxins*, 2015, 7, 5276-5300.
- [2] E-O Ganbold, et al, *Analytical Methods*, 2014, 6, 3573-3577.
- [3] JA. Cruz-Aguado, G. Penner, *Journal of Agricultural Food Chemistry*, 2008, 56, 10456-10461.
- [4] W. E. Smith, *Chemical Society Reviews*, 2008, 37, 955–964.
- [5] K. Wu, et al, *Journal of Physical Chemistry C*, 2015, 119, 2053-2062.